

**SCOTTSDALE POLICE DEPARTMENT CRIME LABORATORY  
BLOOD ALCOHOL FACE SHEET**

ANALYSIS DATE 09/22/2026

SEQUENCE NAME 22Sep26

**EQUIPMENT**

Pipettor

☐

Hamilton ML600GA23163

☒

Hamilton ML600EA20851

Gas Chromatograph

☐

Agilent US14173023

☒

Agilent US2320A024

**INSTRUMENT CALIBRATION**

Vial 1 0.02 calibrator Lot FN03062417

Coefficient of determination ( $r^2$ )

Vial 2 0.10 calibrator Lot FN04042408

1.00000

Vial 3 0.20 calibrator Lot FN07012401

Vial 4 0.40 calibrator Lot FN04092401

**CALIBRATION VERIFICATION AND RESOLUTION TEST**

| Vial | Sample            | Expected result | Measured result | Manufacturer/lot  |
|------|-------------------|-----------------|-----------------|-------------------|
| 5    | Blank             | Not detected    | Not detected    | SPD lab 061026WB  |
| 6    | Volatiles mixture | 6 compounds     | 6 compounds     | SPD lab 031726MIX |
| 7    | Aqueous control   | 0.400 g/dL      | 0.402 g/dL      | Lipomed C1475825  |
| 8    | Aqueous control   | 0.040 g/dL      | 0.040 g/dL      | Lipomed 529862    |
| 9    | Blood control     | 0.201 d/dL      | 0.202 g/dL      | ACQ 4190425129/6  |
| 20   | Aqueous control   | 0.100 g/dL      | 0.101 g/dL      | Lipomed C1455120  |
| 31   | Blood control     | 0.201 g/dL      | 0.204 g/dL      | ACQ 4190425129/6  |
| 42   | Aqueous control   | 0.100 g/dL      | 0.101 g/dL      | Lipomed C1455120  |
| 53   | Blood control     | 0.201 g/dL      | 0.204 g/dL      | ACQ 4190425129/6  |
| 64   | Aqueous control   | 0.400 g/dL      | 0.408 g/dL      | Lipomed C1475825  |
| 65   | Aqueous control   | 0.040 g/dL      | 0.040 g/dL      | Lipomed 529862    |
| 66   | Blood control     | 0.201 g/dL      | 0.201 g/dL      | ACQ 4190425129/6  |
| 67   | Blank             | Not detected    | Not detected    | SPD lab 073126AQ  |
|      |                   |                 |                 |                   |
|      |                   |                 |                 |                   |
|      |                   |                 |                 |                   |
|      |                   |                 |                 |                   |


**SUBJECT SAMPLES**

Subjects in the sequence 25

Subjects requiring reanalysis 3

**ADDITIONAL NOTES:** Subject samples in vials #25, #34, and #37 failed to inject.

Run valid ☒  
Run invalid ☐

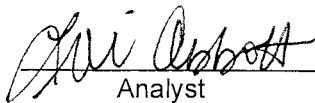  
Analyst

Run valid ☒  
Run invalid ☐

  
Technical Reviewer

# Scottsdale Police Department Crime Laboratory


## Sequence Quality Assurance Summary

Sequence: 22Sep26

Analyst: Abbott *du*

| Sample Name        | Vial | Measured Value (g/dL) | Expected Value (g/dL) | Percent Difference | Absolute Difference (g/dL) |
|--------------------|------|-----------------------|-----------------------|--------------------|----------------------------|
| blank 061026WB     | 5    | negative              | negative              | -                  | -                          |
| 0.400 C1475825     | 7    | 0.402                 | 0.400                 | 0.50               | 0.002                      |
| 0.040 529862       | 8    | 0.040                 | 0.040                 | 0.00               | 0.000                      |
| 0.201 4190425129/6 | 9    | 0.202                 | 0.201                 | 0.50               | 0.001                      |
| 0.100 C1455120     | 20   | 0.101                 | 0.100                 | 1.00               | 0.001                      |
| 0.201 4190425129/6 | 31   | 0.204                 | 0.201                 | 1.49               | 0.003                      |
| 0.100 C1455120     | 42   | 0.101                 | 0.100                 | 1.00               | 0.001                      |
| 0.201 4190425129/6 | 53   | 0.204                 | 0.201                 | 1.49               | 0.003                      |
| 0.400 C1475825     | 64   | 0.408                 | 0.400                 | 2.00               | 0.008                      |
| 0.040 529862       | 65   | 0.040                 | 0.040                 | 0.00               | 0.000                      |
| 0.201 4190425129/6 | 66   | 0.201                 | 0.201                 | 0.00               | 0.000                      |
| blank 073126AQ     | 67   | negative              | negative              | -                  | -                          |

| Calibrator | Ethanol Area | N-propanol Area | Ratio |
|------------|--------------|-----------------|-------|
| 0.020      | 13.736       | 189.071         | 0.073 |
| 0.100      | 71.233       | 193.037         | 0.369 |
| 0.200      | 141.240      | 190.809         | 0.740 |
| 0.400      | 280.038      | 188.684         | 1.484 |

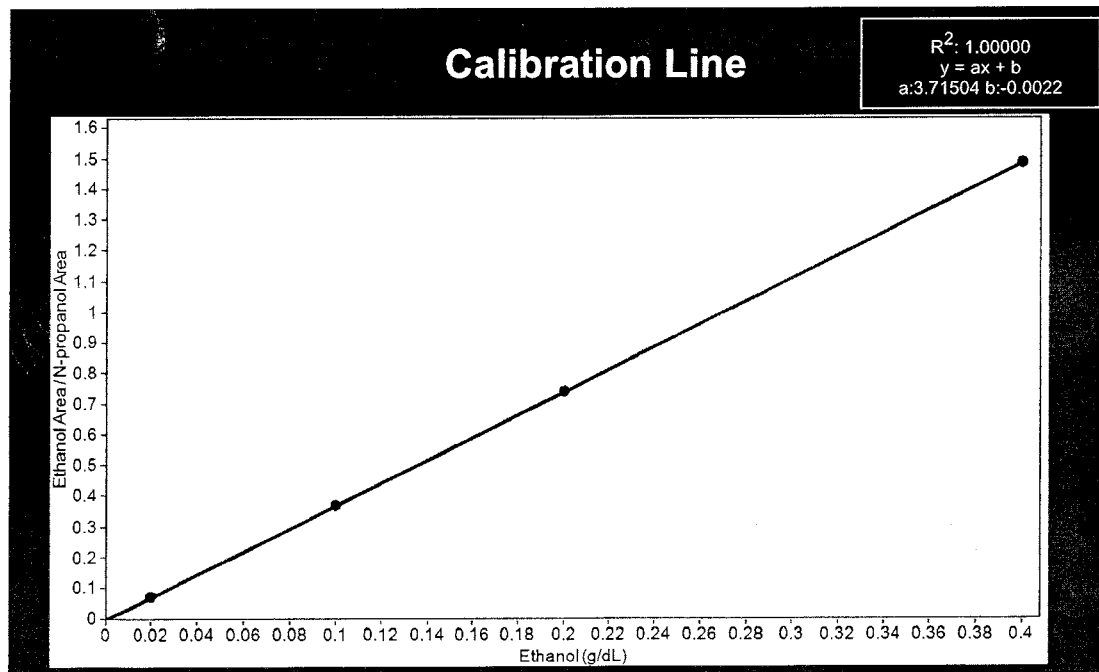

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | FN03062417                   | <b>Vial:</b>           | 1                    |
| <b>Description:</b> | 0.020 Calibrator             | <b>Sequence:</b>       | 22Sep26              |
| <b>Method:</b>      | ethanol quant.M              | <b>Injection date:</b> | 9/22/2026 3:13:35 PM |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Analyst:</b>        | Abbott <i>du</i>     |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\1.D |                        |                      |

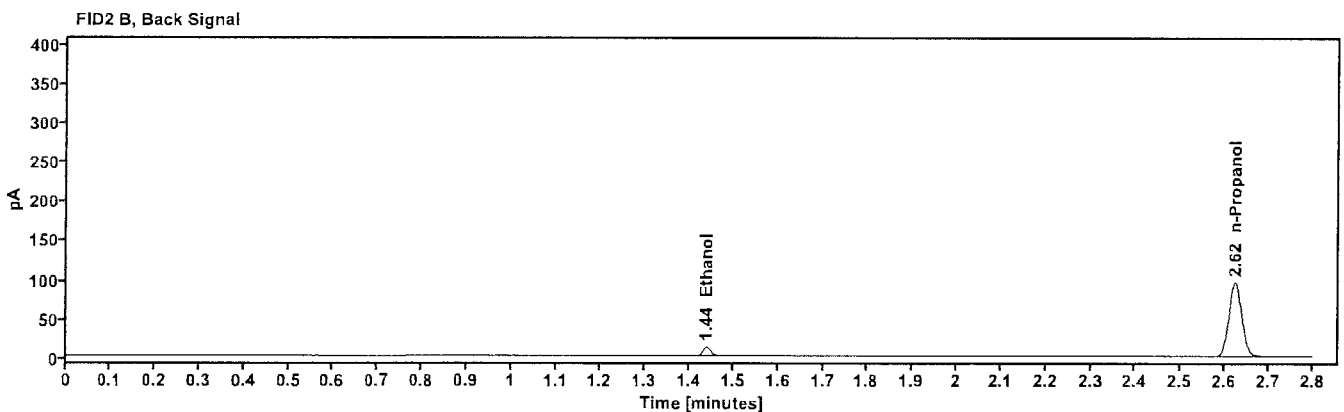
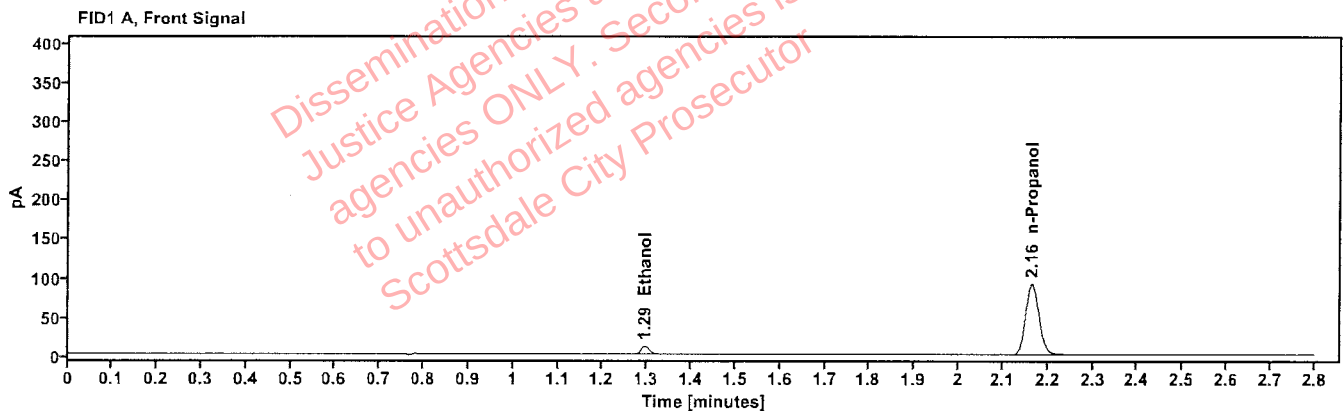


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.295      | 13.736    |
| n-Propanol | 2.164      | 189.071   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 14.276    |
| n-Propanol | 2.624      | 196.254   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | FN04042408                   | <b>Vial:</b>           | 2                    |
| <b>Description:</b> | 0.100 Calibrator             | <b>Sequence:</b>       | 22Sep26              |
| <b>Method:</b>      | ethanol quant.M              | <b>Injection date:</b> | 9/22/2026 3:17:23 PM |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Analyst:</b>        | Abbott <i>DA</i>     |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\2.D |                        |                      |

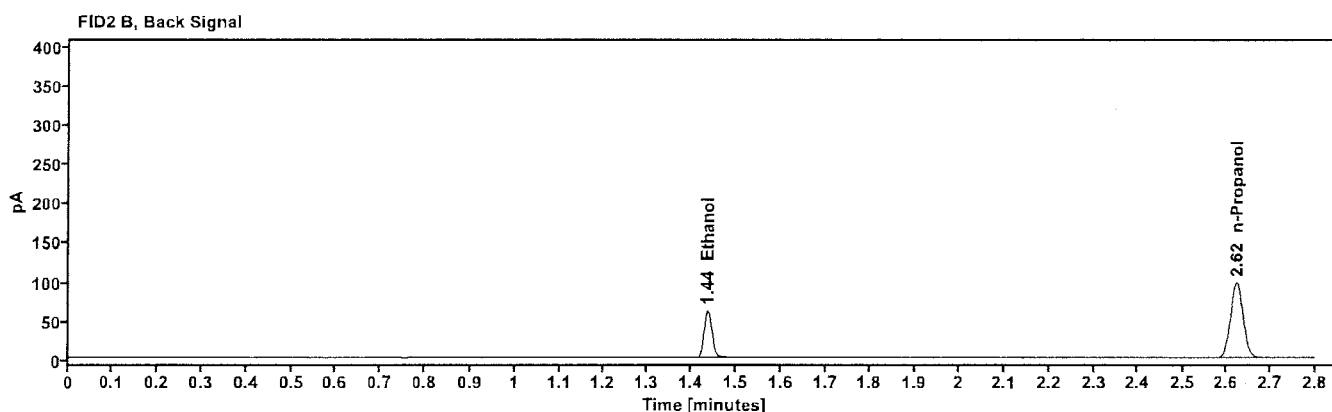
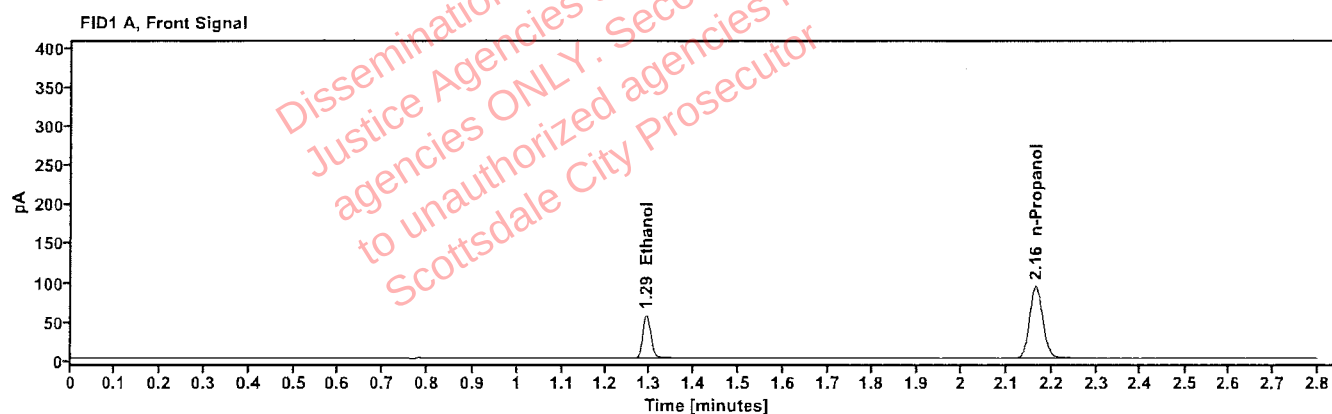


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.292      | 71.233    |
| n-Propanol | 2.164      | 193.037   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 73.771    |
| n-Propanol | 2.624      | 200.305   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                              |                 |                      |
|--------------|------------------------------|-----------------|----------------------|
| Sample:      | FN07012401                   | Vial:           | 3                    |
| Description: | 0.200 Calibrator             | Sequence:       | 22Sep26              |
| Method:      | ethanol quant.M              | Injection date: | 9/22/2026 3:21:11 PM |
| Instrument:  | US2320A024 US23180007        | Analyst:        | Abbott <i>lu</i>     |
| Data file:   | C:\Chem32\1\Data\22Sep26\3.D |                 |                      |

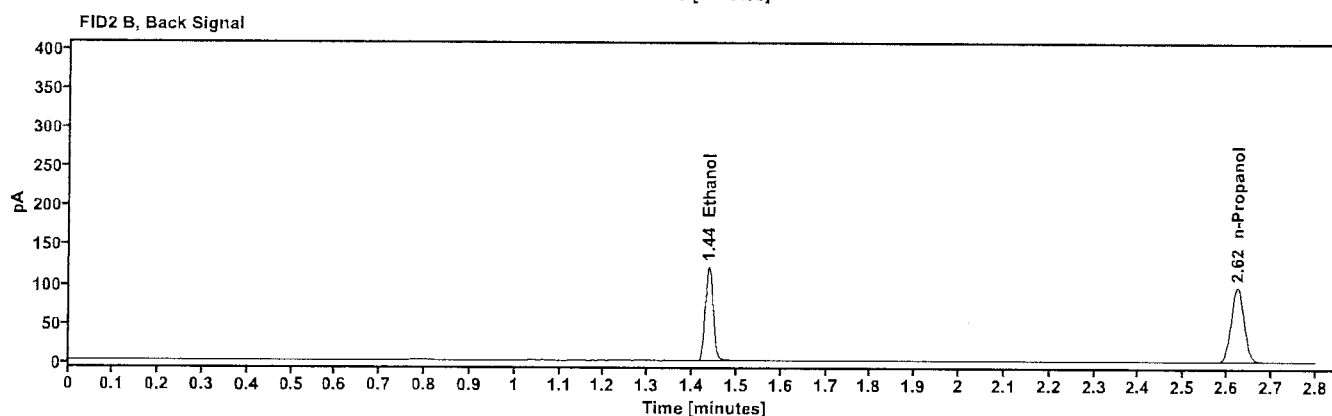
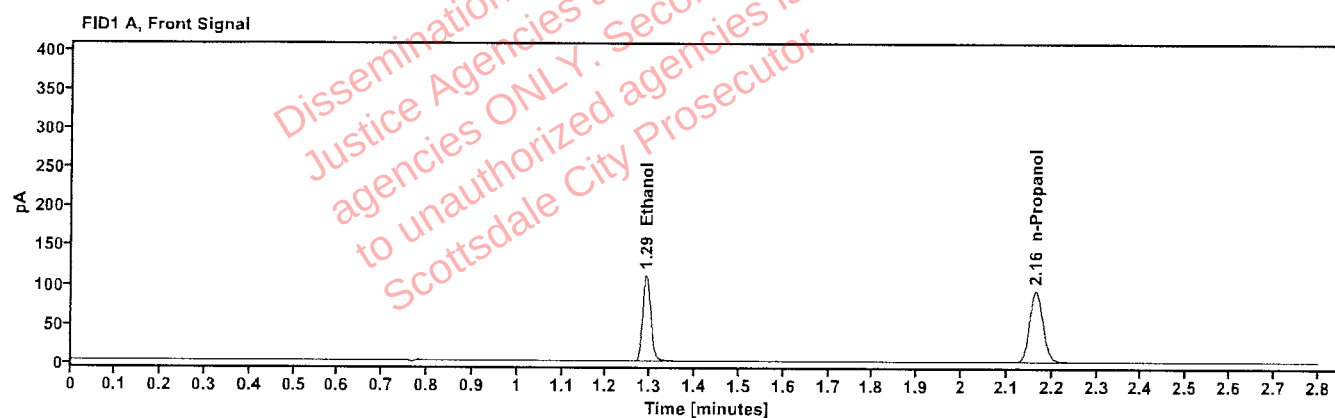


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.291      | 141.240   |
| n-Propanol | 2.164      | 190.809   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 146.418   |
| n-Propanol | 2.624      | 197.959   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | FN04092401                   | <b>Vial:</b>           | 4                    |
| <b>Description:</b> | 0.400 Calibrator             | <b>Sequence:</b>       | 22Sep26              |
| <b>Method:</b>      | ethanol quant.M              | <b>Injection date:</b> | 9/22/2026 3:24:59 PM |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Analyst:</b>        | Abbott <i>dn</i>     |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\4.D |                        |                      |

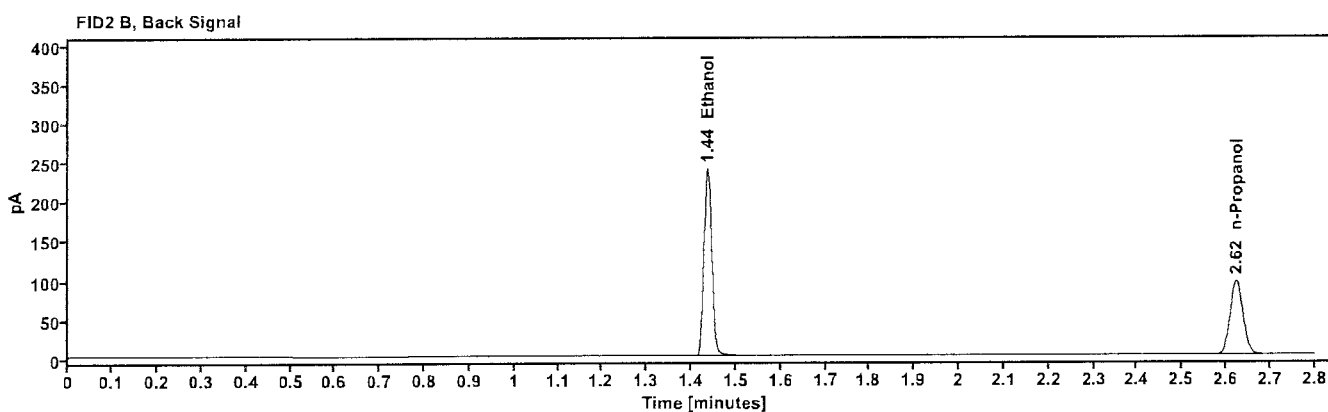
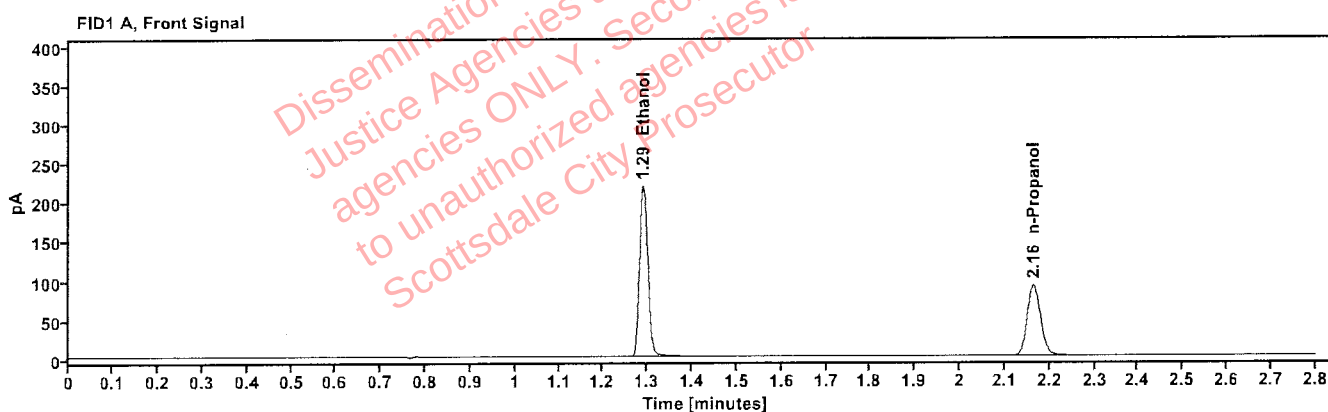


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.290      | 280.038   |
| n-Propanol | 2.164      | 188.684   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.436      | 289.946   |
| n-Propanol | 2.624      | 195.717   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 061026WB                     | <b>Vial:</b>           | 5                    |
| <b>Description:</b> | negative                     | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M              | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Injection date:</b> | 9/22/2026 3:28:47 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\5.D | <b>Analyst:</b>        | Abbott               |

*[Signature]*

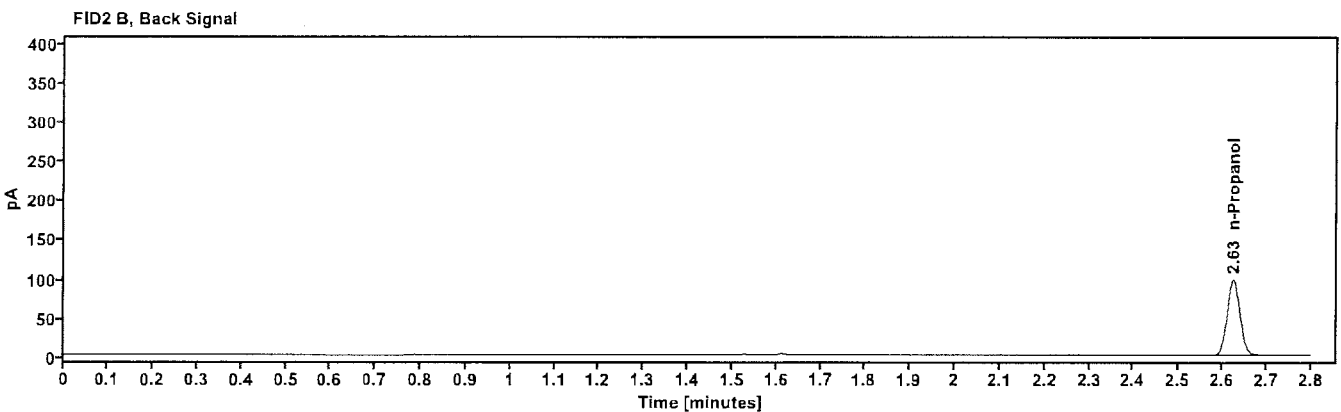
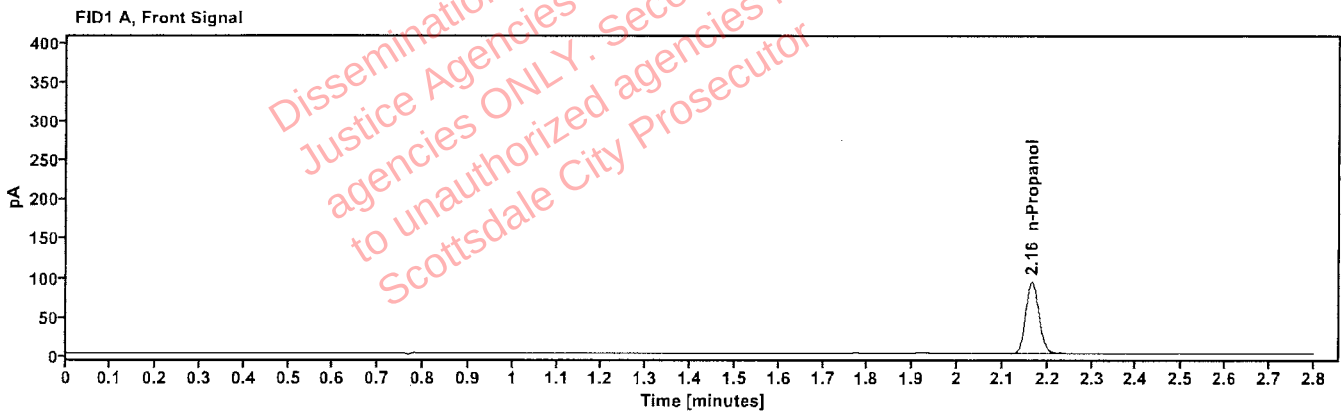


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Amount (g/100mL) | Time (min) | Peak Area |
|------------|------------------|------------|-----------|
| n-Propanol | -----            | 2.165      | 193.123   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| n-Propanol | 2.625      | 200.423   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 031726MIX                    | <b>Vial:</b>           | 6                    |
| <b>Description:</b> | Volatile Mix                 | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M              | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Injection date:</b> | 9/22/2026 3:32:35 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\6.D | <b>Analyst:</b>        | Abbott               |

*[Signature]*

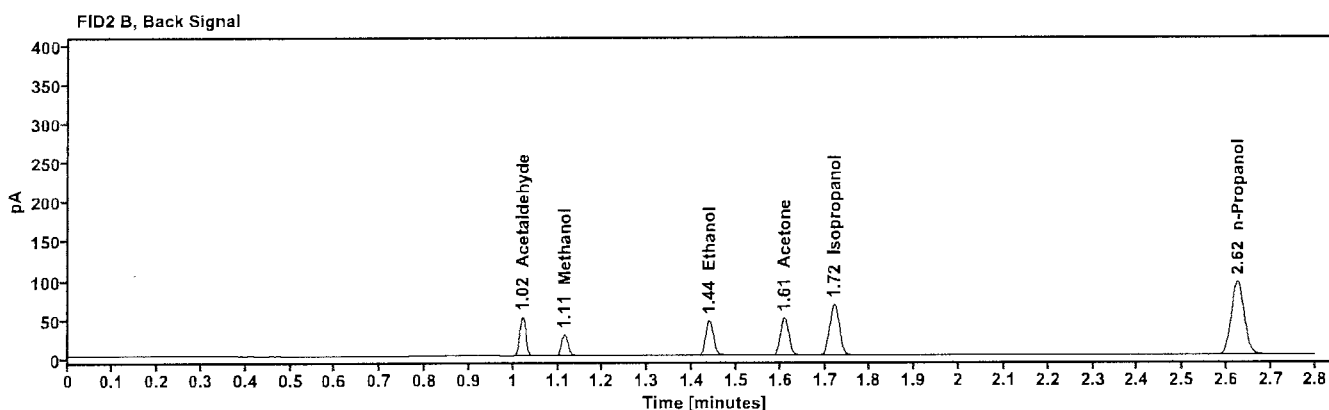
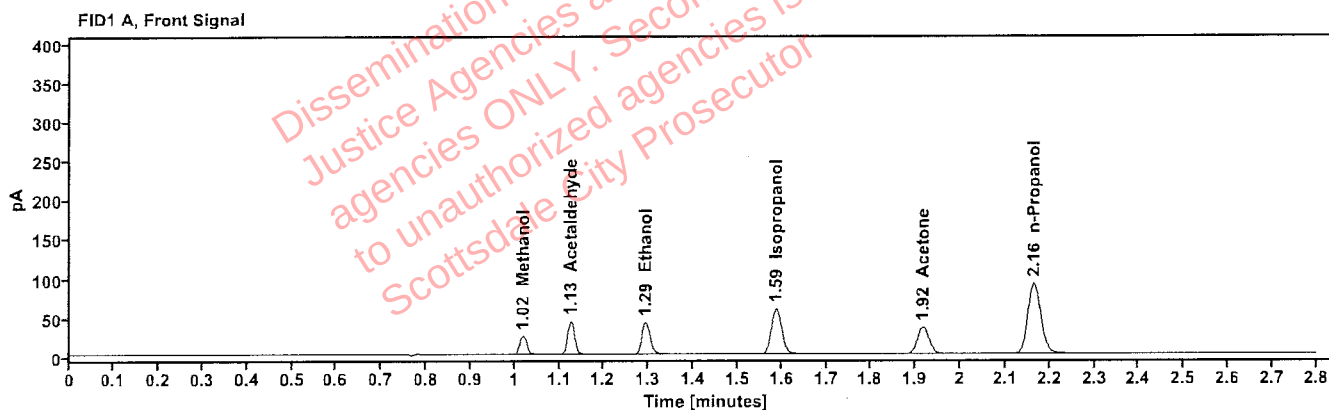


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| Methanol           | -----            | 1.019      | 26.883    |
| Acetaldehyde       | -----            | 1.126      | 46.174    |
| quantified Ethanol | 0.0765           | 1.293      | 53.115    |
| Isopropanol        | -----            | 1.588      | 96.178    |
| Acetone            | -----            | 1.916      | 62.365    |
| n-Propanol         | -----            | 2.164      | 188.234   |

Table 2: FID 2 B (column DB-ALC2)

| Compound     | Time (min) | Peak Area |
|--------------|------------|-----------|
| Acetaldehyde | 1.021      | 47.866    |
| Methanol     | 1.115      | 27.632    |
| Ethanol      | 1.438      | 55.094    |
| Acetone      | 1.609      | 64.608    |
| Isopropanol  | 1.721      | 99.754    |
| n-Propanol   | 2.625      | 195.167   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                              |                 |                      |
|--------------|------------------------------|-----------------|----------------------|
| Sample:      | C1475825                     | Vial:           | 7                    |
| Description: | 0.400                        | LIMS ID:        |                      |
| Method:      | ethanol quant.M              | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007        | Injection date: | 9/22/2026 3:37:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\7.D | Analyst:        | Abbott               |

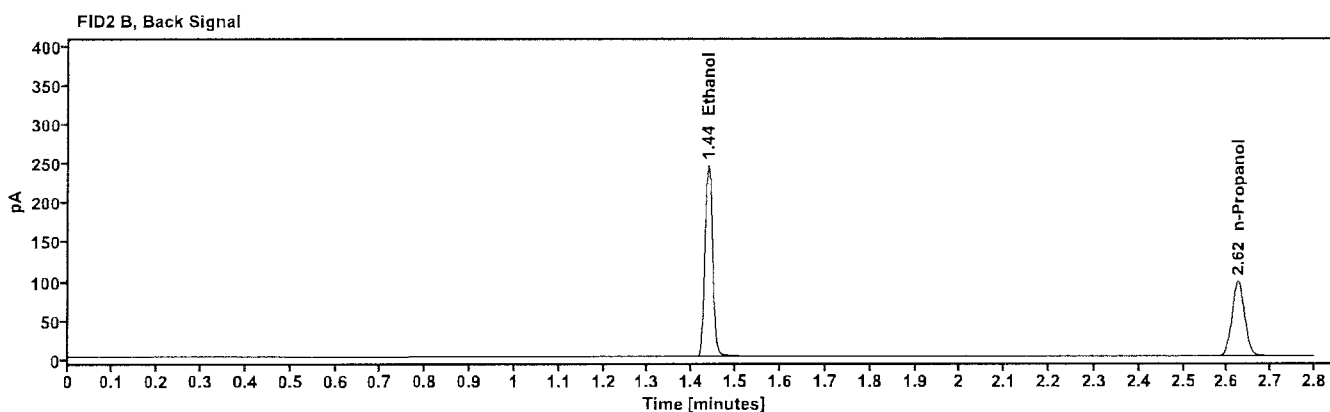
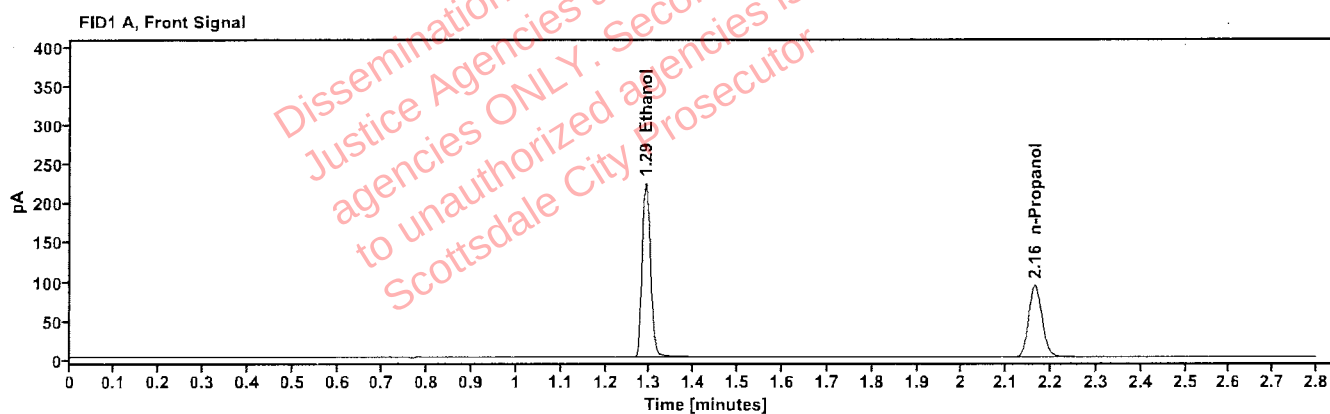


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.4027           | 1.291      | 286.019   |
| n-Propanol         | -----            | 2.164      | 191.486   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 296.153   |
| n-Propanol | 2.624      | 198.612   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                              |                 |                      |
|--------------|------------------------------|-----------------|----------------------|
| Sample:      | 529862                       | Vial:           | 8                    |
| Description: | 0.040                        | LIMS ID:        |                      |
| Method:      | ethanol quant.M              | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007        | Injection date: | 9/22/2026 3:41:22 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\8.D | Analyst:        | Abbott               |

*Handwritten signature/initials*

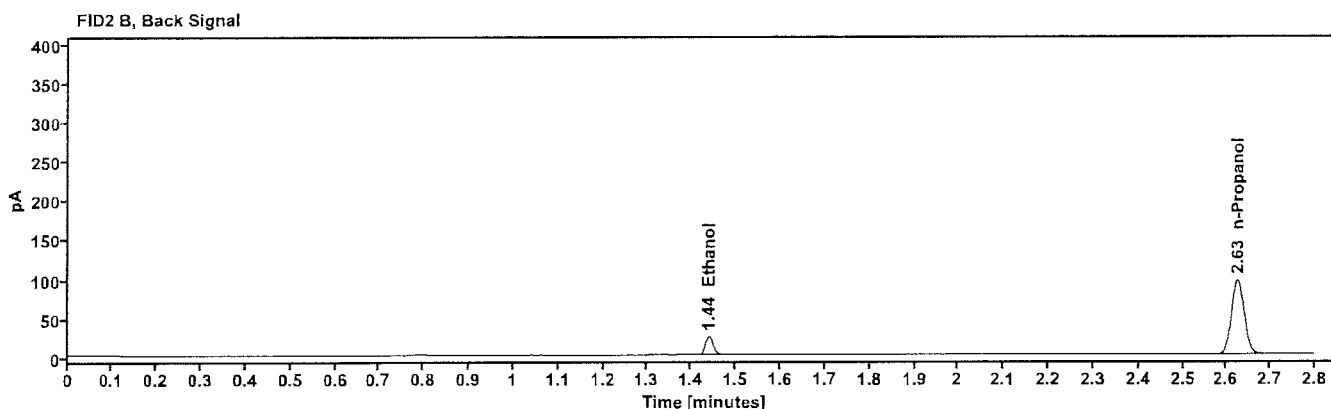
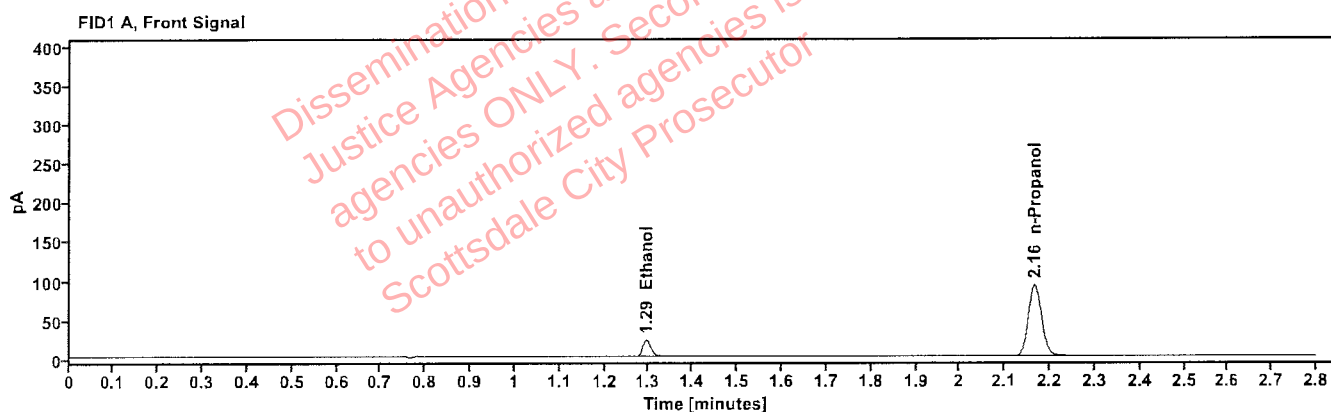


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0401           | 1.294      | 27.915    |
| n-Propanol         | ----             | 2.165      | 190.317   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 29.017    |
| n-Propanol | 2.626      | 197.394   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                              |                        |                      |
|---------------------|------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 4190425129/6                 | <b>Vial:</b>           | 9                    |
| <b>Description:</b> | 0.201                        | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M              | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007        | <b>Injection date:</b> | 9/22/2026 3:45:11 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\9.D | <b>Analyst:</b>        | Abbott               |

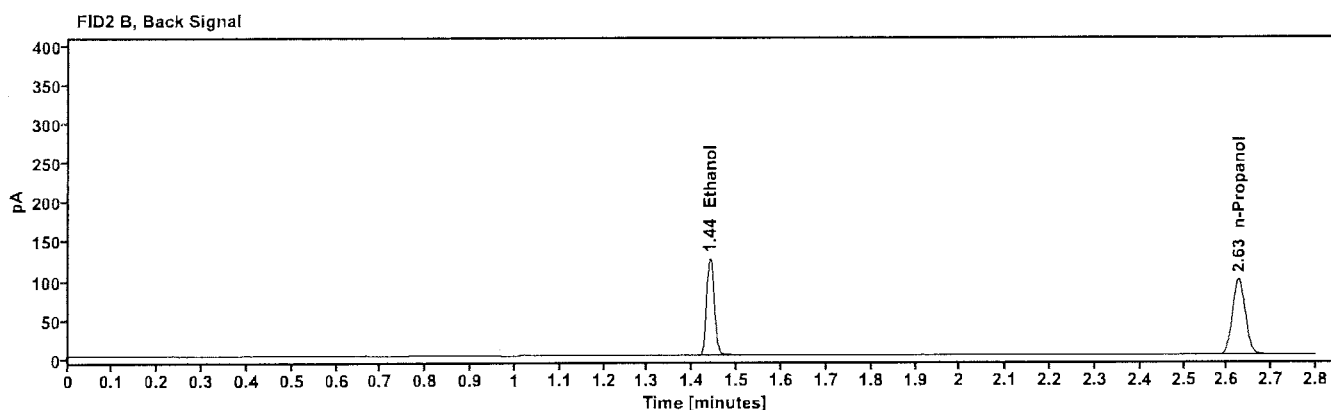
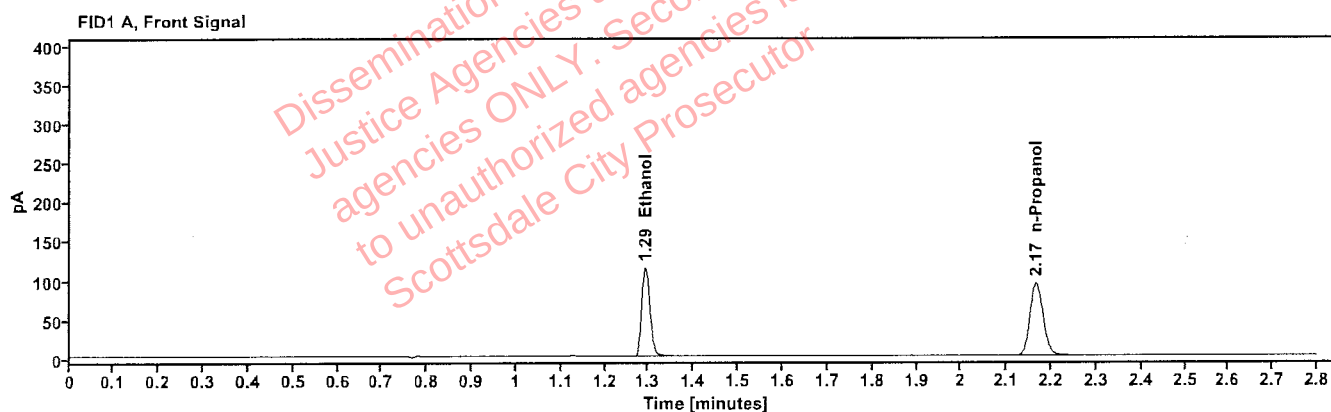


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2029           | 1.292      | 145.610   |
| n-Propanol         | -----            | 2.165      | 193.721   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 150.675   |
| n-Propanol | 2.626      | 200.953   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | C1455120                      | <b>Vial:</b>           | 20                   |
| <b>Description:</b> | 0.100                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 4:29:22 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\20.D | <b>Analyst:</b>        | Abbott               |

*[Handwritten signature]*

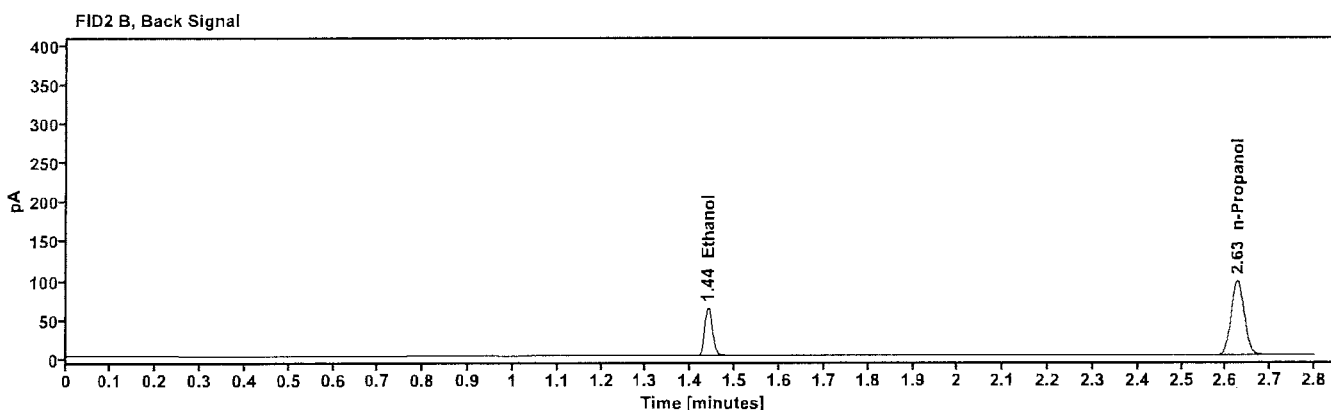
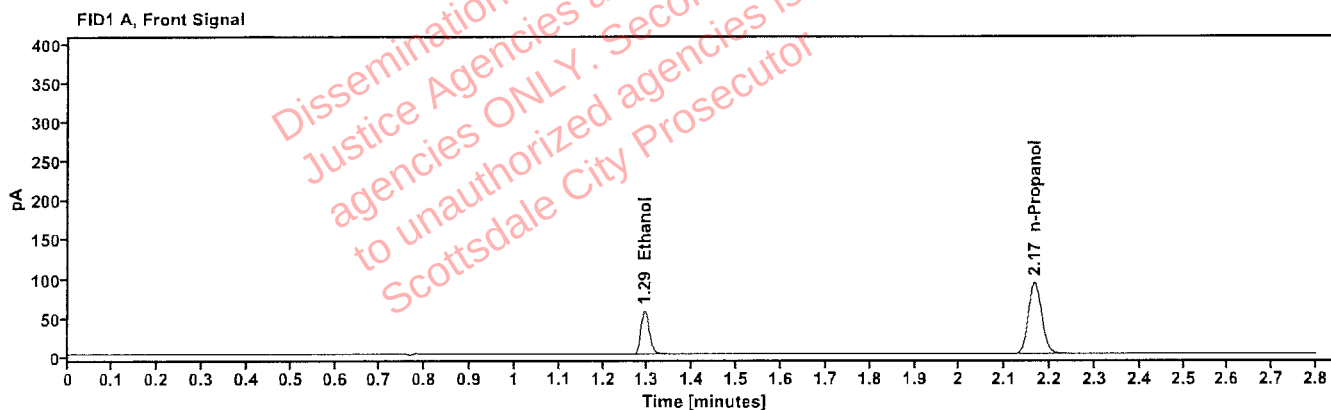


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1019           | 1.293      | 71.977    |
| n-Propanol         | -----            | 2.165      | 191.240   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 74.590    |
| n-Propanol | 2.626      | 198.197   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 4190425129/6                  | <b>Vial:</b>           | 31                   |
| <b>Description:</b> | 0.201                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 5:13:35 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\31.D | <b>Analyst:</b>        | Abbott               |

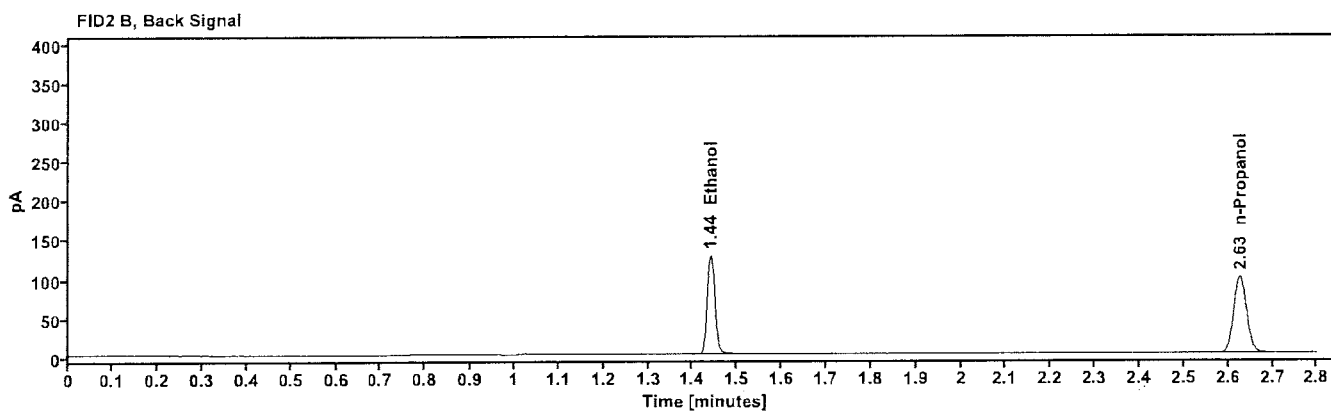
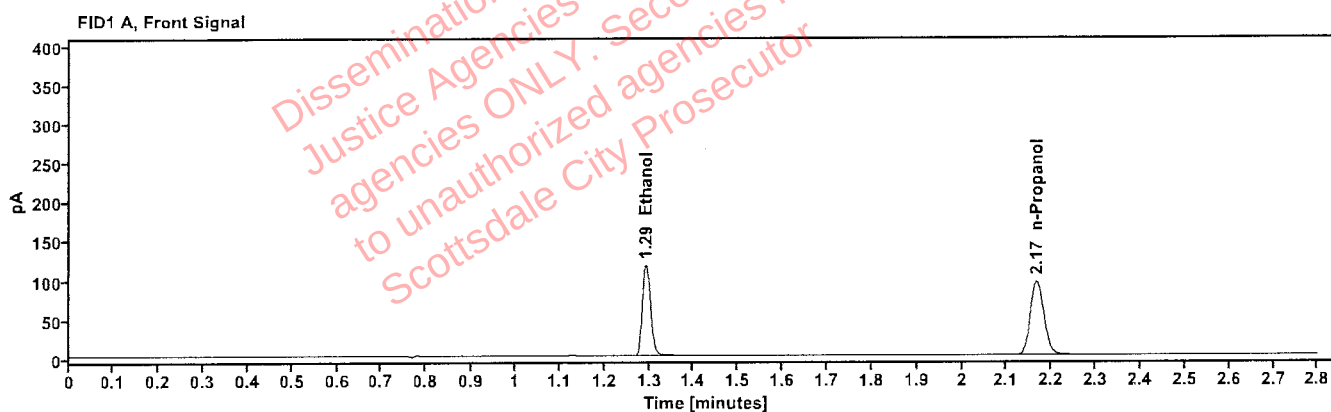


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2042           | 1.293      | 148.823   |
| n-Propanol         | ----             | 2.166      | 196.780   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 154.272   |
| n-Propanol | 2.627      | 204.074   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | C1455120                      | <b>Vial:</b>           | 42                   |
| <b>Description:</b> | 0.100                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 5:56:35 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\42.D | <b>Analyst:</b>        | Abbott               |

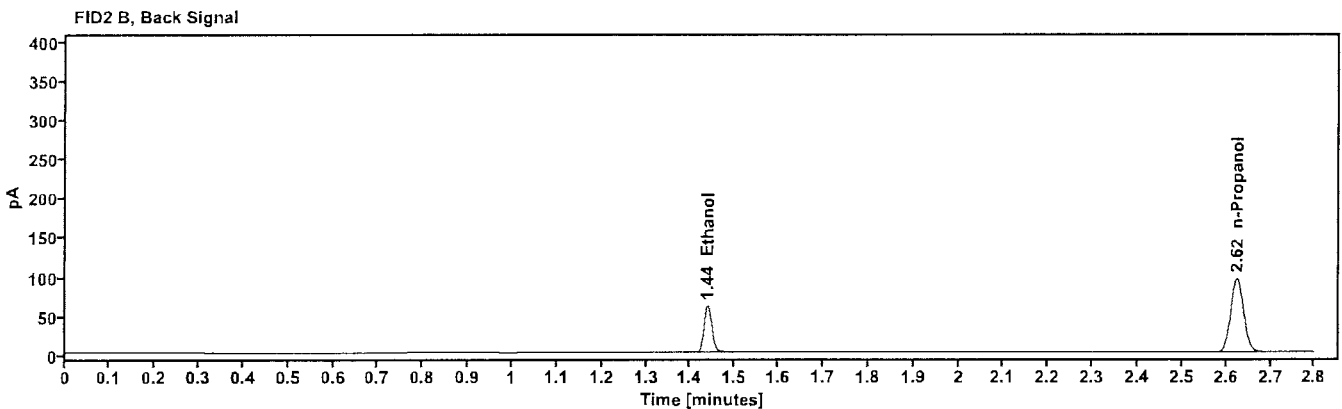
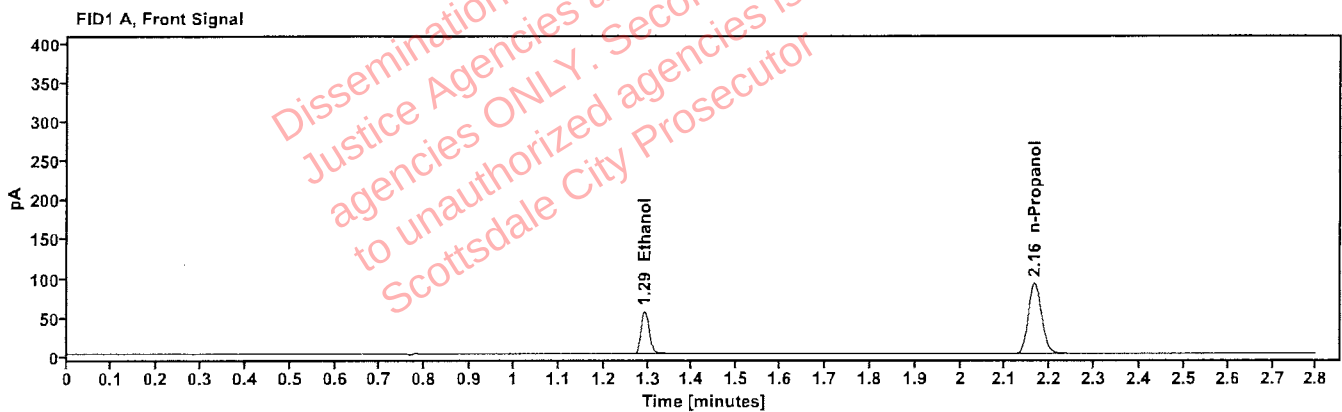


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1014           | 1.292      | 70.606    |
| n-Propanol         | -----            | 2.164      | 188.476   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.438      | 73.075    |
| n-Propanol | 2.624      | 195.452   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 4190425129/6                  | <b>Vial:</b>           | 53                   |
| <b>Description:</b> | 0.201                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 6:40:46 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\53.D | <b>Analyst:</b>        | Abbott               |

*[Handwritten signature]*

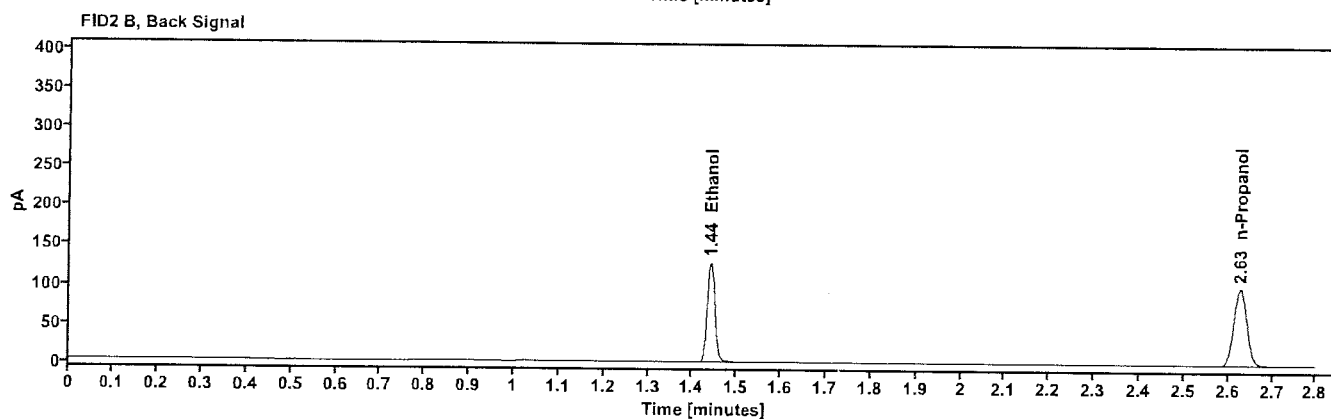
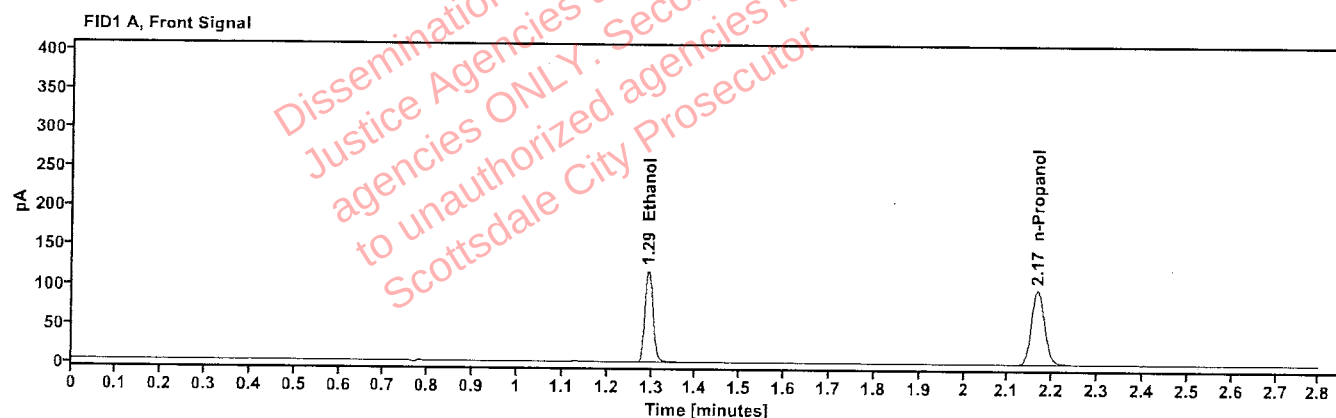


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2043           | 1.293      | 149.839   |
| n-Propanol         | -----            | 2.166      | 197.952   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 154.823   |
| n-Propanol | 2.627      | 205.208   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | C1475825                      | <b>Vial:</b>           | 64                   |
| <b>Description:</b> | 0.400                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 7:24:59 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\64.D | <b>Analyst:</b>        | Abbott               |

*[Signature]*

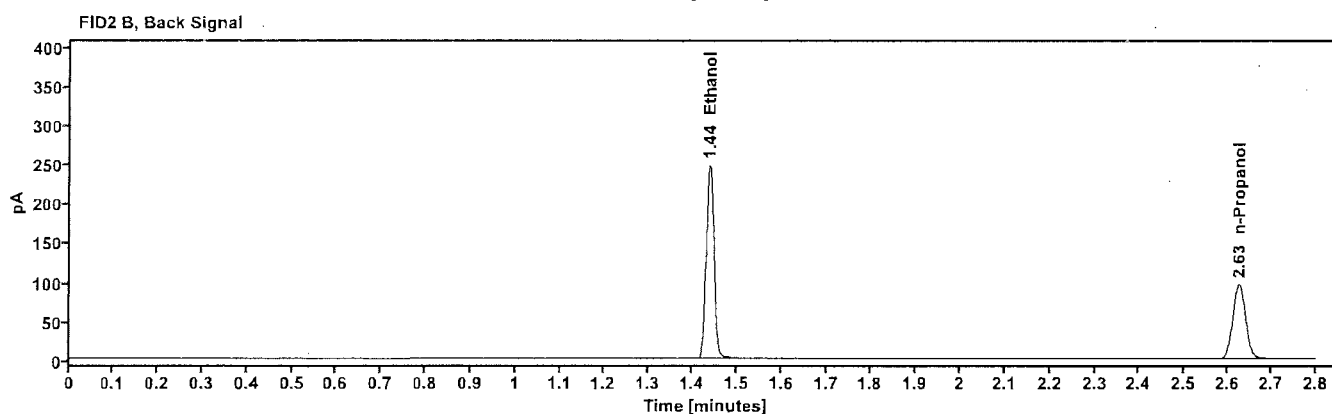
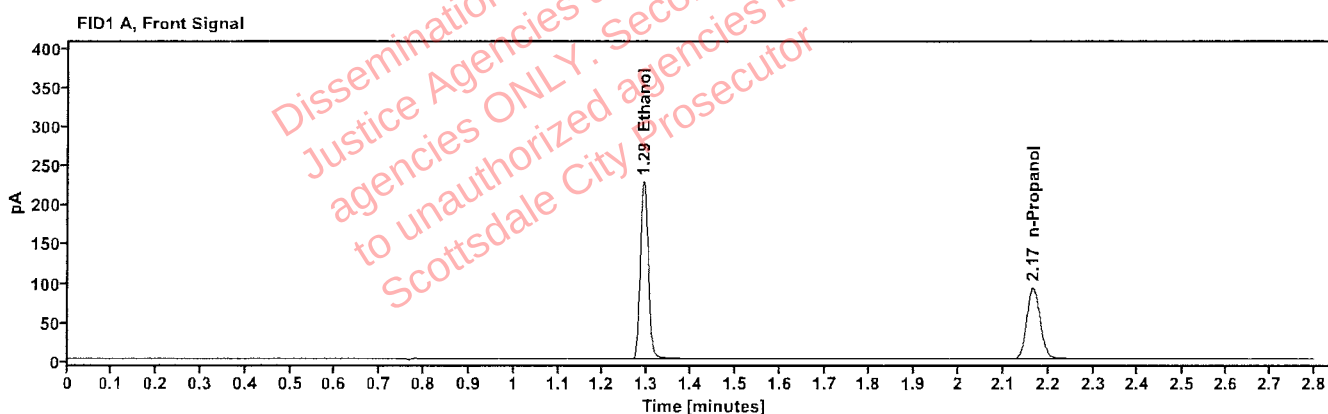


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.4084           | 1.292      | 289.840   |
| n-Propanol         | -----            | 2.165      | 191.320   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 299.961   |
| n-Propanol | 2.626      | 198.105   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|                     |                               |                        |                      |
|---------------------|-------------------------------|------------------------|----------------------|
| <b>Sample:</b>      | 529862                        | <b>Vial:</b>           | 65                   |
| <b>Description:</b> | 0.040                         | <b>LIMS ID:</b>        |                      |
| <b>Method:</b>      | ethanol quant.M               | <b>Sequence:</b>       | 22Sep26              |
| <b>Instrument:</b>  | US2320A024 US23180007         | <b>Injection date:</b> | 9/22/2026 7:28:47 PM |
| <b>Data file:</b>   | C:\Chem32\1\Data\22Sep26\65.D | <b>Analyst:</b>        | Abbott <i>JA</i>     |

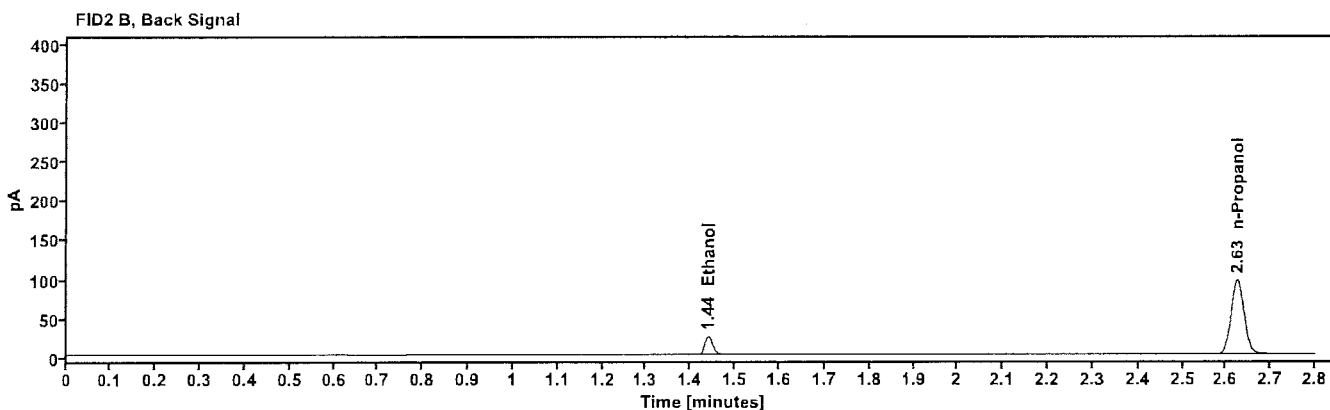
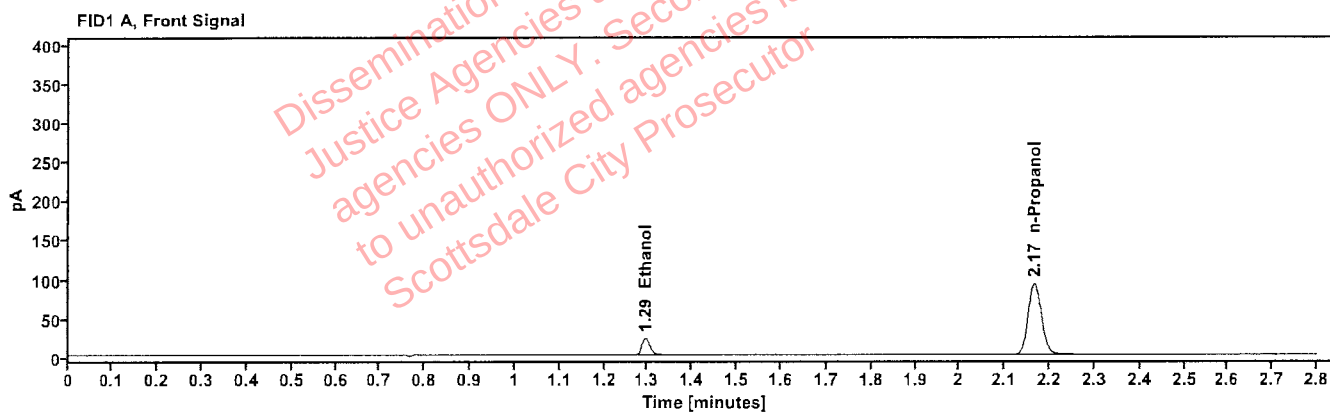


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0403           | 1.295      | 28.120    |
| n-Propanol         | ----             | 2.165      | 190.439   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 29.243    |
| n-Propanol | 2.626      | 197.359   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 4190425129/6                  | Vial:           | 66                   |
| Description: | 0.201                         | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:32:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\66.D | Analyst:        | Abbott               |

sh

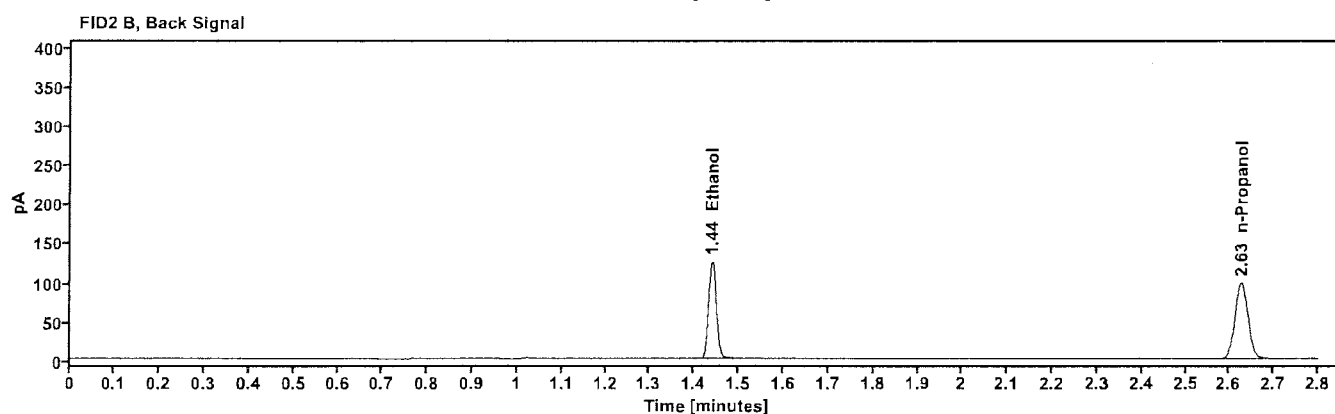
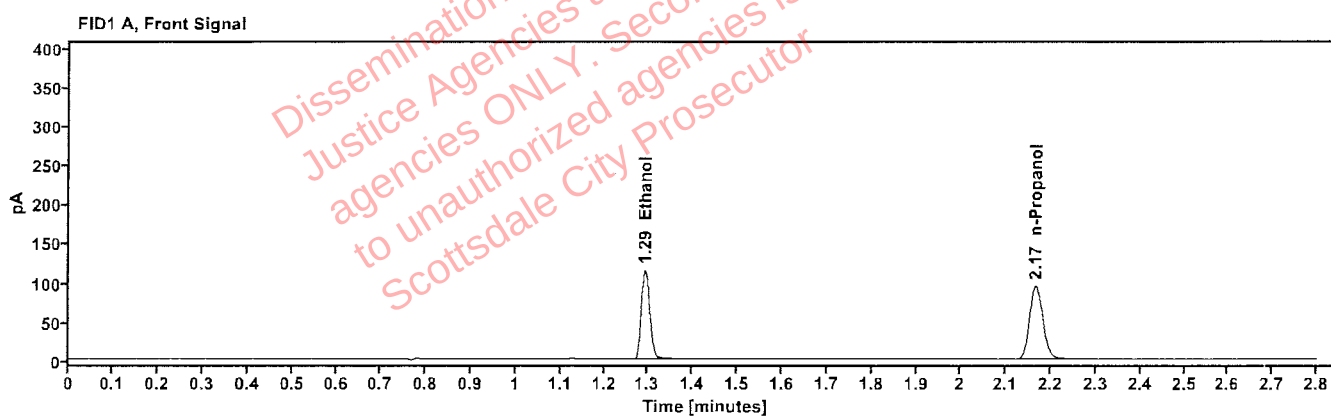


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2019           | 1.293      | 146.451   |
| n-Propanol         | —                | 2.166      | 195.849   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 151.453   |
| n-Propanol | 2.627      | 202.940   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 073126AQ                      | Vial:           | 67                   |
| Description: | negative                      | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:37:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\67.D | Analyst:        | Abbott               |

*Handwritten signature/initials*

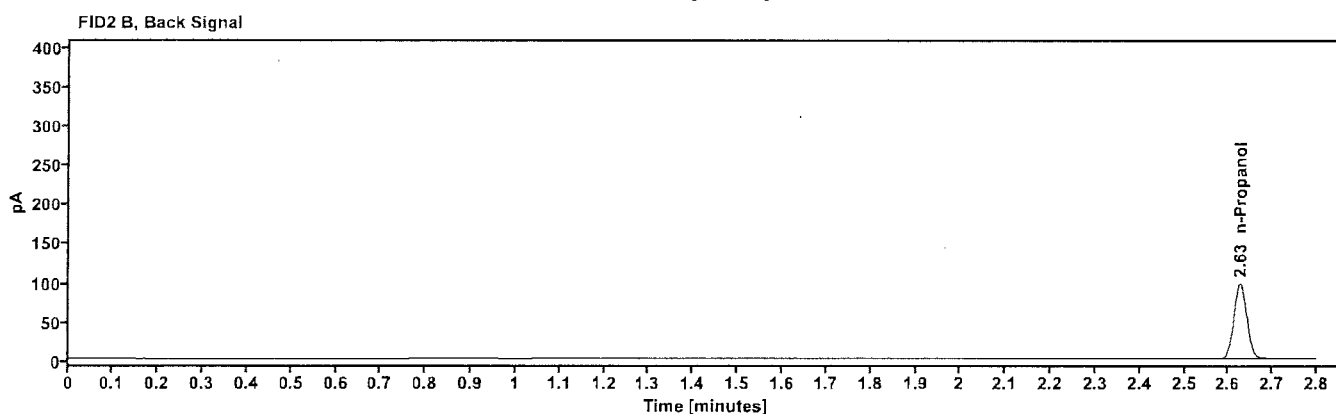
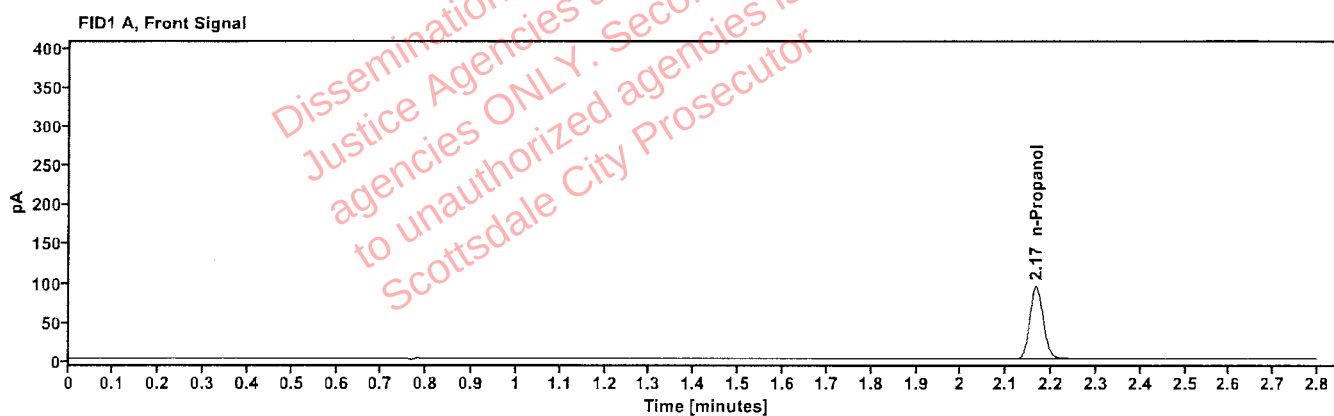


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Amount (g/100mL) | Time (min) | Peak Area |
|------------|------------------|------------|-----------|
| n-Propanol | -----            | 2.166      | 193.784   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| n-Propanol | 2.626      | 200.818   |

# Sequence Summary


Page 1 of 2

Sequence name: 22Sep26

Instrument: US2320A024 US23180007

Analyst: Abbott

| Vial | Sample       | Description      | Type        | LIMS ID | Method          |
|------|--------------|------------------|-------------|---------|-----------------|
| 1    | FN03062417   | 0.020 Calibrator | Calibration |         | ethanol quant.M |
| 2    | FN04042408   | 0.100 Calibrator | Calibration |         | ethanol quant.M |
| 3    | FN07012401   | 0.200 Calibrator | Calibration |         | ethanol quant.M |
| 4    | FN04092401   | 0.400 Calibrator | Calibration |         | ethanol quant.M |
| 5    | 061026VWB    | negative         | Control     |         | ethanol quant.M |
| 6    | 031726MIX    | Volatile Mix     | Control     |         | ethanol quant.M |
| 7    | C1475825     | 0.400            | Control     |         | ethanol quant.M |
| 8    | 529862       | 0.040            | Control     |         | ethanol quant.M |
| 9    | 4190425129/6 | 0.201            | Control     |         | ethanol quant.M |
| 10   | 436-01S.C.   |                  | Sample      |         | ethanol quant.M |
| 11   | 436-01S.C.   |                  | Sample      |         | ethanol quant.M |
| 12   | 347-01B.G.   |                  | Sample      |         | ethanol quant.M |
| 13   | 347-01B.G.   |                  | Sample      |         | ethanol quant.M |
| 14   | 438-01P.A.   |                  | Sample      |         | ethanol quant.M |
| 15   | 438-01P.A.   |                  | Sample      |         | ethanol quant.M |
| 16   | 434-01I.Q.   |                  | Sample      |         | ethanol quant.M |
| 17   | 434-01I.Q.   |                  | Sample      |         | ethanol quant.M |
| 18   | 327-01C.P.   |                  | Sample      |         | ethanol quant.M |
| 19   | 327-01C.P.   |                  | Sample      |         | ethanol quant.M |
| 20   | C1455120     | 0.100            | Control     |         | ethanol quant.M |
| 21   | 327-01F.N.   |                  | Sample      |         | ethanol quant.M |
| 22   | 327-01F.N.   |                  | Sample      |         | ethanol quant.M |
| 23   | 327-01C.C.   |                  | Sample      |         | ethanol quant.M |
| 24   | 327-01C.C.   |                  | Sample      |         | ethanol quant.M |
| 25   | 2618391-2    |                  | Sample      |         | ethanol quant.M |
| 26   | 2618391-2    |                  | Sample      |         | ethanol quant.M |
| 27   | 2617901-2    |                  | Sample      |         | ethanol quant.M |
| 28   | 2617901-2    |                  | Sample      |         | ethanol quant.M |
| 29   | 2618383-2    |                  | Sample      |         | ethanol quant.M |
| 30   | 2618383-2    |                  | Sample      |         | ethanol quant.M |
| 31   | 4190425129/6 | 0.201            | Control     |         | ethanol quant.M |
| 32   | 2618498-1    |                  | Sample      |         | ethanol quant.M |
| 33   | 2618498-1    |                  | Sample      |         | ethanol quant.M |
| 34   | 2618512-1    |                  | Sample      |         | ethanol quant.M |
| 35   | 2618512-1    |                  | Sample      |         | ethanol quant.M |
| 36   | 2618522-1    |                  | Sample      |         | ethanol quant.M |
| 37   | 2618522-1    |                  | Sample      |         | ethanol quant.M |
| 38   | 2618555-3    |                  | Sample      |         | ethanol quant.M |
| 39   | 2618555-3    |                  | Sample      |         | ethanol quant.M |
| 40   | 2618554-1    |                  | Sample      |         | ethanol quant.M |
| 41   | 2618554-1    |                  | Sample      |         | ethanol quant.M |
| 42   | C1455120     | 0.100            | Control     |         | ethanol quant.M |
| 43   | 2618544-2    |                  | Sample      |         | ethanol quant.M |
| 44   | 2618544-2    |                  | Sample      |         | ethanol quant.M |
| 45   | 2618053-2    |                  | Sample      |         | ethanol quant.M |
| 46   | 2618053-2    |                  | Sample      |         | ethanol quant.M |
| 47   | 2618598-1    |                  | Sample      |         | ethanol quant.M |

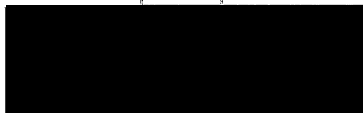

## Sequence Summary

Page 2 of 2

|    |              |          |         |  |                 |
|----|--------------|----------|---------|--|-----------------|
| 48 | 2618598-1    |          | Sample  |  | ethanol quant.M |
| 49 | 2618621-2    |          | Sample  |  | ethanol quant.M |
| 50 | 2618621-2    |          | Sample  |  | ethanol quant.M |
| 51 | 2618610-1    |          | Sample  |  | ethanol quant.M |
| 52 | 2618610-1    |          | Sample  |  | ethanol quant.M |
| 53 | 4190425129/6 | 0.201    | Control |  | ethanol quant.M |
| 54 | 2618625-2    |          | Sample  |  | ethanol quant.M |
| 55 | 2618625-2    |          | Sample  |  | ethanol quant.M |
| 56 | 2618602-2    |          | Sample  |  | ethanol quant.M |
| 57 | 2618602-2    |          | Sample  |  | ethanol quant.M |
| 58 | 2609690-2    |          | Sample  |  | ethanol quant.M |
| 59 | 2609690-2    |          | Sample  |  | ethanol quant.M |
| 60 | 2605726-2    |          | Sample  |  | ethanol quant.M |
| 61 | 2605726-2    |          | Sample  |  | ethanol quant.M |
| 62 | 2618613-2    |          | Sample  |  | ethanol quant.M |
| 63 | 2618613-2    |          | Sample  |  | ethanol quant.M |
| 64 | C1475825     | 0.400    | Control |  | ethanol quant.M |
| 65 | 529862       | 0.040    | Control |  | ethanol quant.M |
| 66 | 4190425129/6 | 0.201    | Control |  | ethanol quant.M |
| 67 | 073126AQ     | negative | Control |  | ethanol quant.M |

# Scottsdale Police Department Crime Laboratory

## Summary of Cases

SEQUENCE NAME: 22Sep26

ANALYST: Abbott

*du*

| Vials | Test 1 (g/dL) | Test 2 (g/dL) | Mean (g/dL) | Percent Difference* | Absolute Difference(g/dL)* |
|-------|---------------|---------------|-------------|---------------------|----------------------------|
| 10 11 | 0.0481        | 0.0478        | 0.04795     | 0.31                | 0.00015                    |
| 12 13 | 0.1688        | 0.1678        | 0.16830     | 0.30                | 0.00050                    |
| 14 15 | 0.1568        | 0.1583        | 0.15755     | 0.48                | 0.00075                    |
| 16 17 | 0.0000        | 0.0000        | 0.00000     | 0.00                | 0.00000                    |
| 18 19 | 0.0367        | 0.0371        | 0.03690     | 0.54                | 0.00020                    |
| 21 22 | 0.1145        | 0.1145        | 0.11450     | 0.00                | 0.00000                    |
| 23 24 | 0.1103        | 0.1112        | 0.11075     | 0.41                | 0.00045                    |
| 25 26 | 0.0000        | 0.1678        | 0.08390     | 100.00              | 0.08390                    |
| 27 28 | 0.2121        | 0.2130        | 0.21255     | 0.21                | 0.00045                    |
| 29 30 | 0.1109        | 0.1109        | 0.11090     | 0.00                | 0.00000                    |
| 32 33 | 0.1248        | 0.1251        | 0.12495     | 0.12                | 0.00015                    |
| 34 35 | 0.0000        | 0.2481        | 0.12405     | 100.00              | 0.12405                    |
| 36 37 | 0.2075        | 0.0000        | 0.10375     | 100.00              | 0.10375                    |
| 38 39 | 0.0754        | 0.0753        | 0.07535     | 0.07                | 0.00005                    |
| 40 41 | 0.1659        | 0.1654        | 0.16565     | 0.15                | 0.00025                    |
| 43 44 | 0.2564        | 0.2578        | 0.25710     | 0.27                | 0.00070                    |
| 45 46 | 0.1298        | 0.1256        | 0.12770     | 1.64                | 0.00210                    |
| 47 48 | 0.0294        | 0.0287        | 0.02905     | 1.20                | 0.00035                    |
| 49 50 | 0.1999        | 0.2015        | 0.20070     | 0.40                | 0.00080                    |
| 51 52 | 0.0821        | 0.0842        | 0.08315     | 1.26                | 0.00105                    |
| 54 55 | 0.1119        | 0.1116        | 0.11175     | 0.13                | 0.00015                    |
| 56 57 | 0.0996        | 0.0997        | 0.09965     | 0.05                | 0.00005                    |
| 58 59 | 0.2621        | 0.2598        | 0.26095     | 0.44                | 0.00115                    |
| 60 61 | 0.2311        | 0.2331        | 0.23210     | 0.43                | 0.00100                    |
| 62 63 | 0.2333        | 0.2339        | 0.23360     | 0.13                | 0.00030                    |

\*Calculated differences are differences from the mean of the two results.

SEQUENCE: 22Sep26

page 1 of 1

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 436-01S.C.                    | Vial:           | 10                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 3:48:59 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\10.D | Analyst:        | Abbott               |

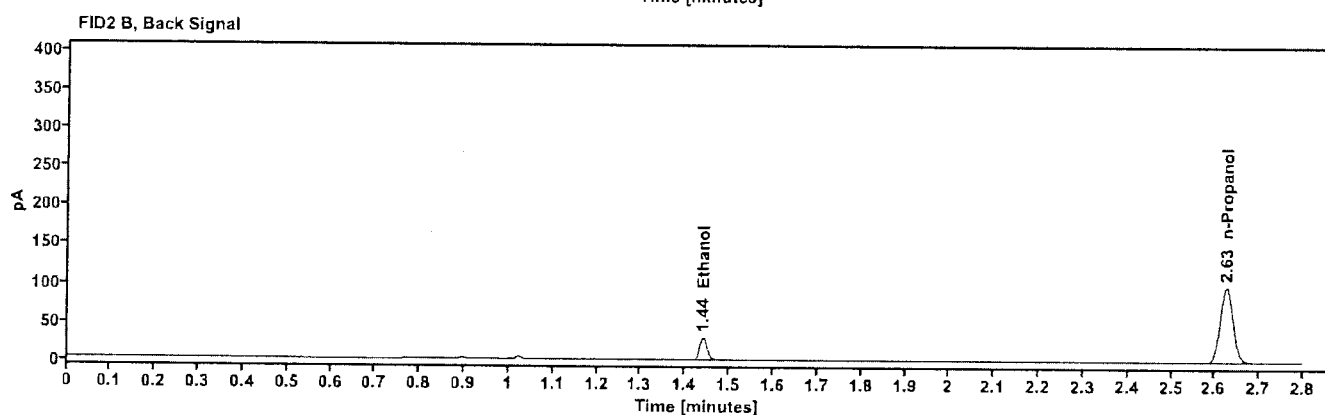
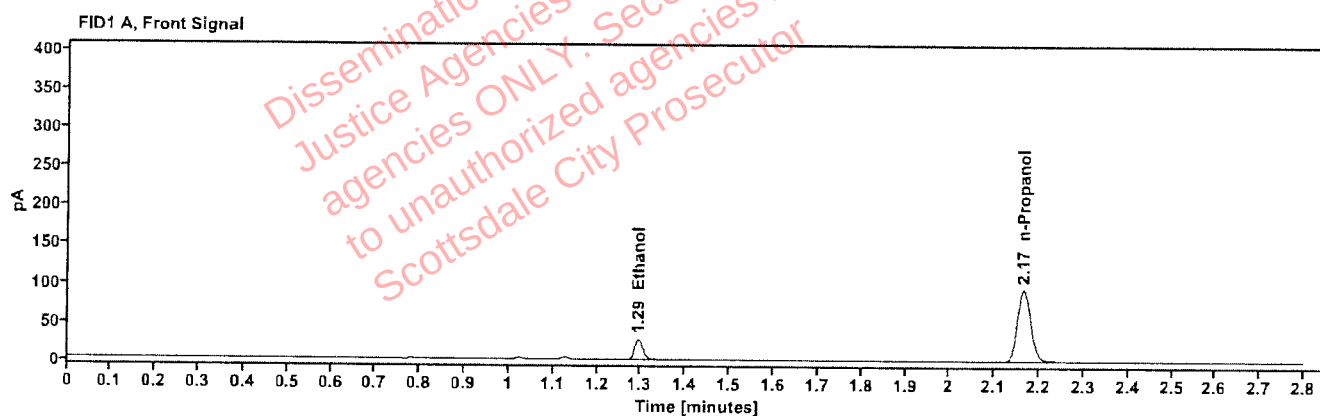


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0481           | 1.294      | 34.232    |
| n-Propanol         | ---              | 2.165      | 194.044   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 35.778    |
| n-Propanol | 2.626      | 201.269   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 436-01S.C.                    | Vial:           | 11                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 3:52:47 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\11.D | Analyst:        | Abbott               |

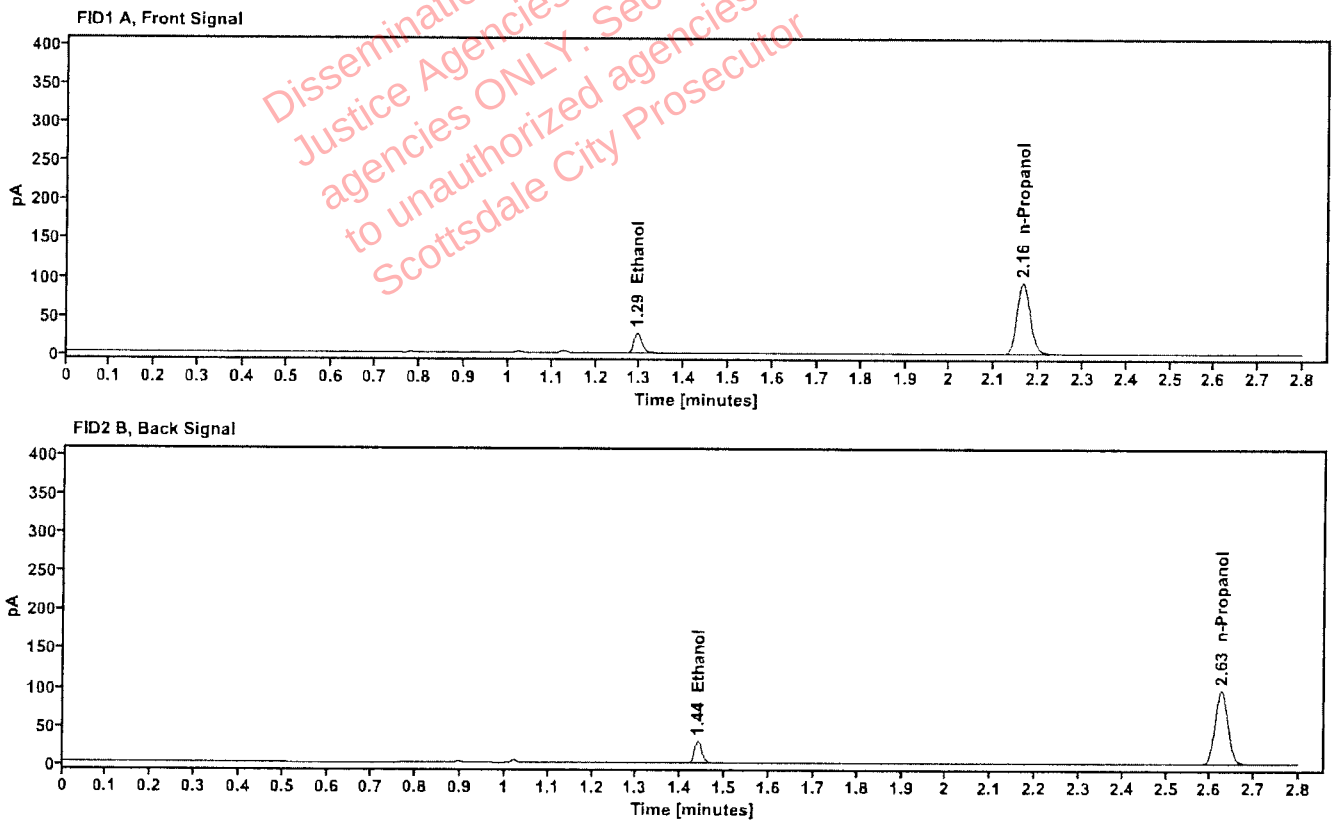


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0478           | 1.294      | 33.652    |
| n-Propanol         | ---              | 2.165      | 192.069   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 35.127    |
| n-Propanol | 2.625      | 199.231   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 347-01B.G.                    | Vial:           | 12                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 3:56:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\12.D | Analyst:        | Abbott               |

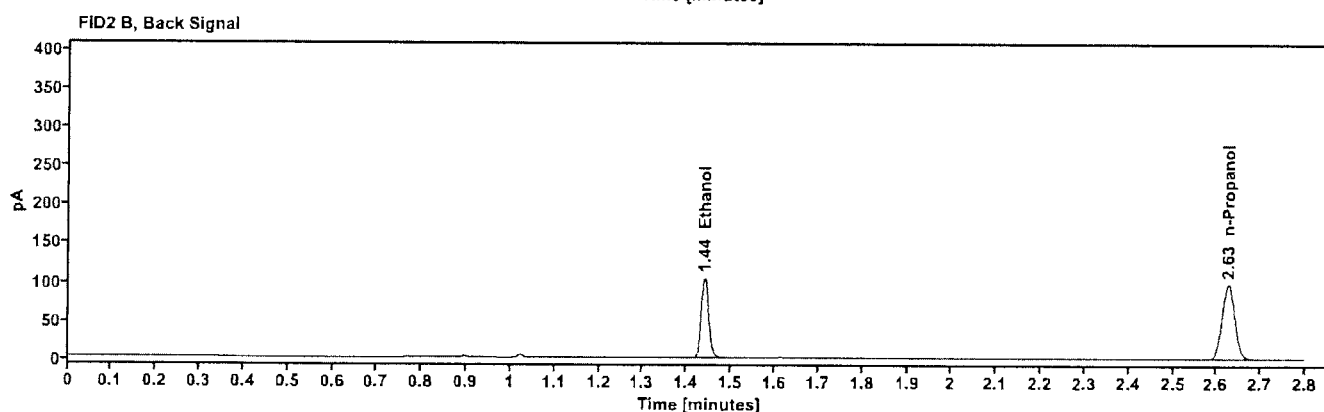
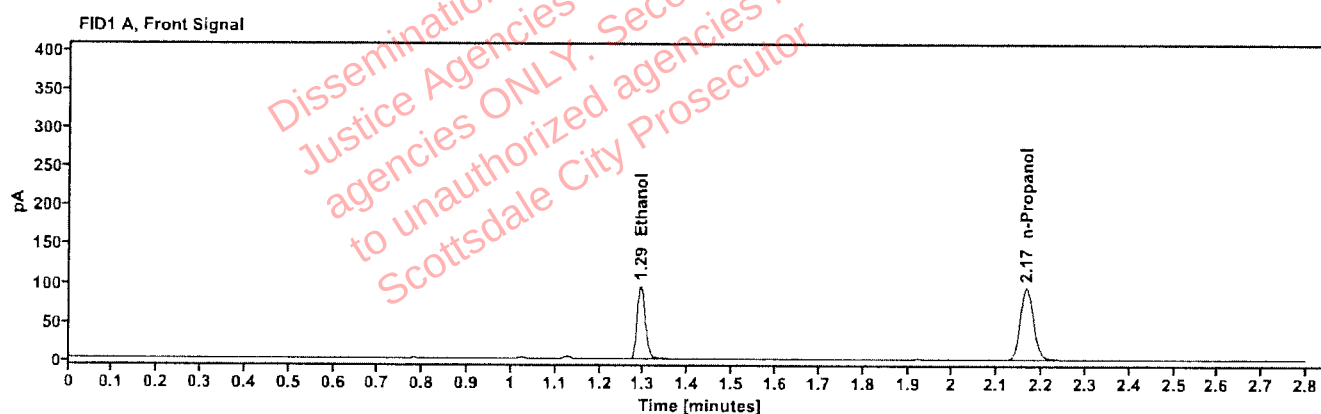


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1688           | 1.293      | 120.659   |
| n-Propanol         | -----            | 2.165      | 193.128   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 124.722   |
| n-Propanol | 2.626      | 199.947   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 347-01B.G.                    | Vial:           | 13                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:01:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\13.D | Analyst:        | Abbott               |

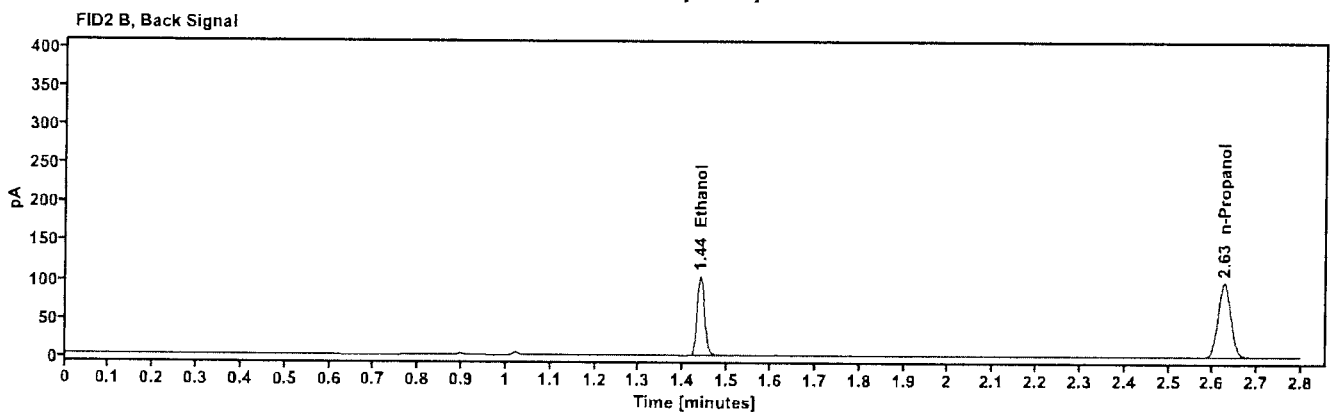
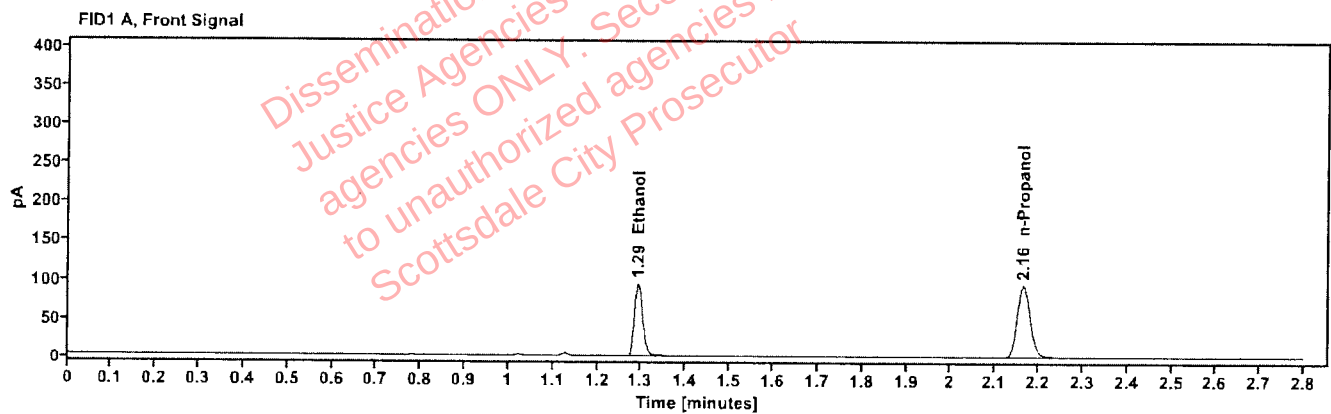


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1678           | 1.292      | 120.410   |
| n-Propanol         | ---              | 2.165      | 193.864   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 124.517   |
| n-Propanol | 2.625      | 200.838   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 438-01P.A.                     | Vial:           | 14                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:05:23 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\14.D | Analyst:        | Abbott               |

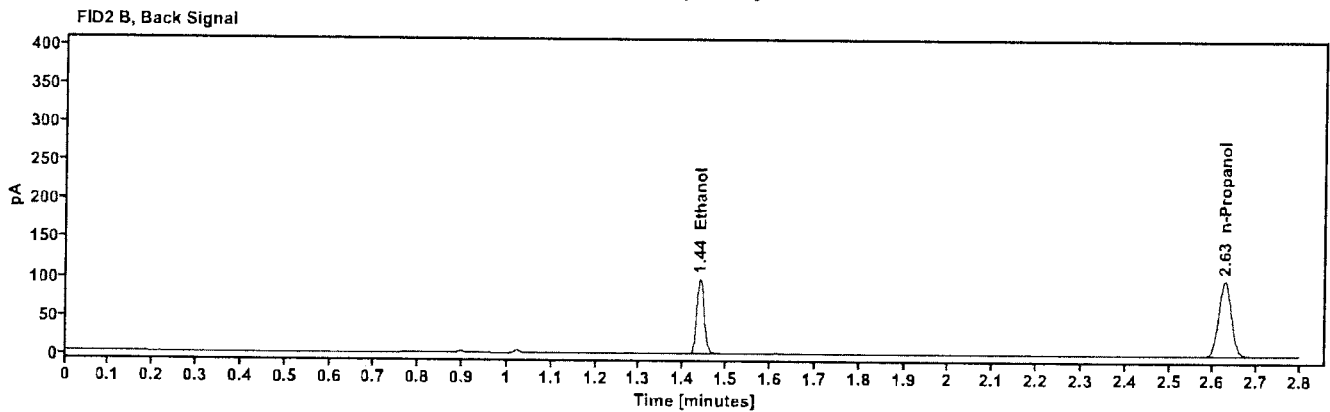
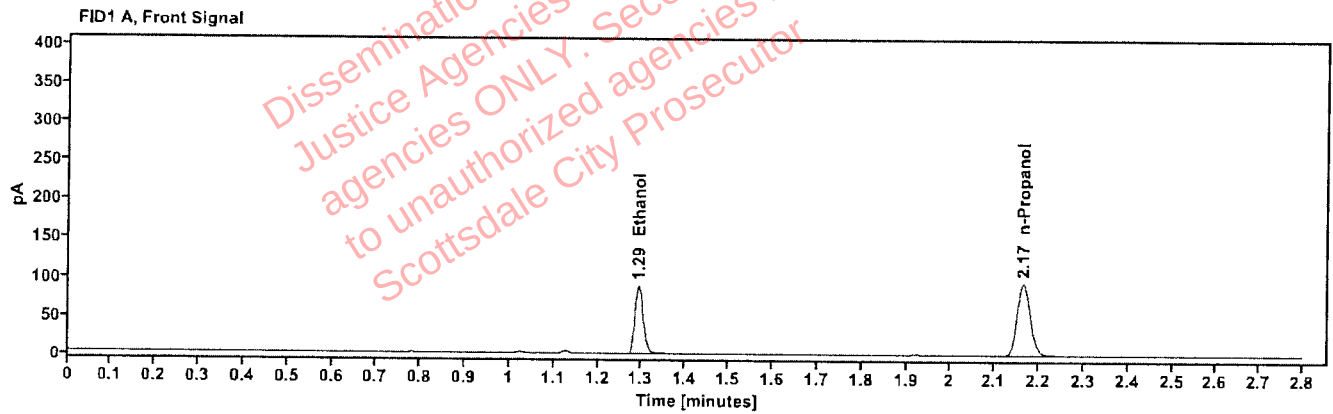


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1568           | 1.293      | 113.183   |
| n-Propanol         | -----            | 2.165      | 195.021   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 116.781   |
| n-Propanol | 2.626      | 201.894   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 438-01P.A.                     | Vial:           | 15                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:09:11 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\15.D | Analyst:        | Abbott               |

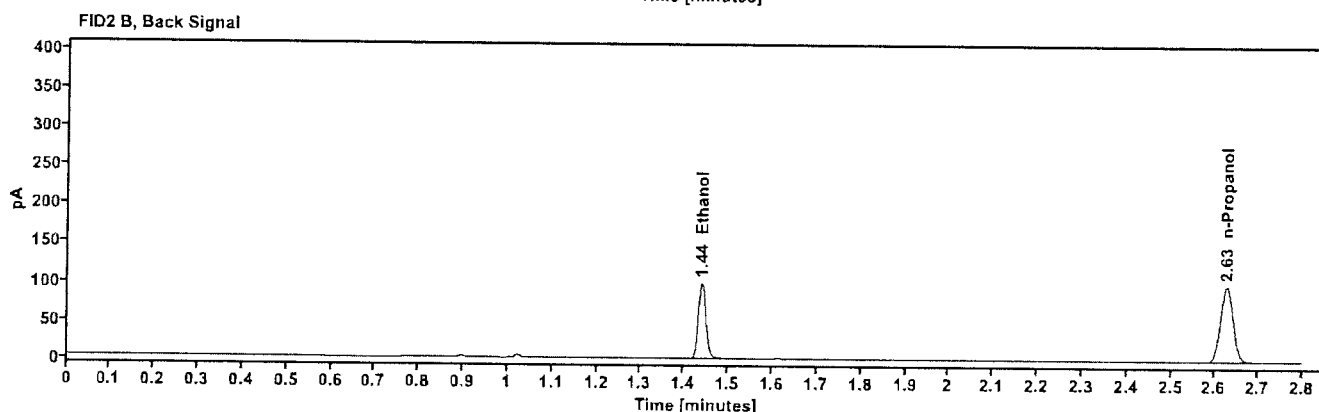
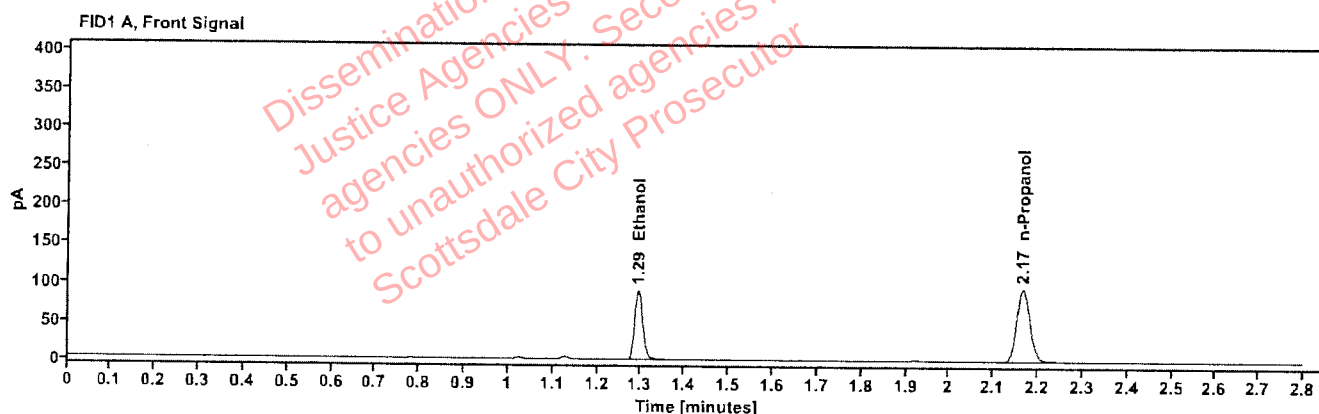


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1583           | 1.293      | 114.543   |
| n-Propanol         |                  | 2.166      | 195.495   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 118.400   |
| n-Propanol | 2.626      | 202.432   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 434-01I.Q.                    | Vial:           | 16                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:12:58 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\16.D | Analyst:        | Abbott               |

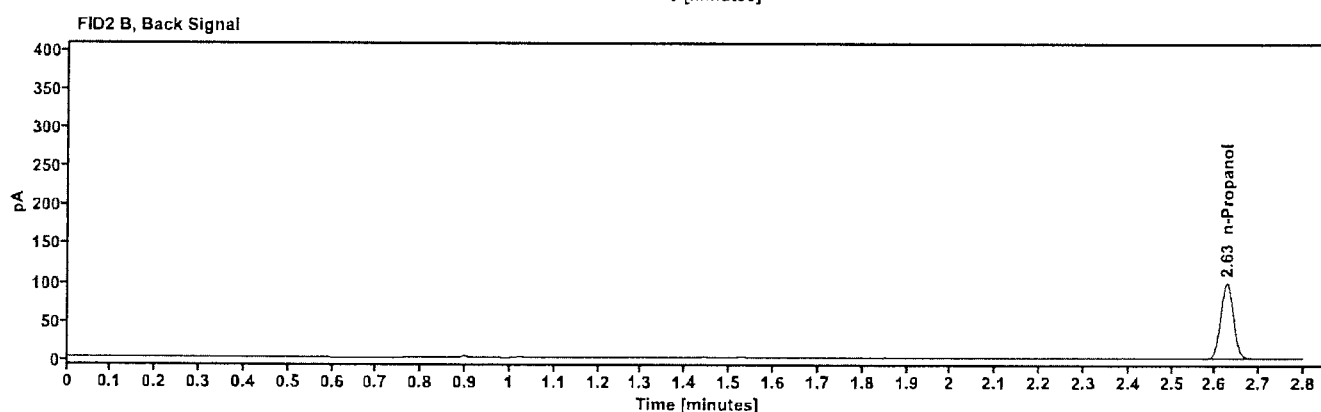
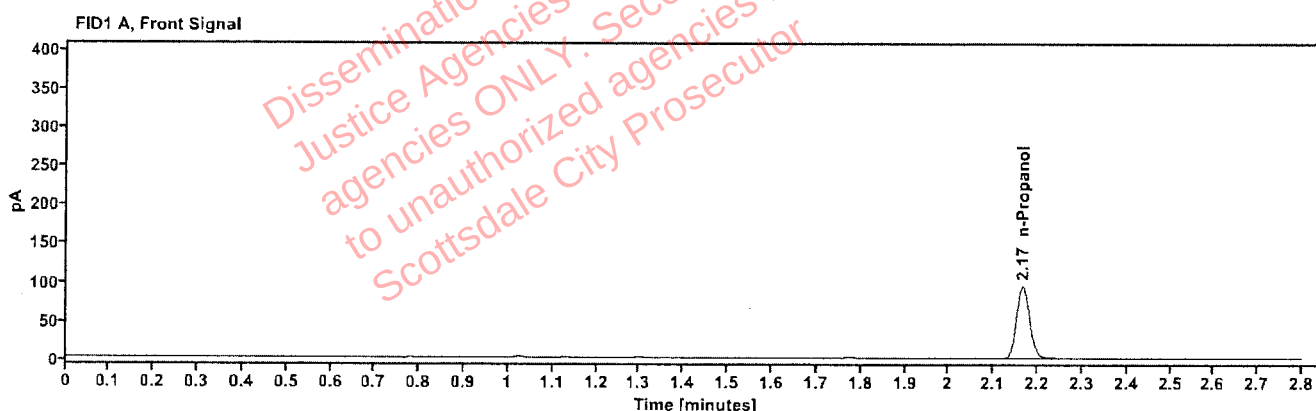


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Amount (g/100mL) | Time (min) | Peak Area |
|------------|------------------|------------|-----------|
| n-Propanol | ---              | 2.166      | 195.002   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| n-Propanol | 2.627      | 202.499   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 434-011.Q.                     | Vial:           | 17                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:16:46 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\17.D | Analyst:        | Abbott               |

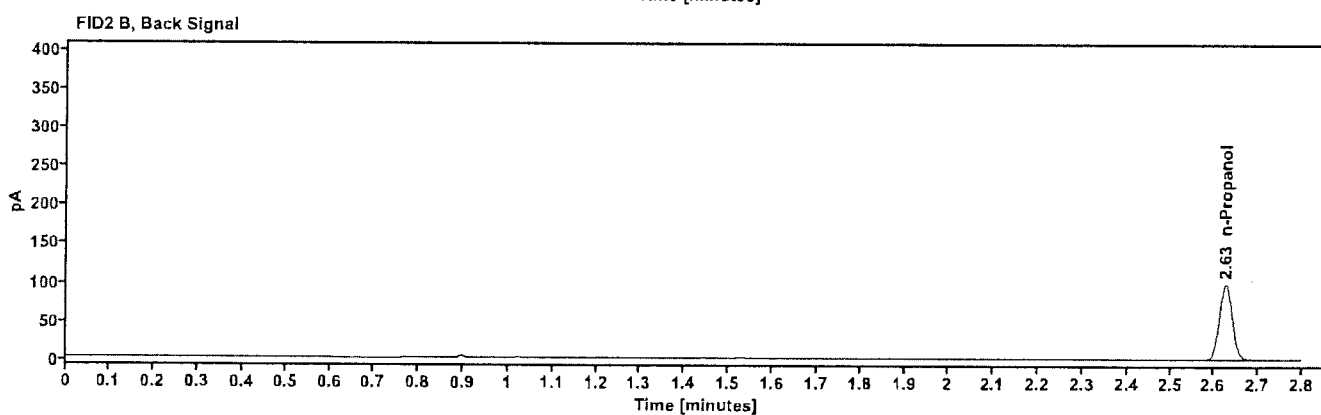
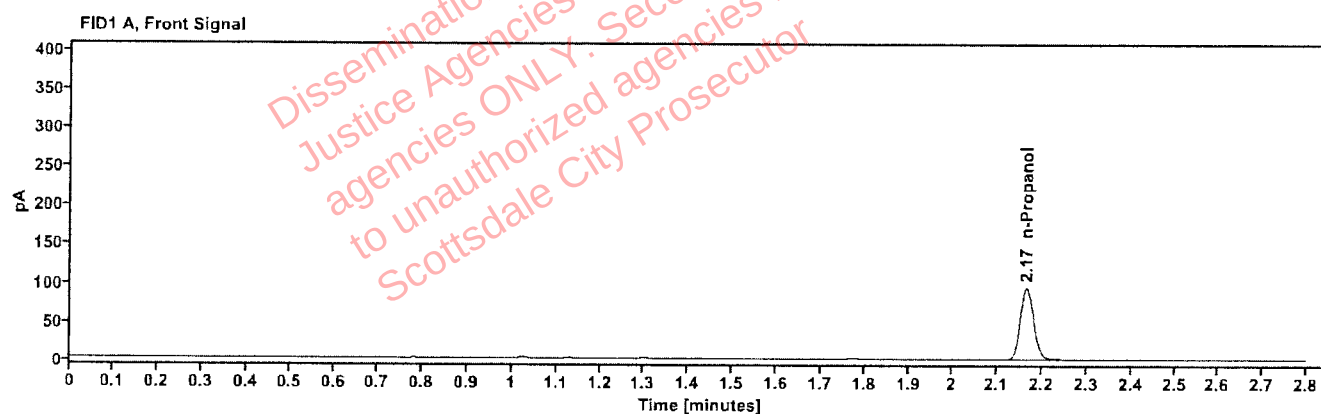


Table 1: FID 1 A (column DB-ALC1)

| Compound   | Amount (g/100mL) | Time (min) | Peak Area |
|------------|------------------|------------|-----------|
| n-Propanol | ---              | 2.165      | 194.352   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| n-Propanol | 2.626      | 201.886   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 327-01C.P.                    | Vial:           | 18                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:20:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\18.D | Analyst:        | Abbott               |

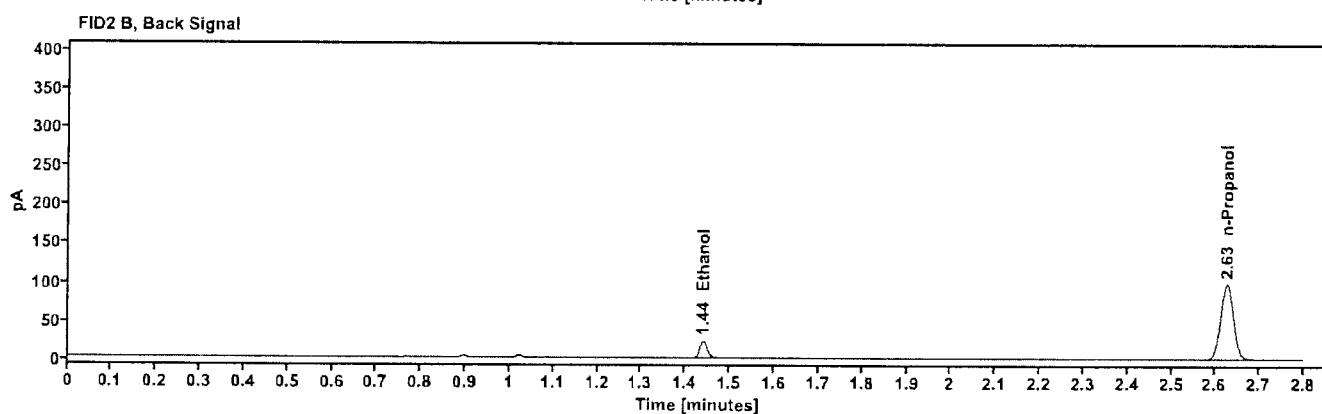
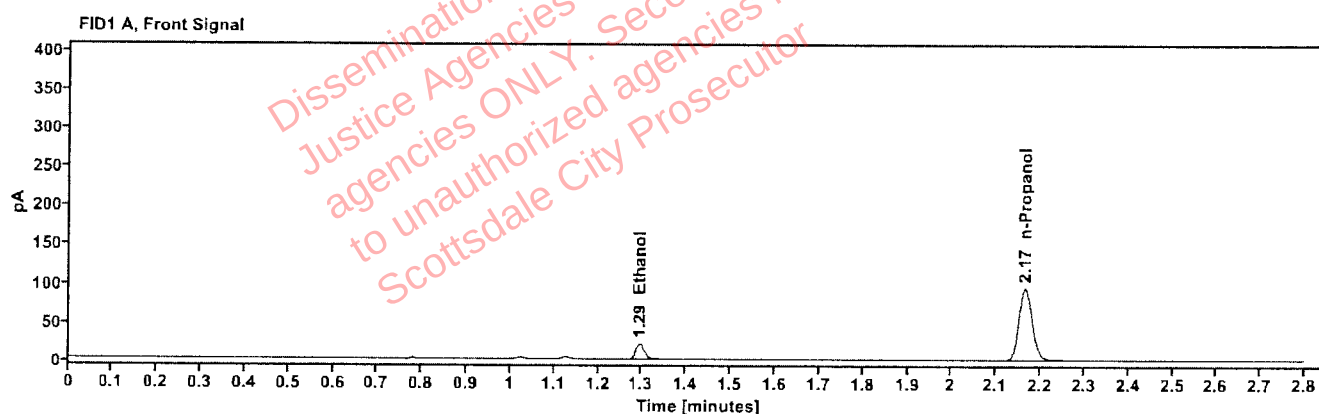


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0367           | 1.295      | 25.988    |
| n-Propanol         | —                | 2.165      | 193.626   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 27.138    |
| n-Propanol | 2.626      | 200.998   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 327-01C.P.                    | Vial:           | 19                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:25:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\19.D | Analyst:        | Abbott               |

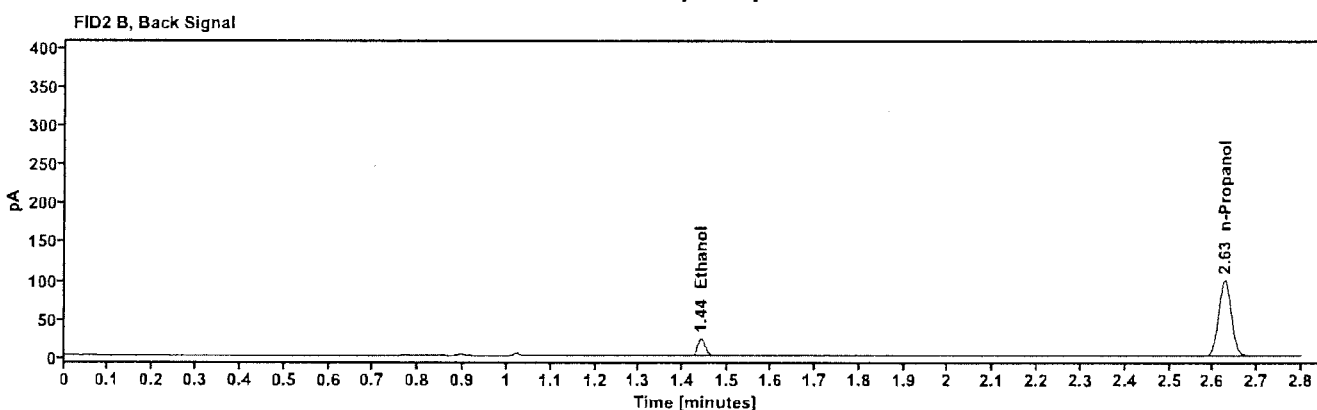
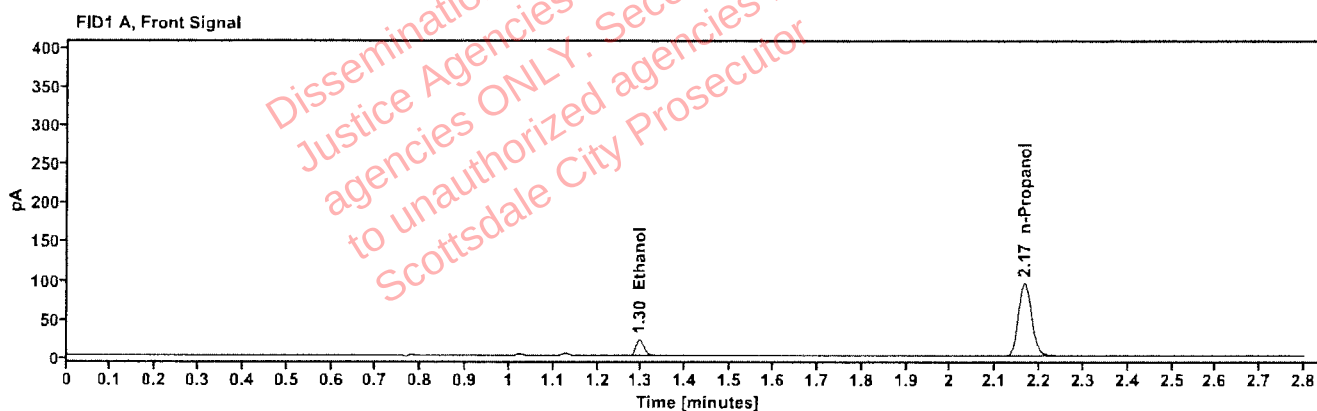


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0371           | 1.295      | 26.523    |
| n-Propanol         | —                | 2.166      | 195.754   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.441      | 27.734    |
| n-Propanol | 2.627      | 203.173   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 327-01F.N.                     | Vial:           | 21                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:33:11 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\21.D | Analyst:        | Abbott               |

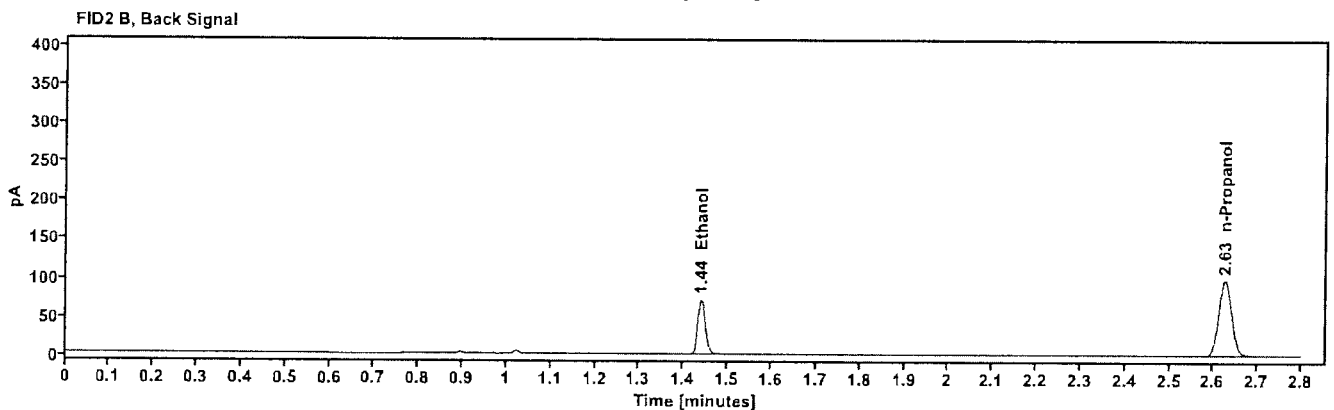
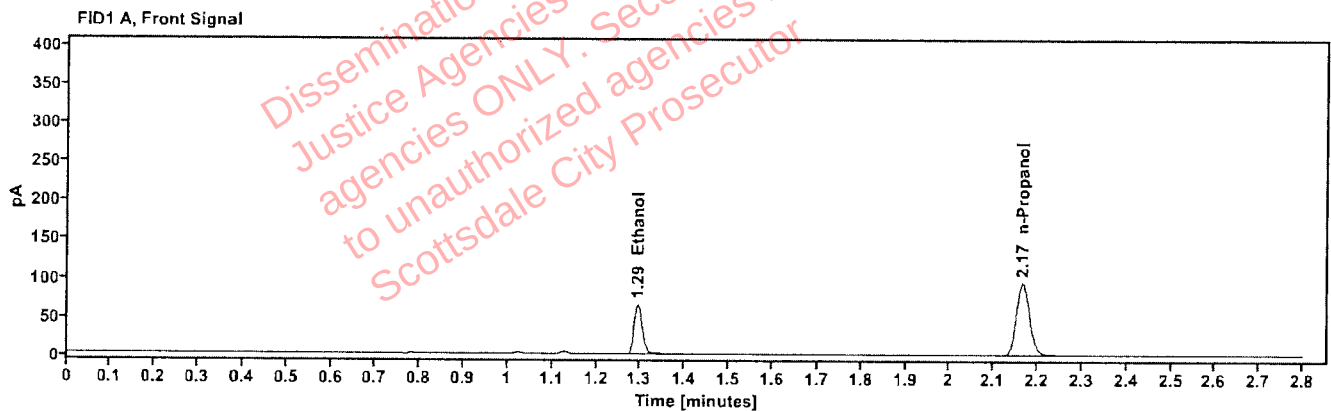


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1145           | 1.293      | 82.443    |
| n-Propanol         | -----            | 2.166      | 194.874   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 85.135    |
| n-Propanol | 2.627      | 201.745   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 327-01F.N.                     | Vial:           | 22                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:36:58 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\22.D | Analyst:        | Abbott               |

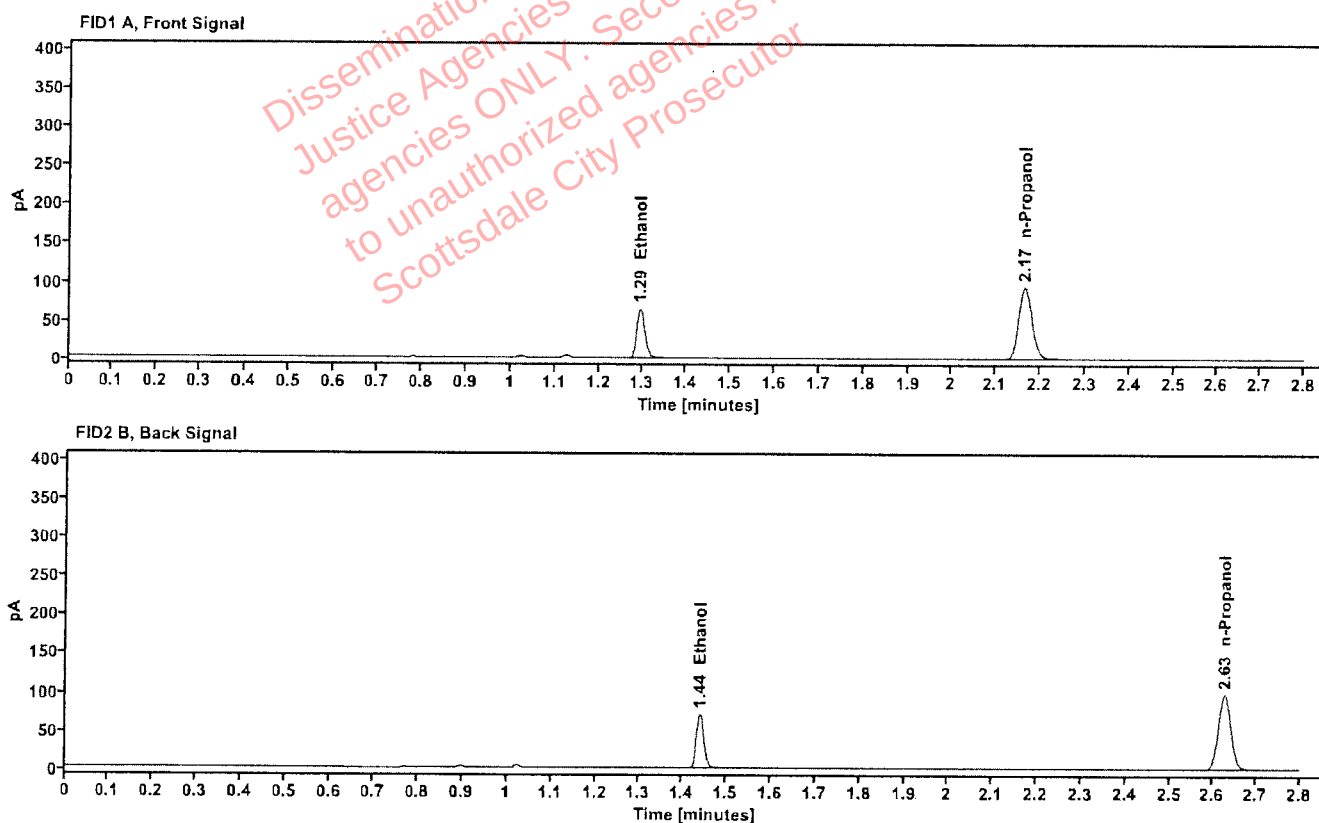


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1145           | 1.293      | 82.239    |
| n-Propanol         | —                | 2.166      | 194.306   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 84.871    |
| n-Propanol | 2.626      | 201.263   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 327-01C.C.                     | Vial:           | 23                   |
| Description: | [REDACTED]                     | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 4:40:47 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\23.D | Analyst:        | Abbott               |

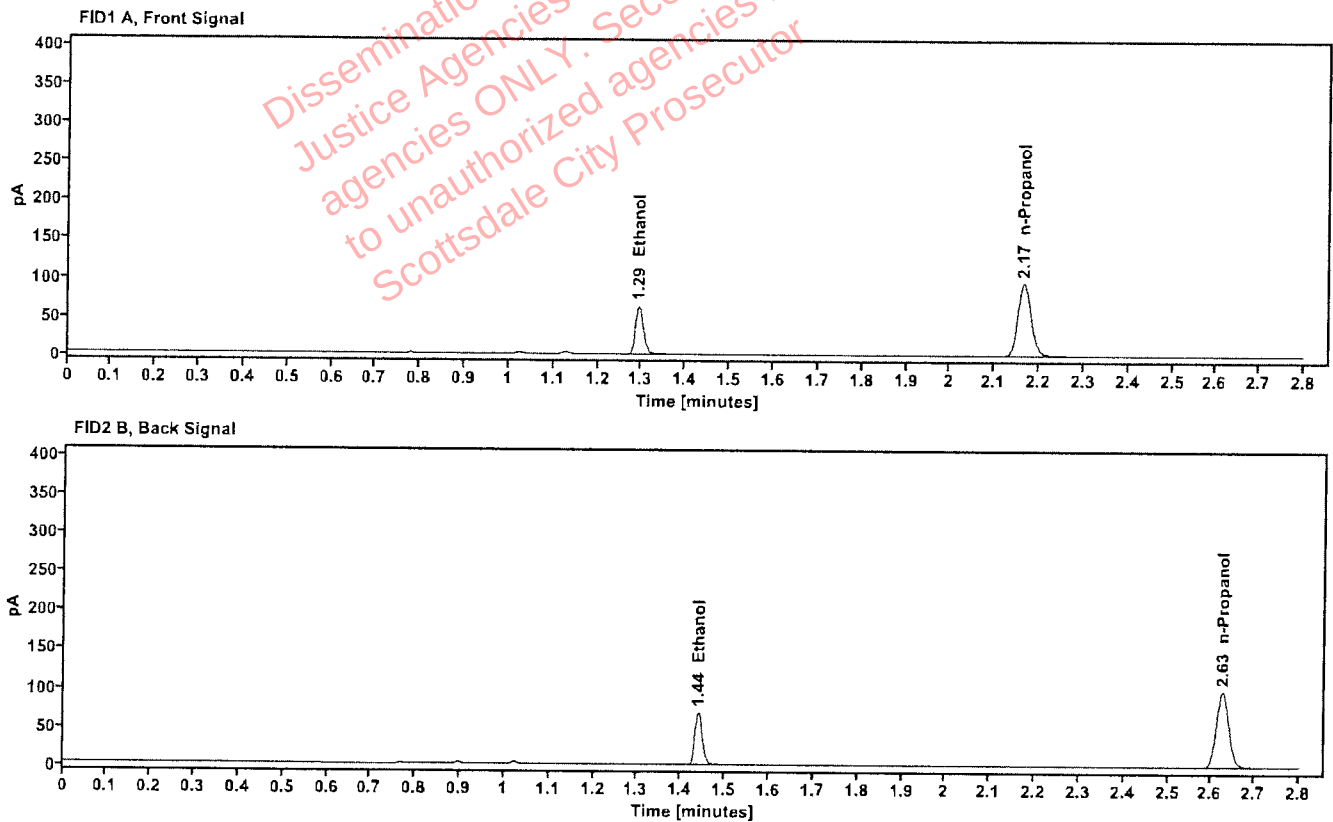


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1103           | 1.294      | 79.712    |
| n-Propanol         | ---              | 2.166      | 195.595   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 82.184    |
| n-Propanol | 2.627      | 202.573   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 327-01C.C.                    | Vial:           | 24                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:44:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\24.D | Analyst:        | Abbott               |

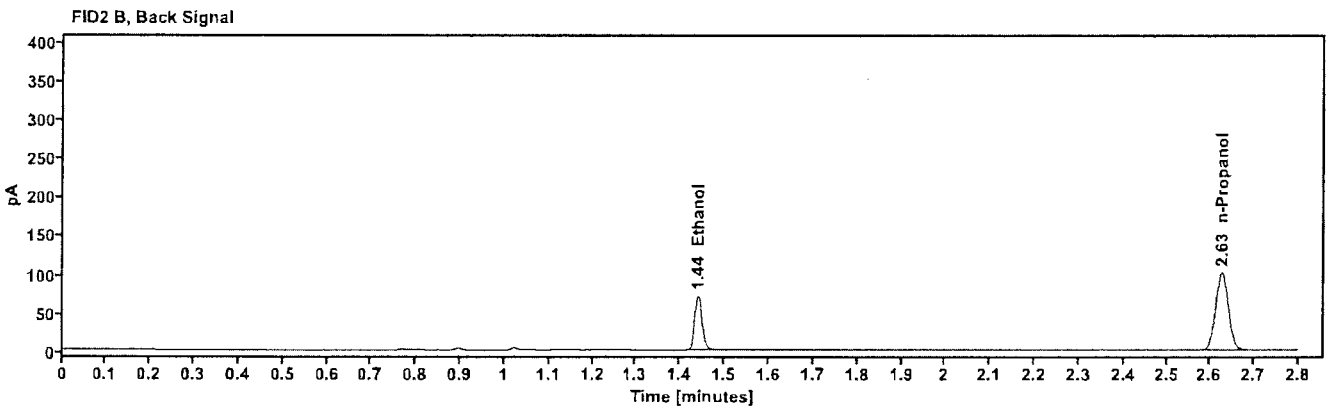
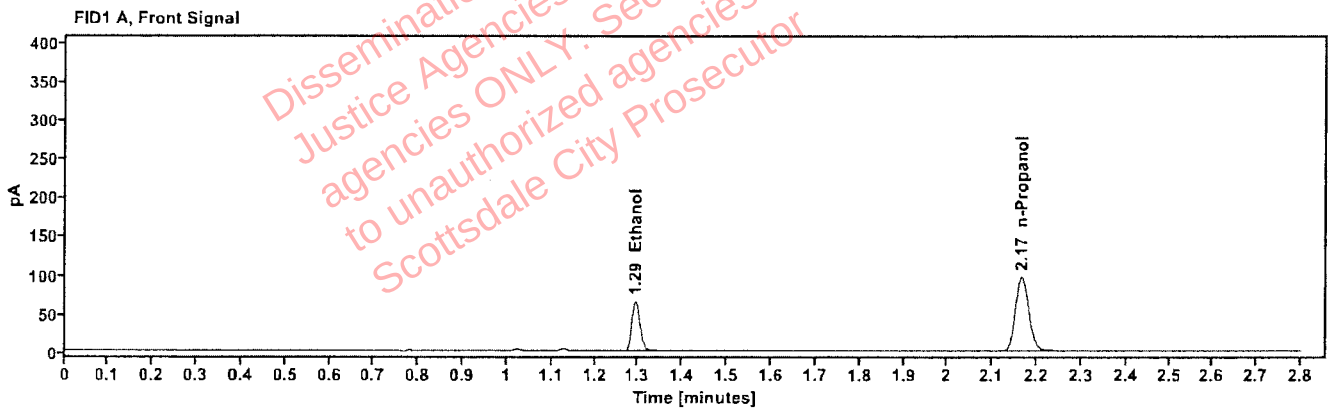


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1112           | 1.294      | 82.384    |
| n-Propanol         | —                | 2.166      | 200.400   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 84.853    |
| n-Propanol | 2.627      | 207.639   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618391-2                     | Vial:           | 25                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:49:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\25.D | Analyst:        | Abbott               |

**VOID**

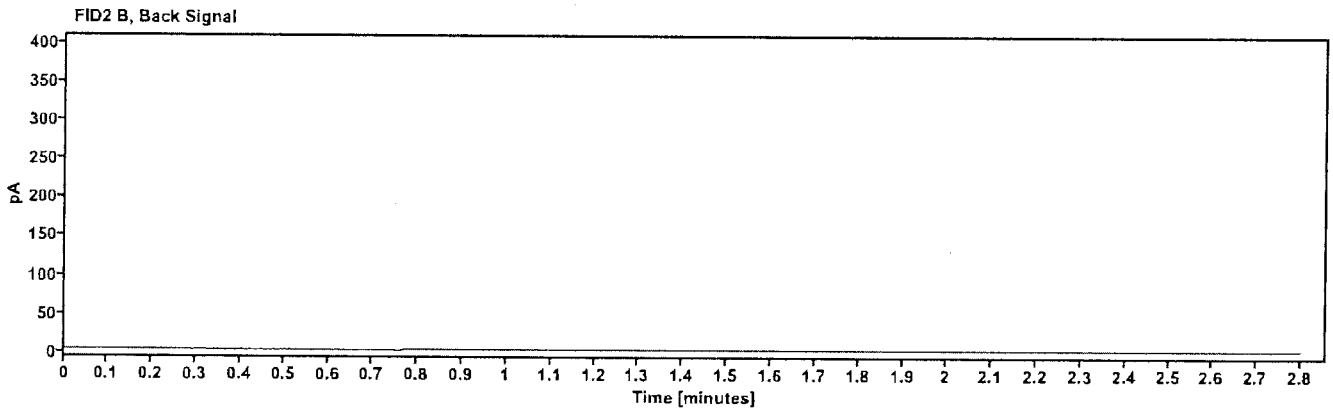
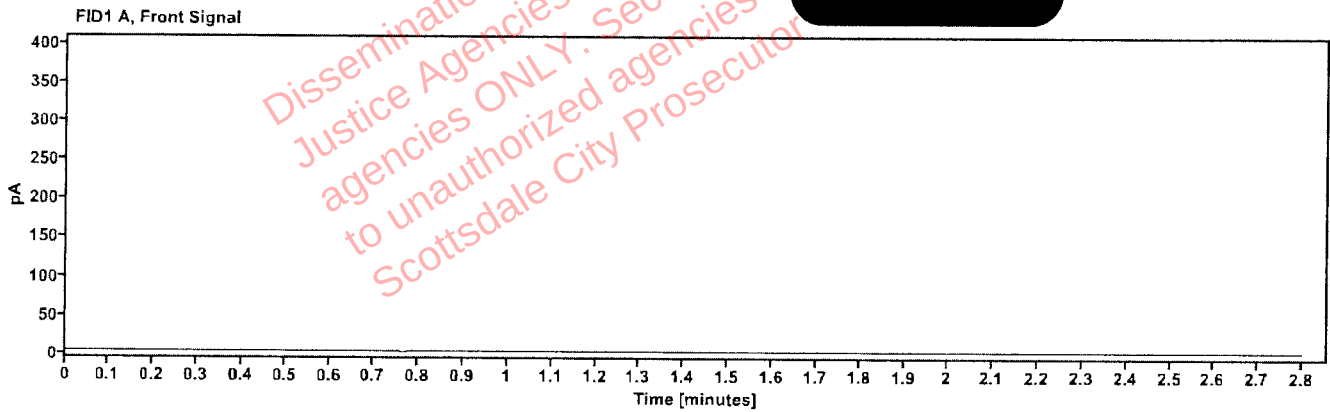


Table 1: FID 1 A (column DB-ALC1)

| Compound | Amount<br>(g/100mL) | Time<br>(min) | Peak<br>Area |
|----------|---------------------|---------------|--------------|
|----------|---------------------|---------------|--------------|

Table 2: FID 2 B (column DB-ALC2)

| Compound | Time<br>(min) | Peak<br>Area |
|----------|---------------|--------------|
|----------|---------------|--------------|

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618391-2                     | Vial:           | 26                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:53:23 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\26.D | Analyst:        | Abbott               |

**VOID**

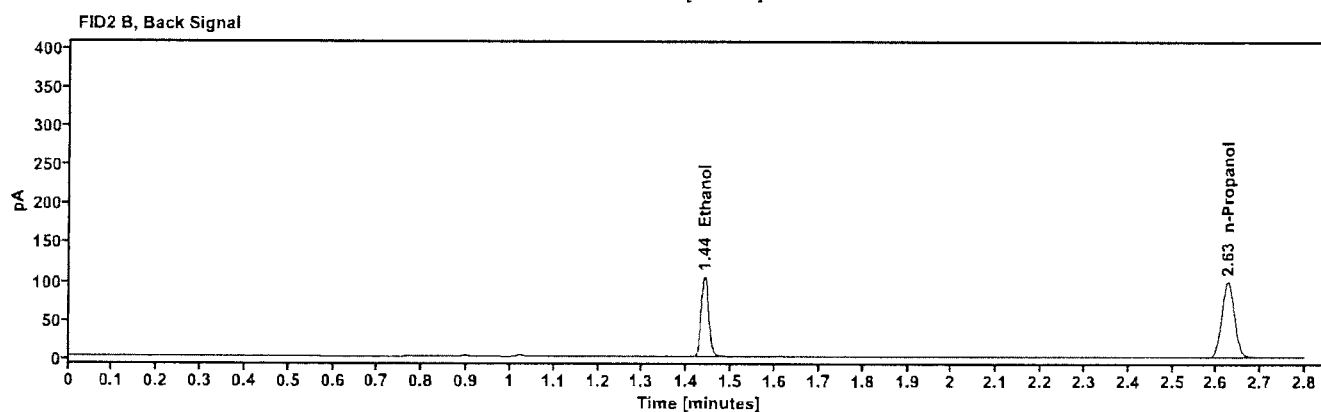
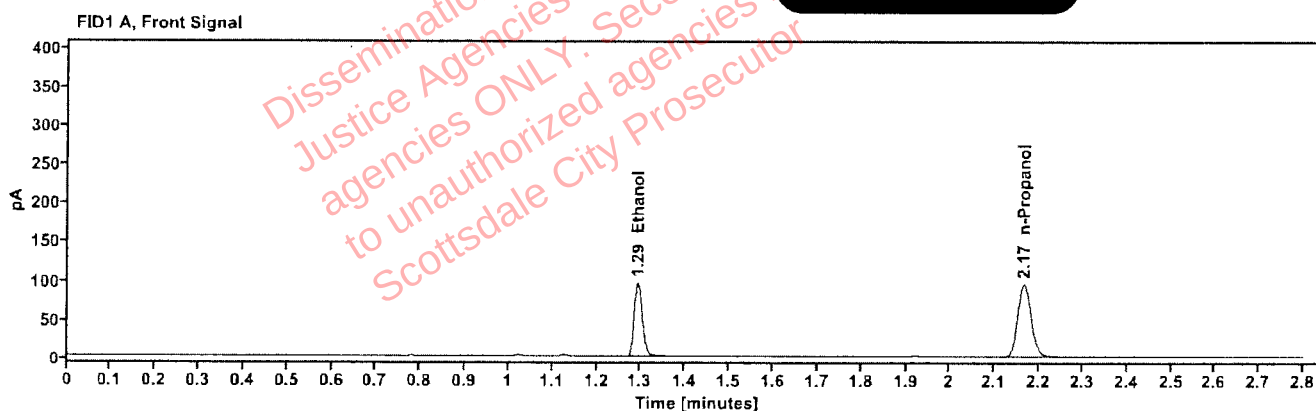


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1678           | 1.293      | 121.595   |
| n-Propanol         | —                | 2.166      | 195.739   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 125.744   |
| n-Propanol | 2.626      | 202.872   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2617901-2                     | Vial:           | 27                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 4:57:10 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\27.D | Analyst:        | Abbott               |

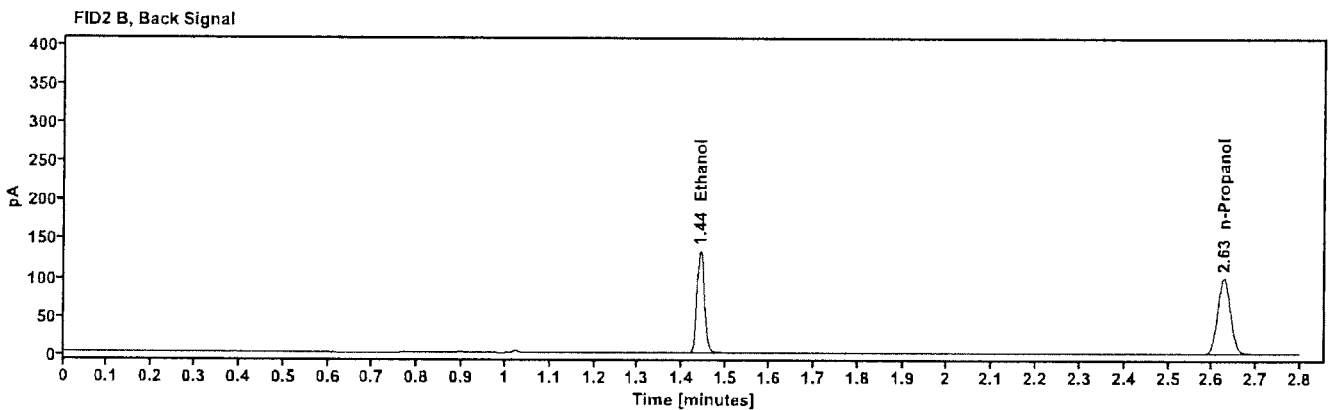
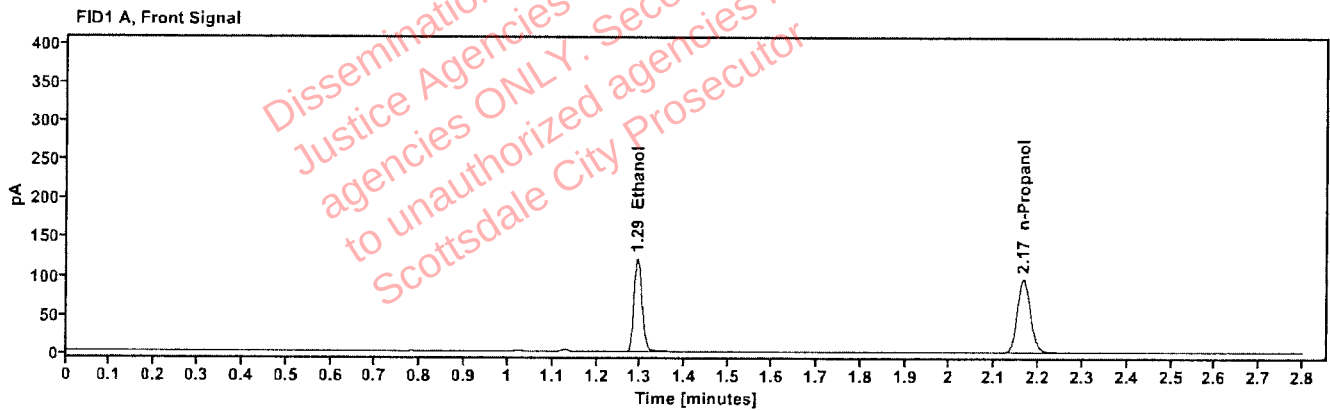


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2121           | 1.293      | 155.069   |
| n-Propanol         | ---              | 2.166      | 197.371   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 160.199   |
| n-Propanol | 2.627      | 204.227   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2617901-2                      | Vial:           | 28                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 5:00:59 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\28.D | Analyst:        | Abbott               |

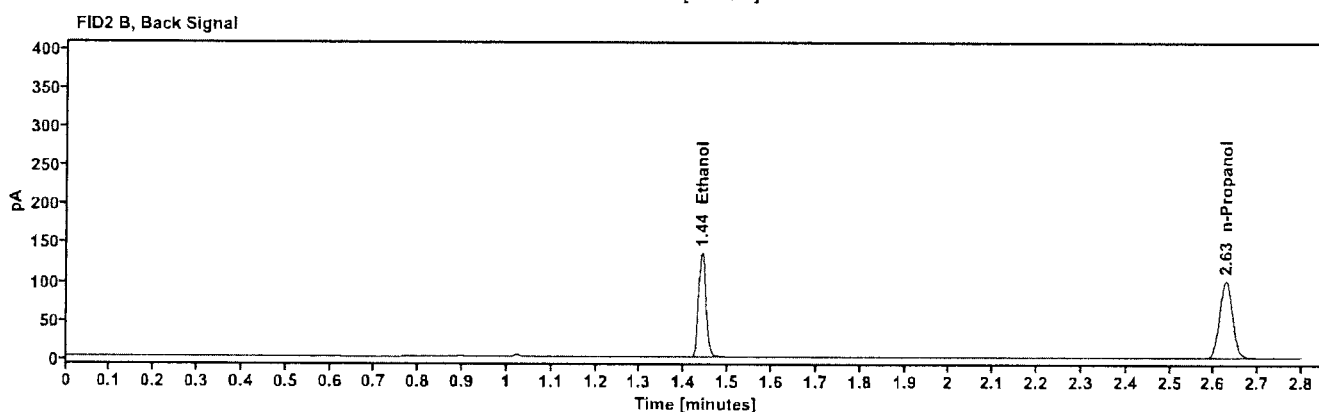
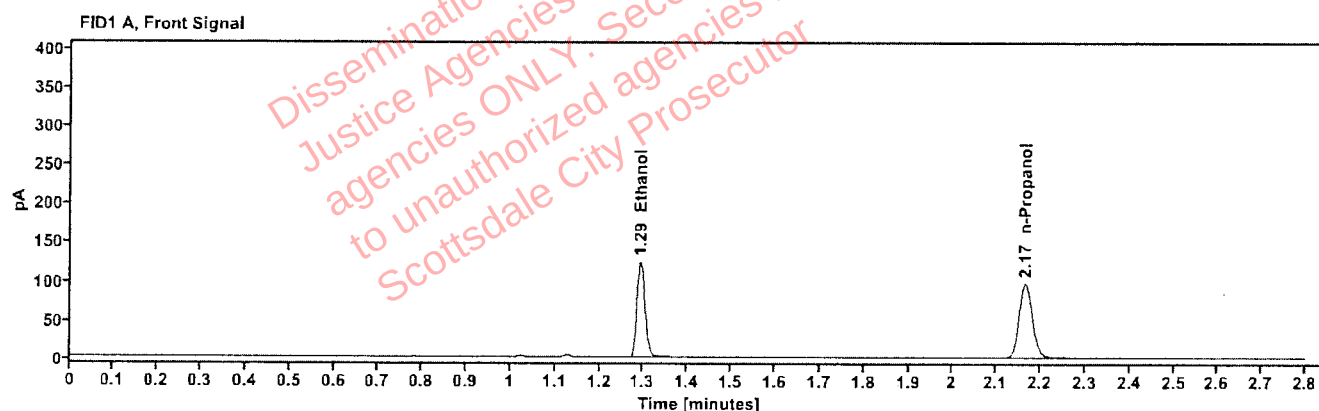


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2130           | 1.293      | 157.603   |
| n-Propanol         | ---              | 2.166      | 199.687   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 162.881   |
| n-Propanol | 2.627      | 206.744   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618383-2                     | Vial:           | 29                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:04:47 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\29.D | Analyst:        | Abbott               |

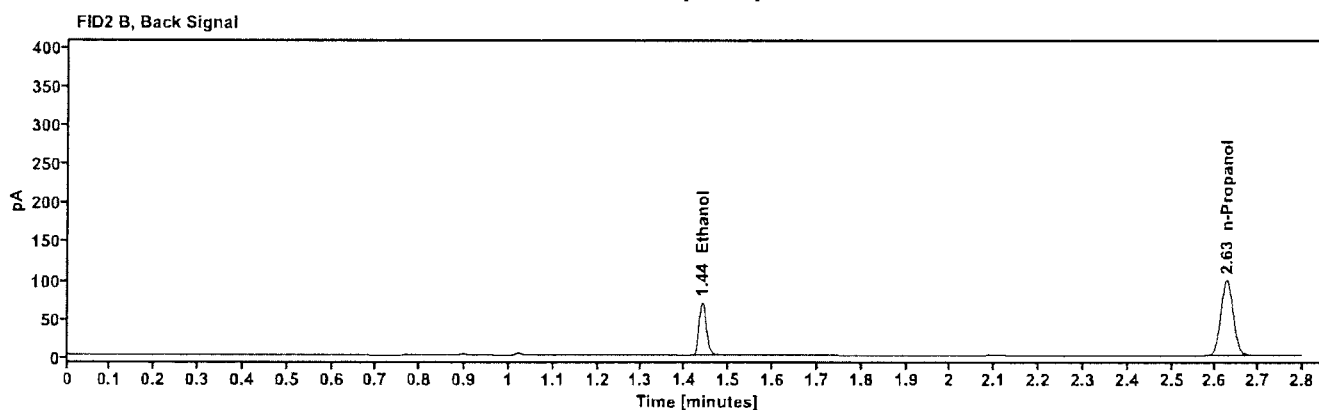
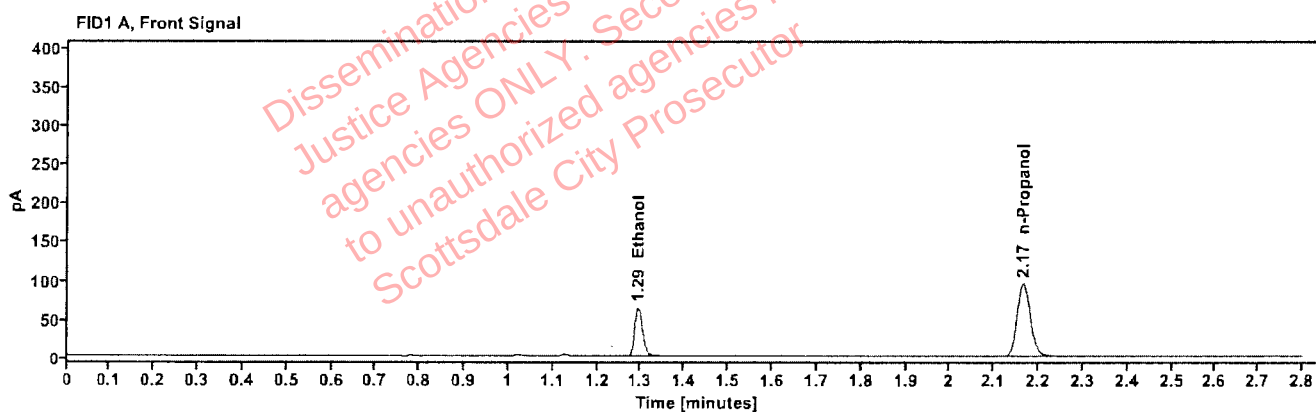


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1109           | 1.294      | 80.273    |
| n-Propanol         | ---              | 2.166      | 195.837   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 82.823    |
| n-Propanol | 2.627      | 202.762   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618383-2                      | Vial:           | 30                   |
| Description: | [REDACTED]                     | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 5:08:34 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\30.D | Analyst:        | Abbott               |

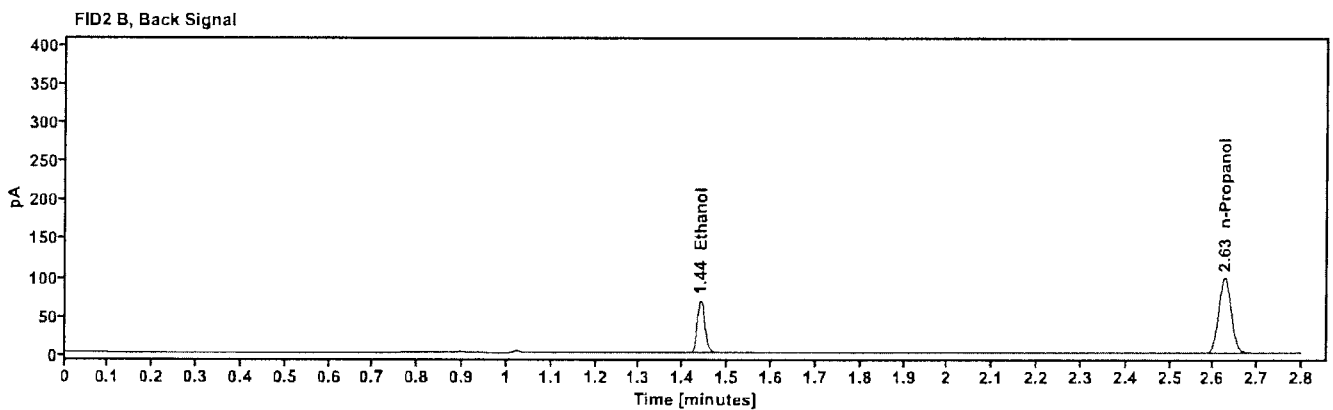
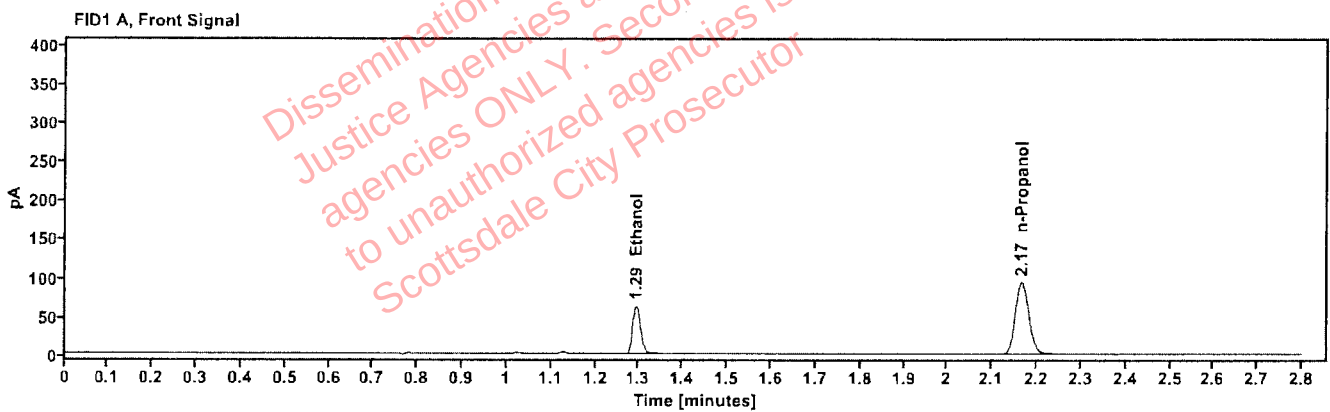


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1109           | 1.293      | 79.843    |
| n-Propanol         | ----             | 2.166      | 194.768   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 82.421    |
| n-Propanol | 2.626      | 201.865   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618498-1                      | Vial:           | 32                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 5:17:23 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\32.D | Analyst:        | Abbott               |

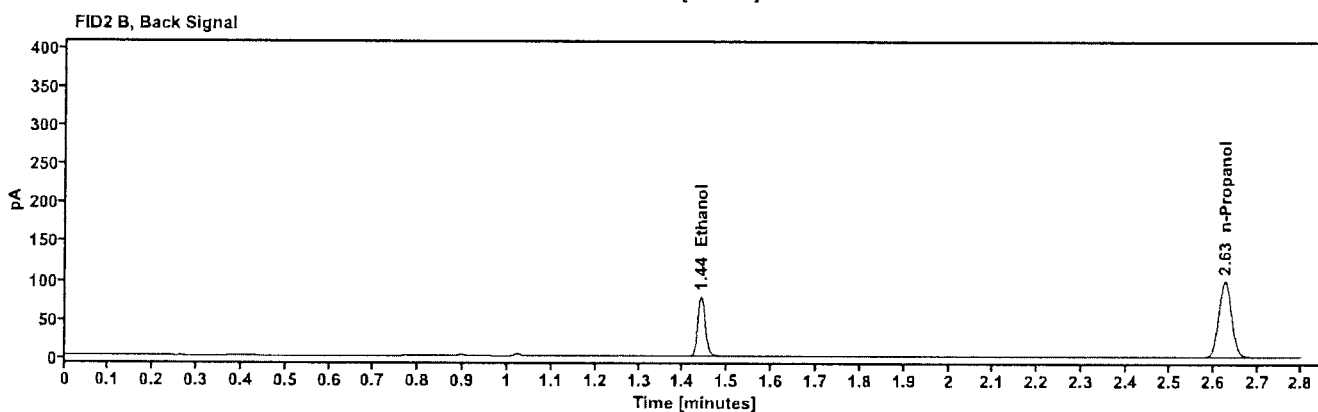
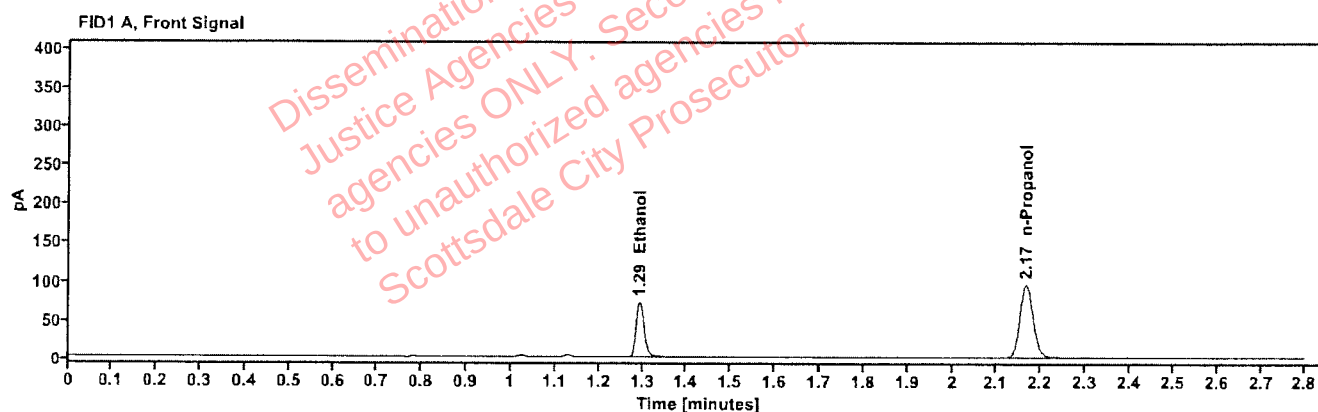


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1248           | 1.293      | 90.655    |
| n-Propanol         | ---              | 2.166      | 196.483   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 93.603    |
| n-Propanol | 2.627      | 203.852   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618498-1                     | Vial:           | 33                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:21:11 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\33.D | Analyst:        | Abbott               |

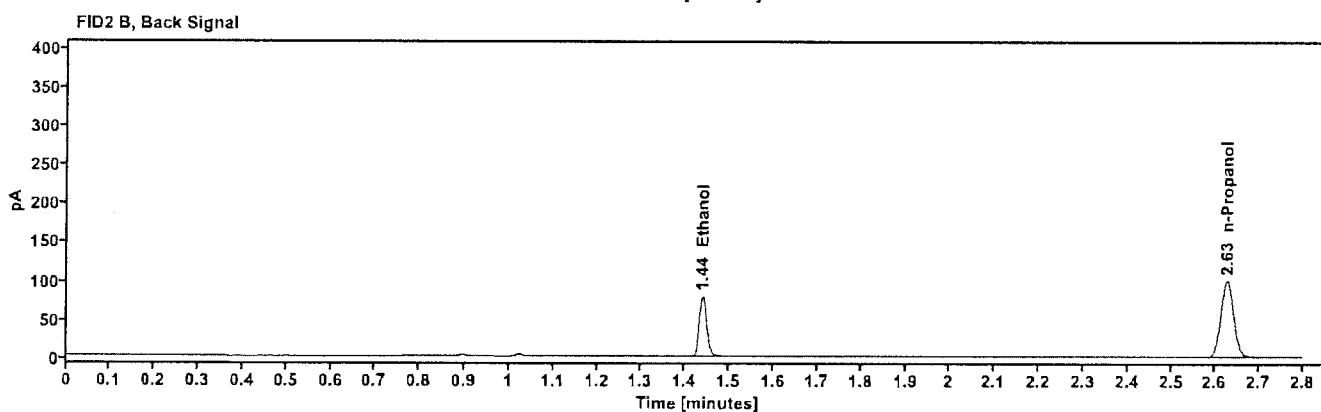
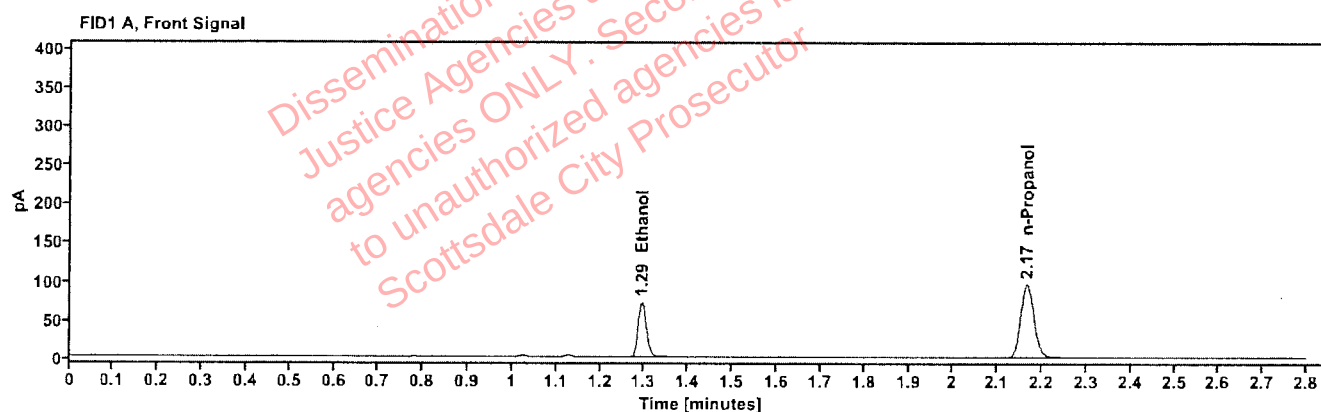


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1251           | 1.294      | 91.141    |
| n-Propanol         | -----            | 2.166      | 197.067   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 94.029    |
| n-Propanol | 2.627      | 204.464   |

Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                                                        |                 |                                                                        |
|--------------|------------------------------------------------------------------------|-----------------|------------------------------------------------------------------------|
| Sample:      | 2618512-1                                                              | Vial:           | 34                                                                     |
| Description: | <span style="background-color: black; color: black;">[REDACTED]</span> | LIMS ID:        | <span style="background-color: black; color: black;">[REDACTED]</span> |
| Method:      | ethanol quant.M                                                        | Sequence:       | 22Sep26                                                                |
| Instrument:  | US2320A024 US23180007                                                  | Injection date: | 9/22/2026 5:24:59 PM                                                   |
| Data file:   | C:\Chem32\1\Data\22Sep26\34.D                                          | Analyst:        | Abbott                                                                 |

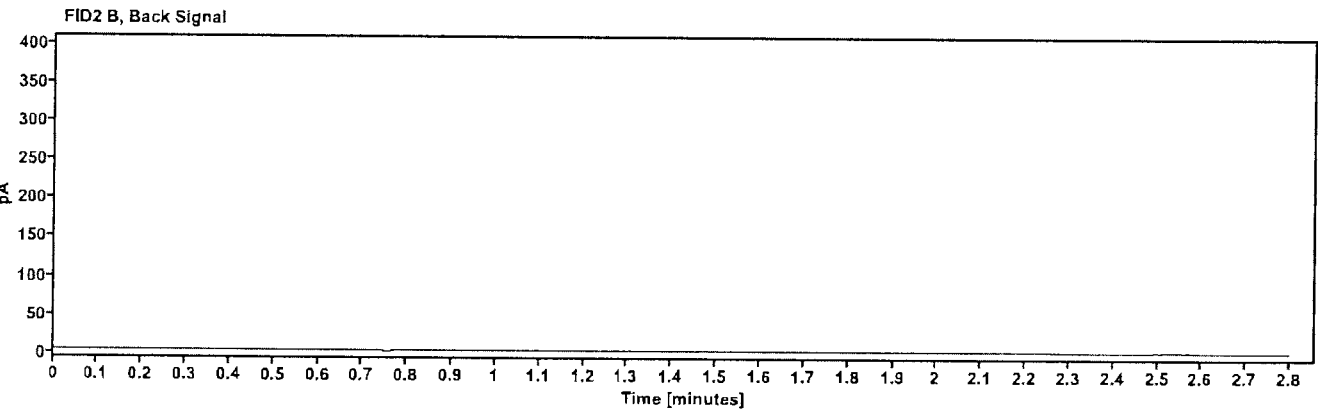
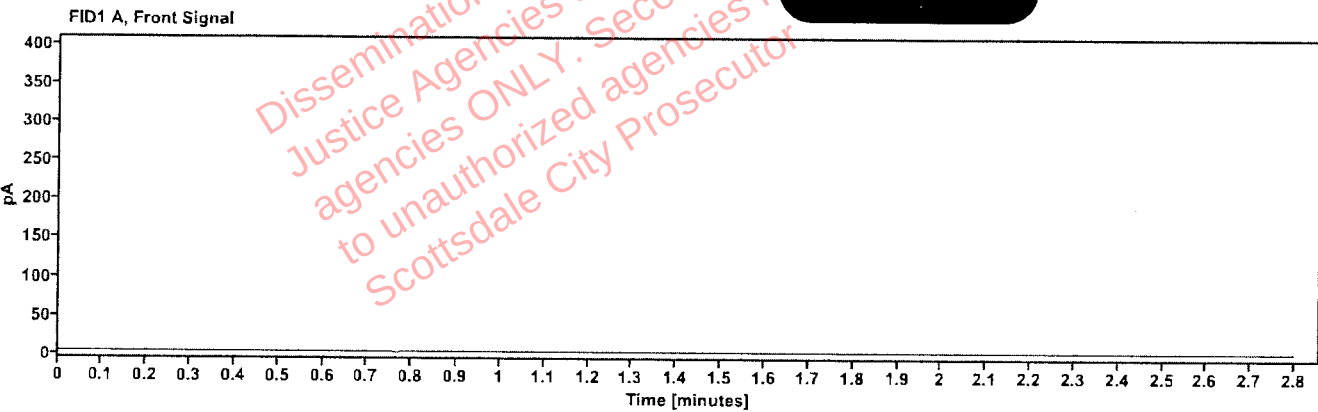


Table 1: FID 1 A (column DB-ALC1)

| Compound | Amount<br>(g/100mL) | Time<br>(min) | Peak<br>Area |
|----------|---------------------|---------------|--------------|
|----------|---------------------|---------------|--------------|

Table 2: FID 2 B (column DB-ALC2)

| Compound | Time<br>(min) | Peak<br>Area |
|----------|---------------|--------------|
|----------|---------------|--------------|

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618512-1                     | Vial:           | 35                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:28:47 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\35.D | Analyst:        | Abbott               |

**VOID**

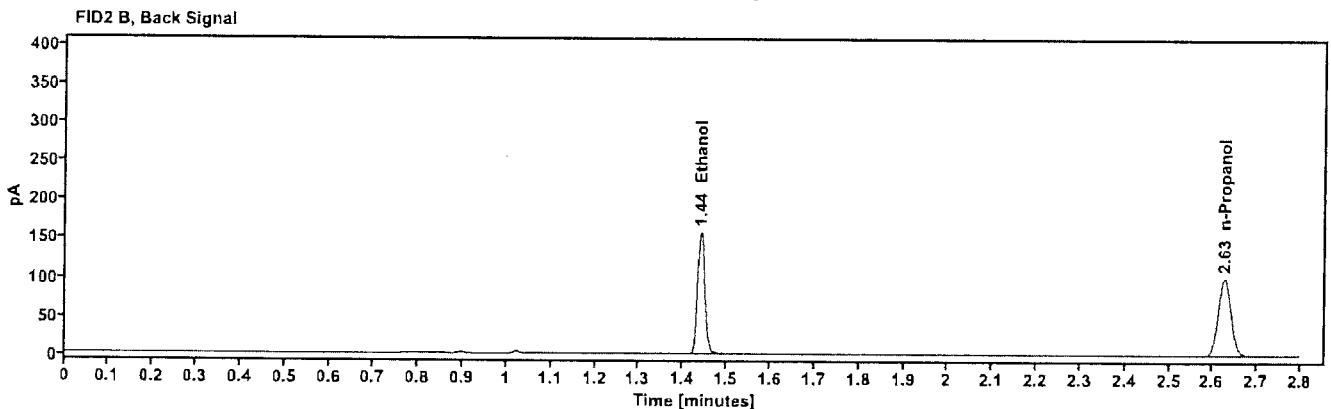
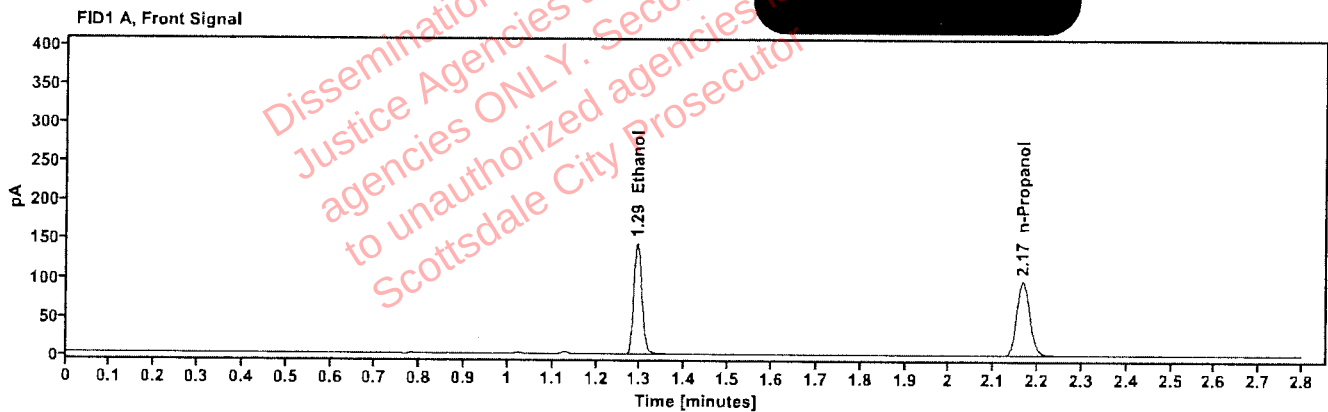


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2481           | 1.293      | 183.804   |
| n-Propanol         | —                | 2.166      | 199.884   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 190.004   |
| n-Propanol | 2.627      | 206.870   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618522-1                     | Vial:           | 36                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:32:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\36.D | Analyst:        | Abbott               |

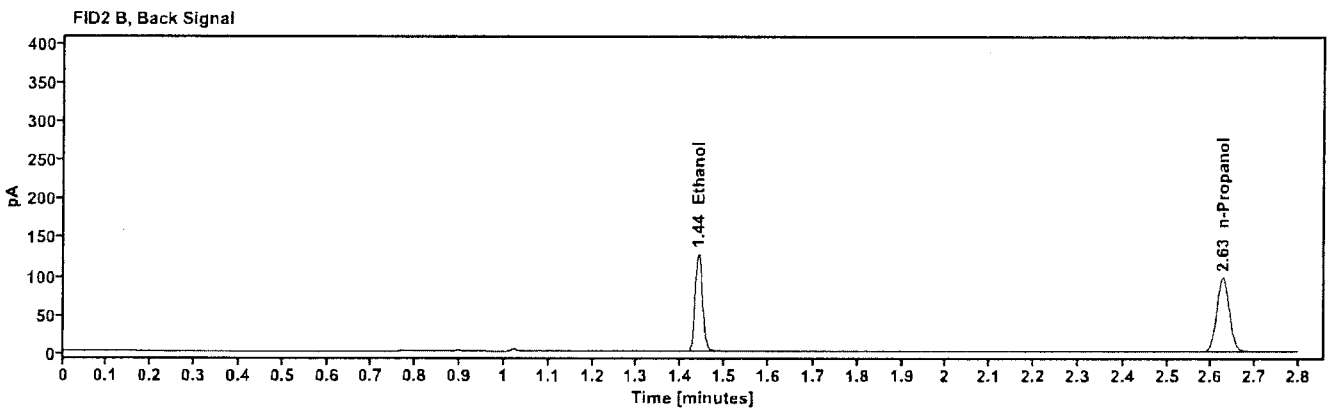
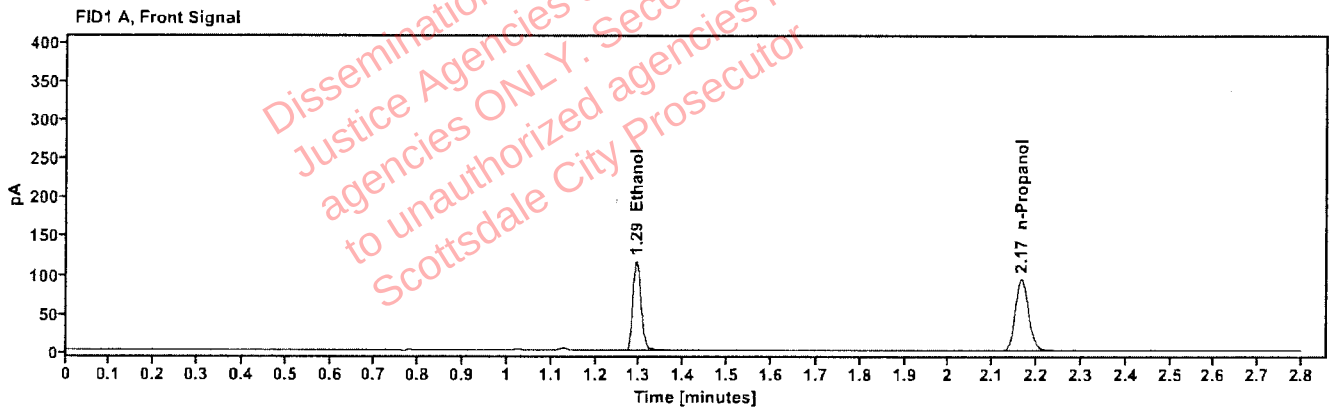


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2075           | 1.293      | 147.982   |
| n-Propanol         | -----            | 2.166      | 192.511   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 152.898   |
| n-Propanol | 2.627      | 199.380   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618522-1                      | Vial:           | 37                   |
| Description: | [REDACTED]                     | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 5:37:34 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\37.D | Analyst:        | Abbott               |

**VOID**

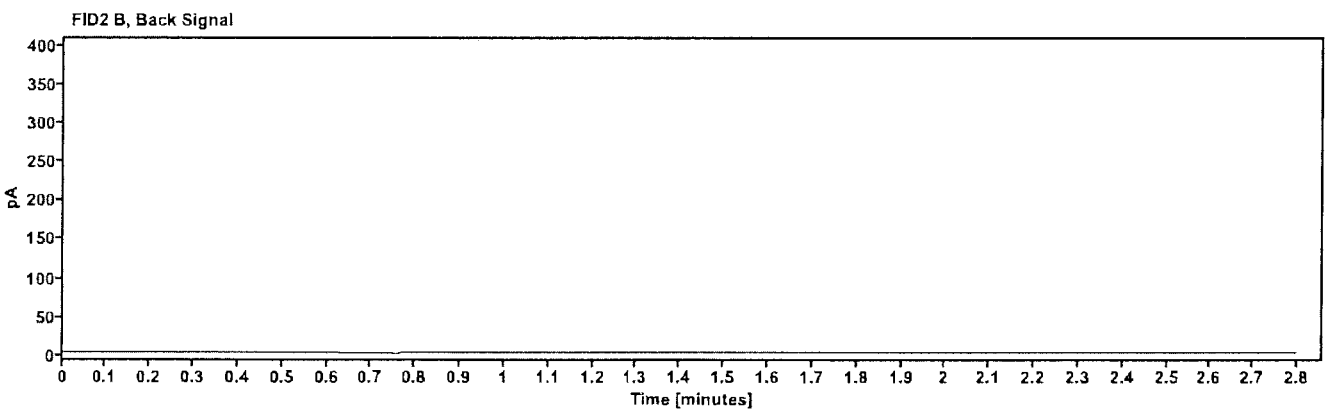
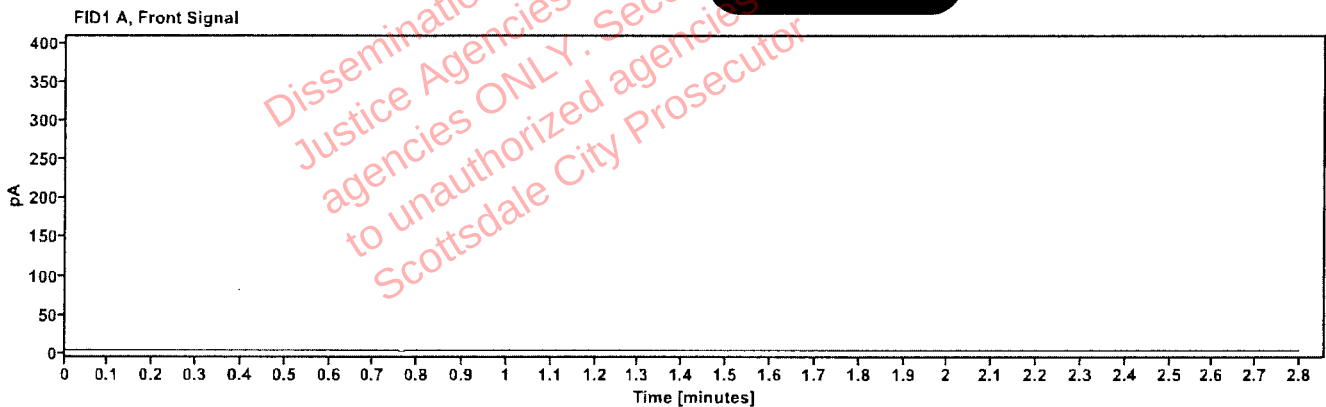


Table 1: FID 1 A (column DB-ALC1)

| Compound | Amount<br>(g/100mL) | Time<br>(min) | Peak<br>Area |
|----------|---------------------|---------------|--------------|
|----------|---------------------|---------------|--------------|

Table 2: FID 2 B (column DB-ALC2)

| Compound | Time<br>(min) | Peak<br>Area |
|----------|---------------|--------------|
|----------|---------------|--------------|

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618555-3                     | Vial:           | 38                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:41:23 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\38.D | Analyst:        | Abbott               |

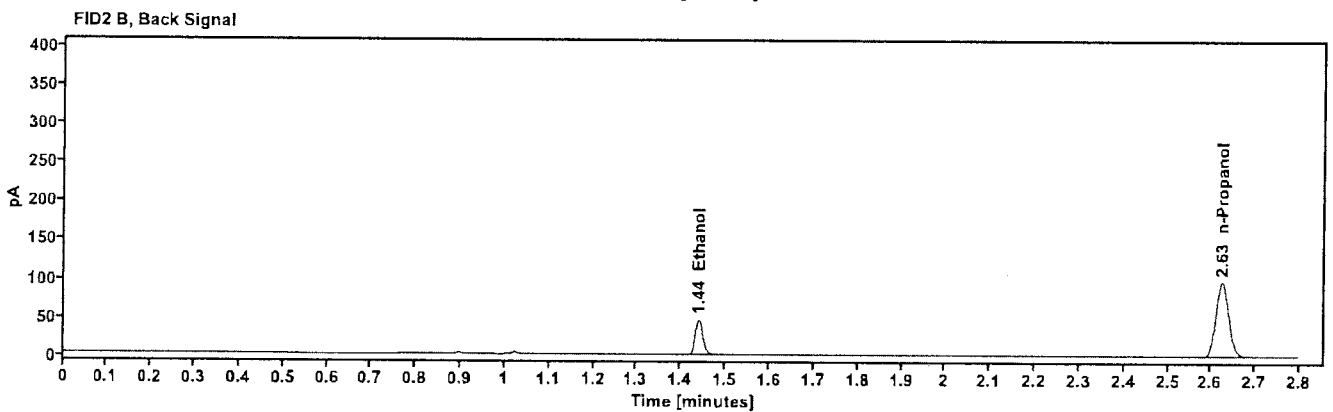
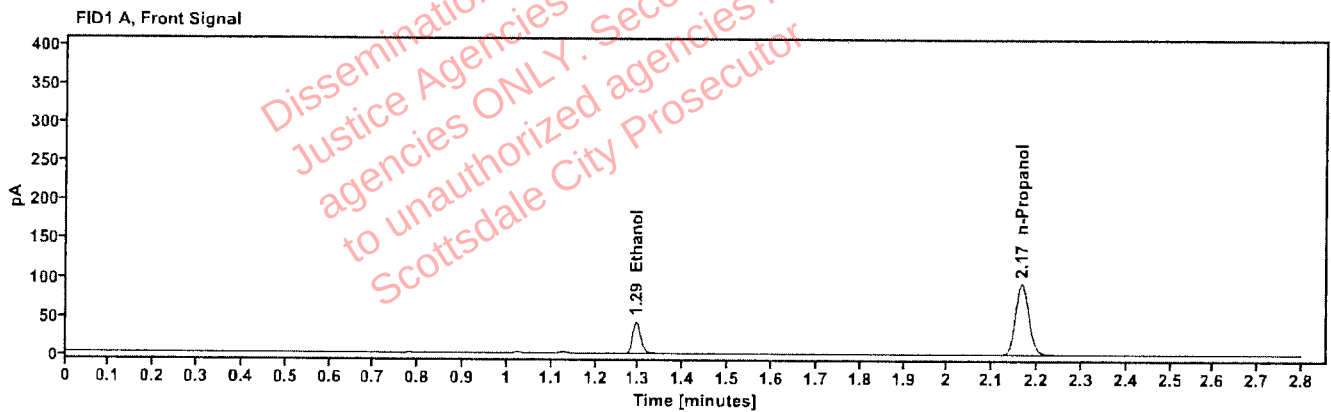


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0754           | 1.294      | 53.357    |
| n-Propanol         | ---              | 2.165      | 191.975   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 55.019    |
| n-Propanol | 2.626      | 198.965   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618555-3                     | Vial:           | 39                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:45:10 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\39.D | Analyst:        | Abbott               |

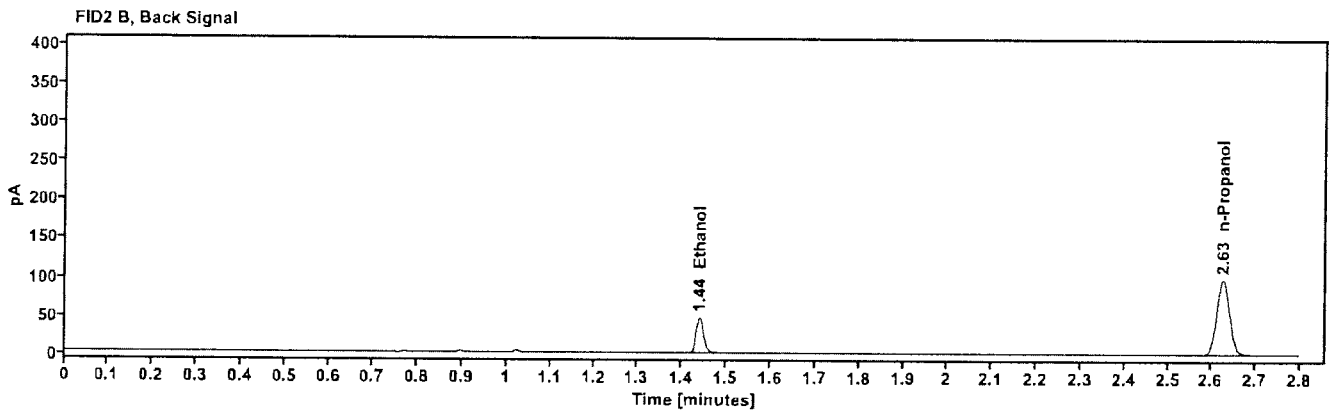
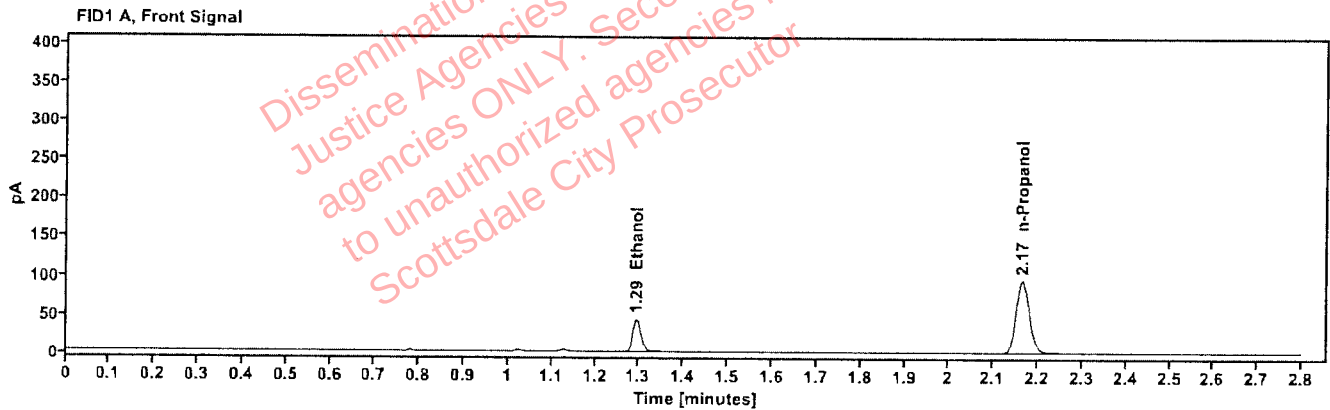


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0753           | 1.294      | 54.182    |
| n-Propanol         | ---              | 2.165      | 195.237   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 55.823    |
| n-Propanol | 2.626      | 202.244   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618554-1                     | Vial:           | 40                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 5:48:59 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\40.D | Analyst:        | Abbott               |

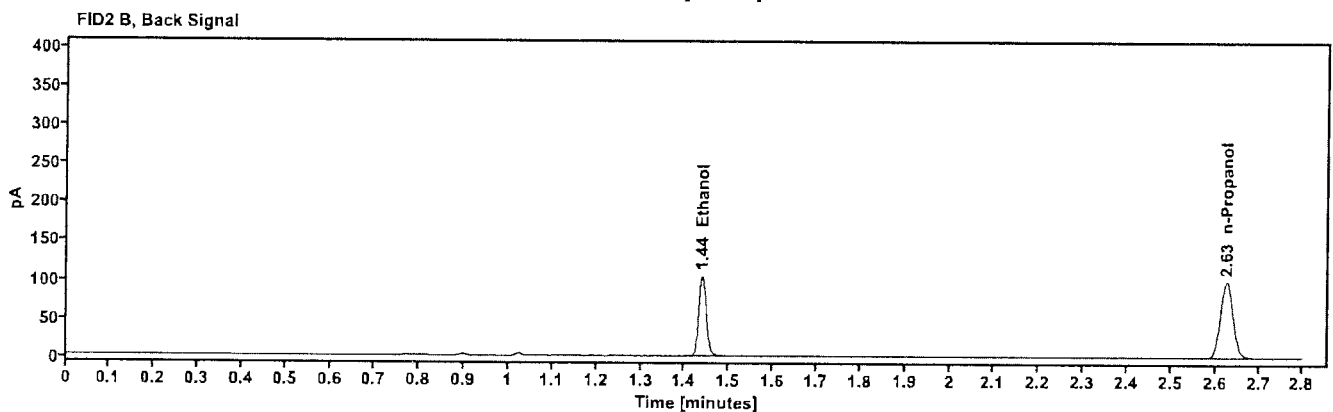
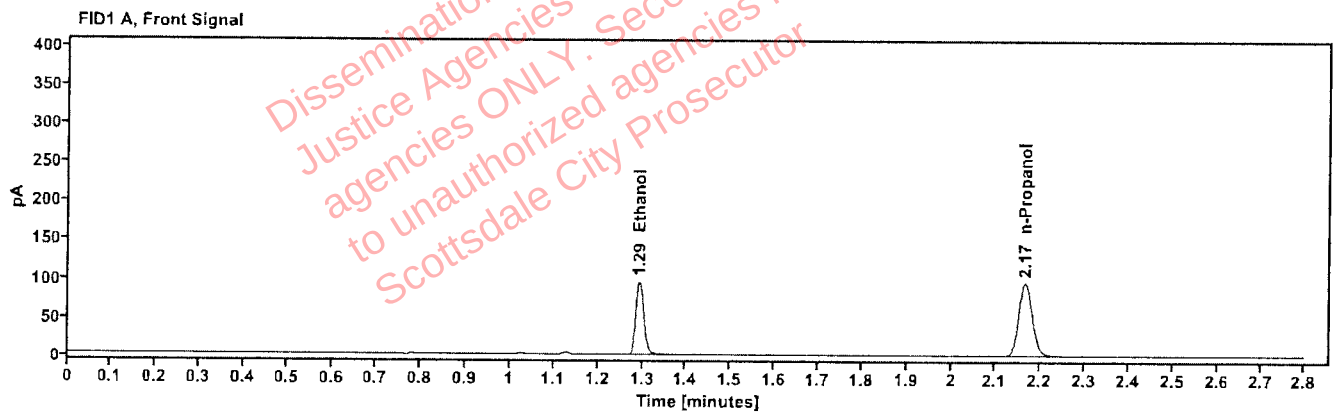


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1659           | 1.293      | 121.115   |
| n-Propanol         | -----            | 2.165      | 197.174   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 125.174   |
| n-Propanol | 2.626      | 203.991   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618554-1                      | Vial:           | 41                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 5:52:46 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\41.D | Analyst:        | Abbott               |

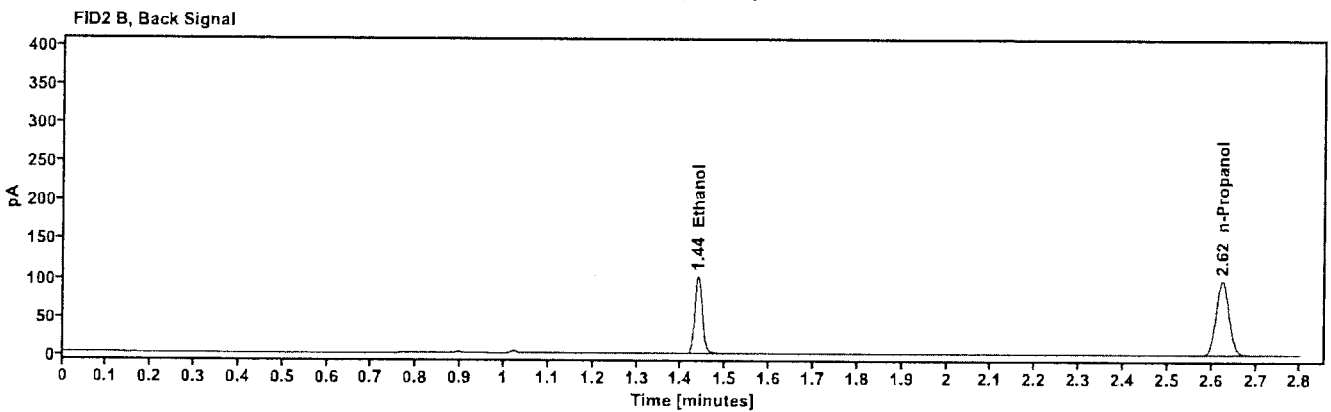
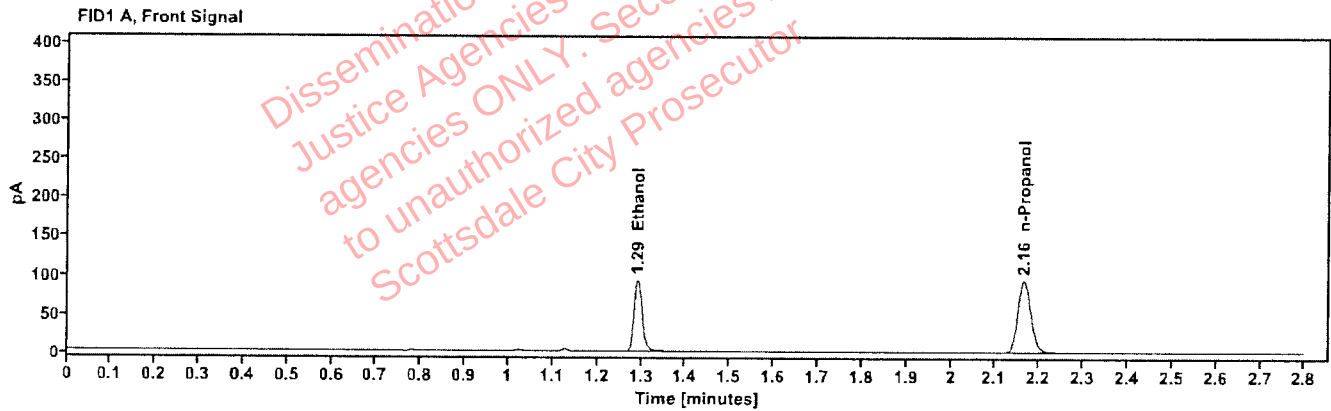


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1654           | 1.292      | 118.634   |
| n-Propanol         |                  | 2.164      | 193.773   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 122.278   |
| n-Propanol | 2.624      | 200.601   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618544-2                     | Vial:           | 43                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:01:35 PM |
| Data file:   | C:\Chem3211\Data\22Sep26\43.D | Analyst:        | Abbott               |

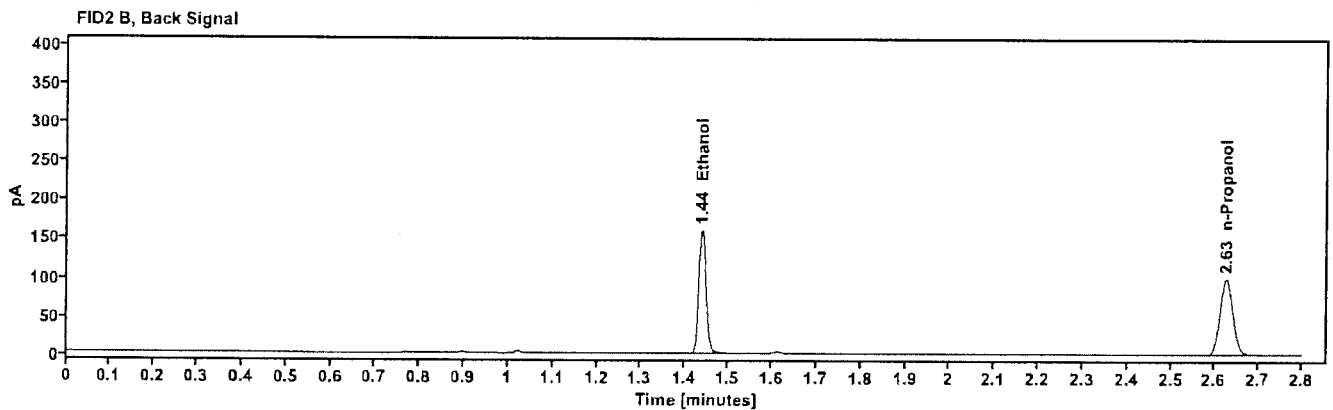
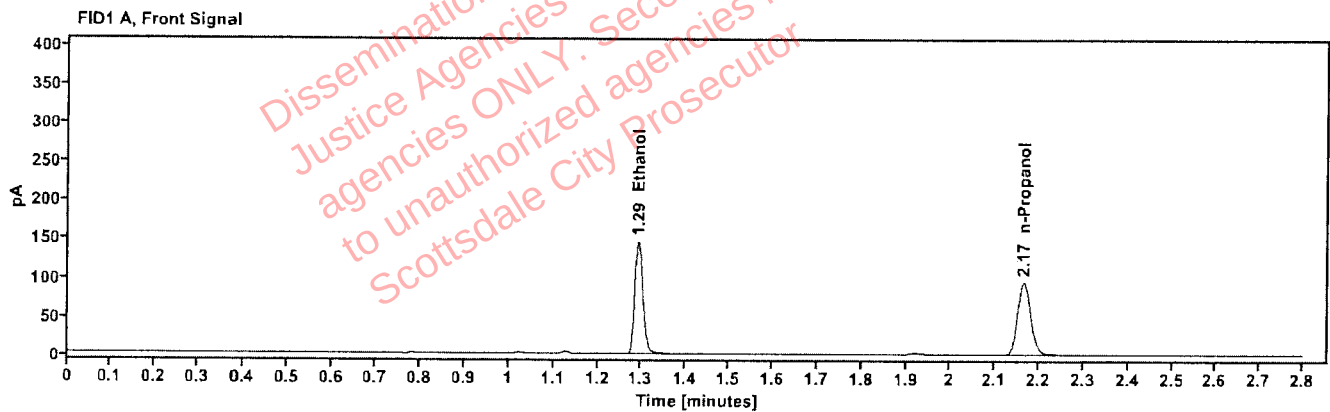


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2564           | 1.292      | 185.915   |
| n-Propanol         | ---              | 2.166      | 195.621   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 192.770   |
| n-Propanol | 2.626      | 202.551   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618544-2                      | Vial:           | 44                   |
| Description: | [REDACTED]                     | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 6:05:23 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\44.D | Analyst:        | Abbott               |

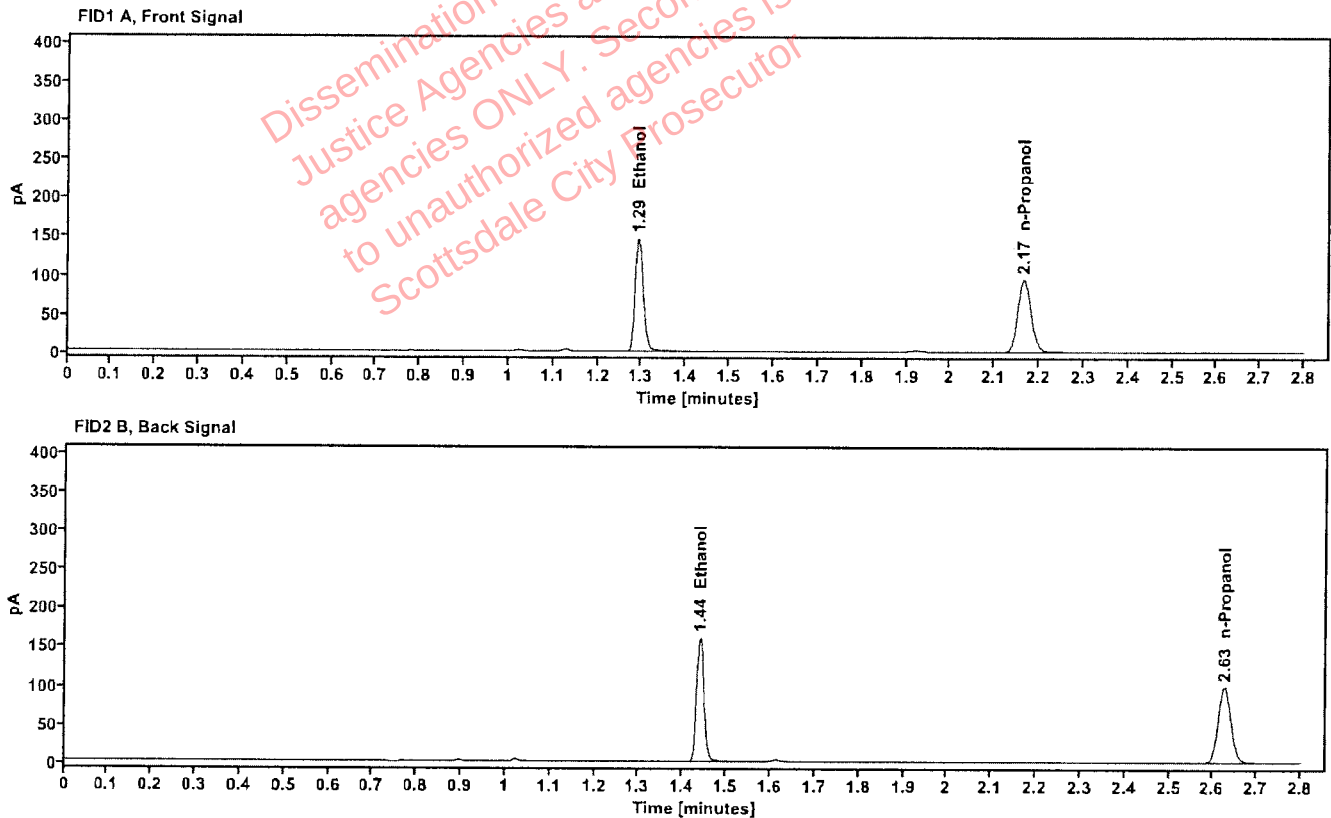


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2578           | 1.293      | 187.033   |
| n-Propanol         | -----            | 2.166      | 195.699   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 193.258   |
| n-Propanol | 2.626      | 202.462   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618053-2                     | Vial:           | 45                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:09:11 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\45.D | Analyst:        | Abbott               |

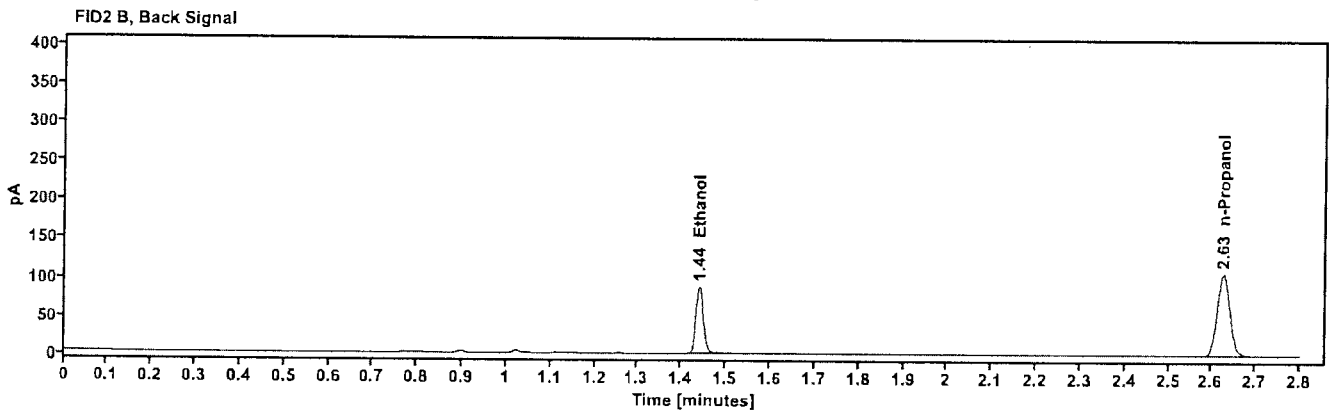
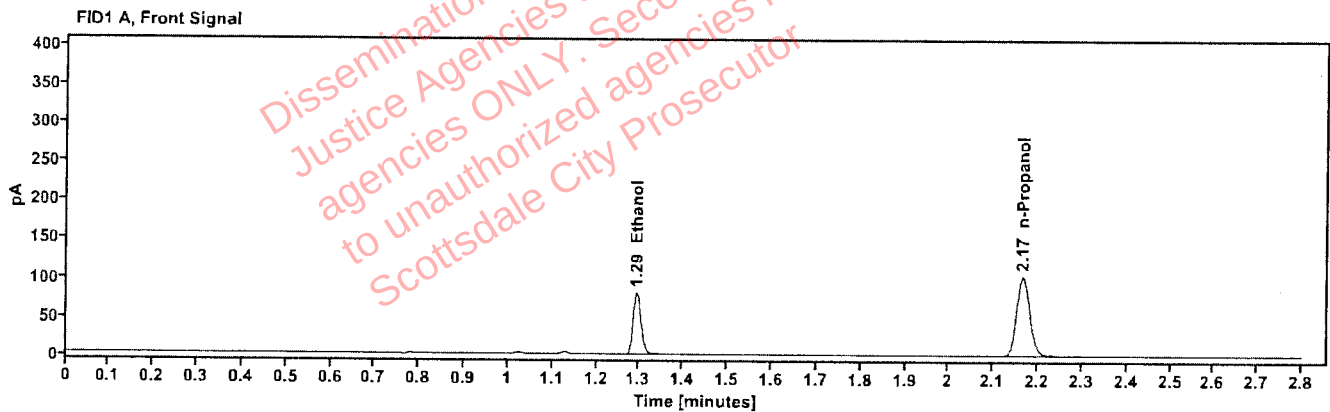


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1298           | 1.293      | 102.015   |
| n-Propanol         | -----            | 2.166      | 212.509   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 105.257   |
| n-Propanol | 2.627      | 220.089   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618053-2                     | Vial:           | 46                   |
| Description: | [REDACTED]                    | LIMS ID:        | [REDACTED]           |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:12:59 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\46.D | Analyst:        | Abbott               |

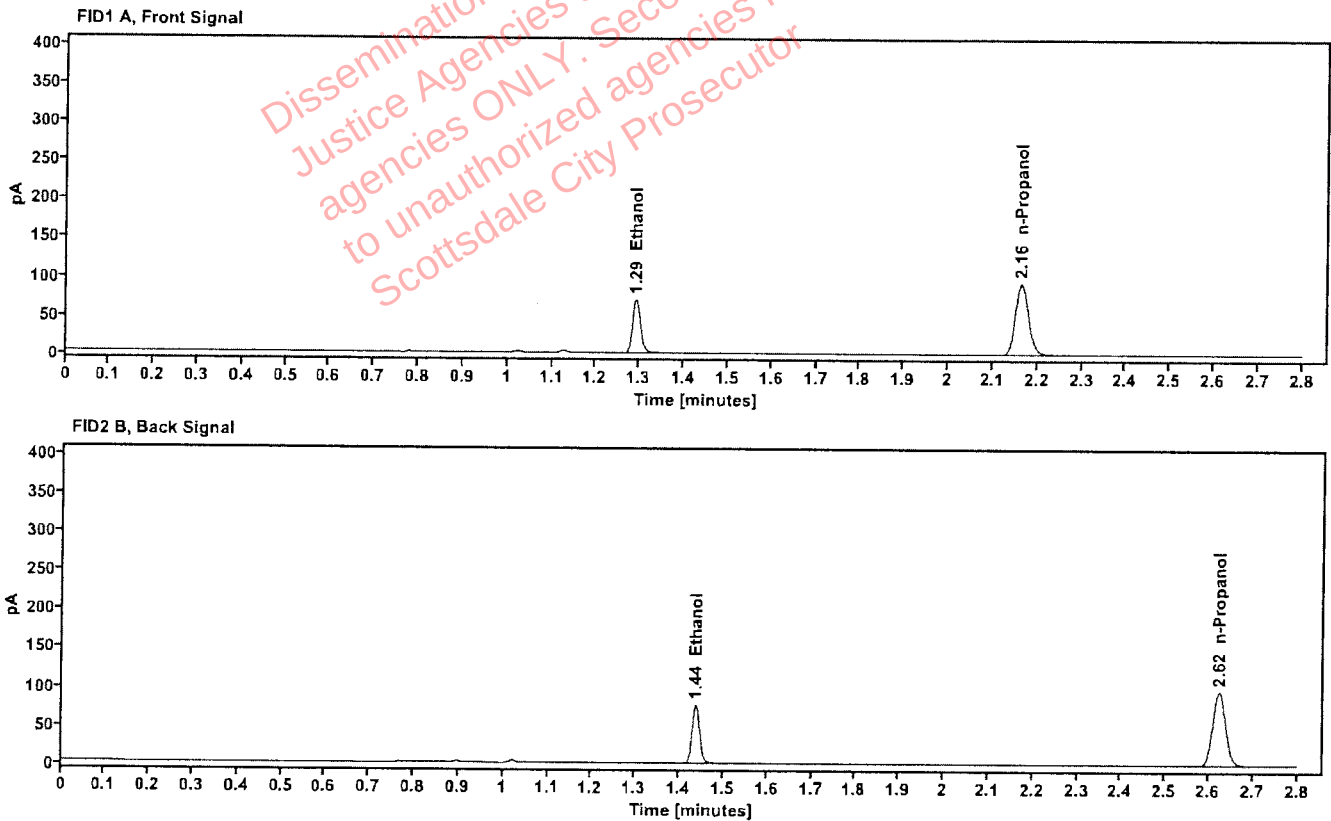


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1256           | 1.292      | 89.418    |
| n-Propanol         | -----            | 2.164      | 192.496   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 92.244    |
| n-Propanol | 2.624      | 199.290   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618598-1                      | Vial:           | 47                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 6:16:47 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\47.D | Analyst:        | Abbott               |

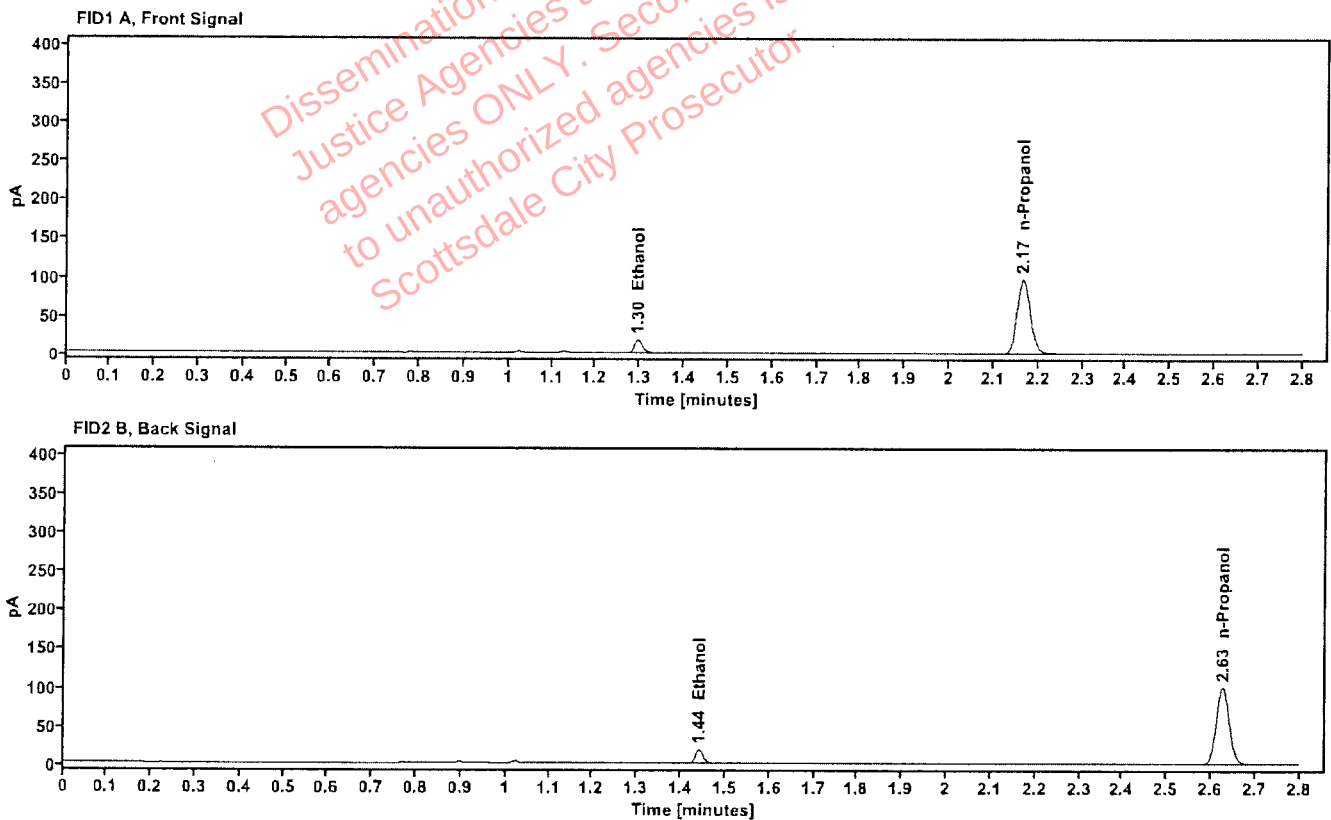


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0294           | 1.296      | 21.443    |
| n-Propanol         | -----            | 2.166      | 200.542   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.441      | 22.417    |
| n-Propanol | 2.627      | 208.007   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618598-1                     | Vial:           | 48                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:20:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\48.D | Analyst:        | Abbott               |

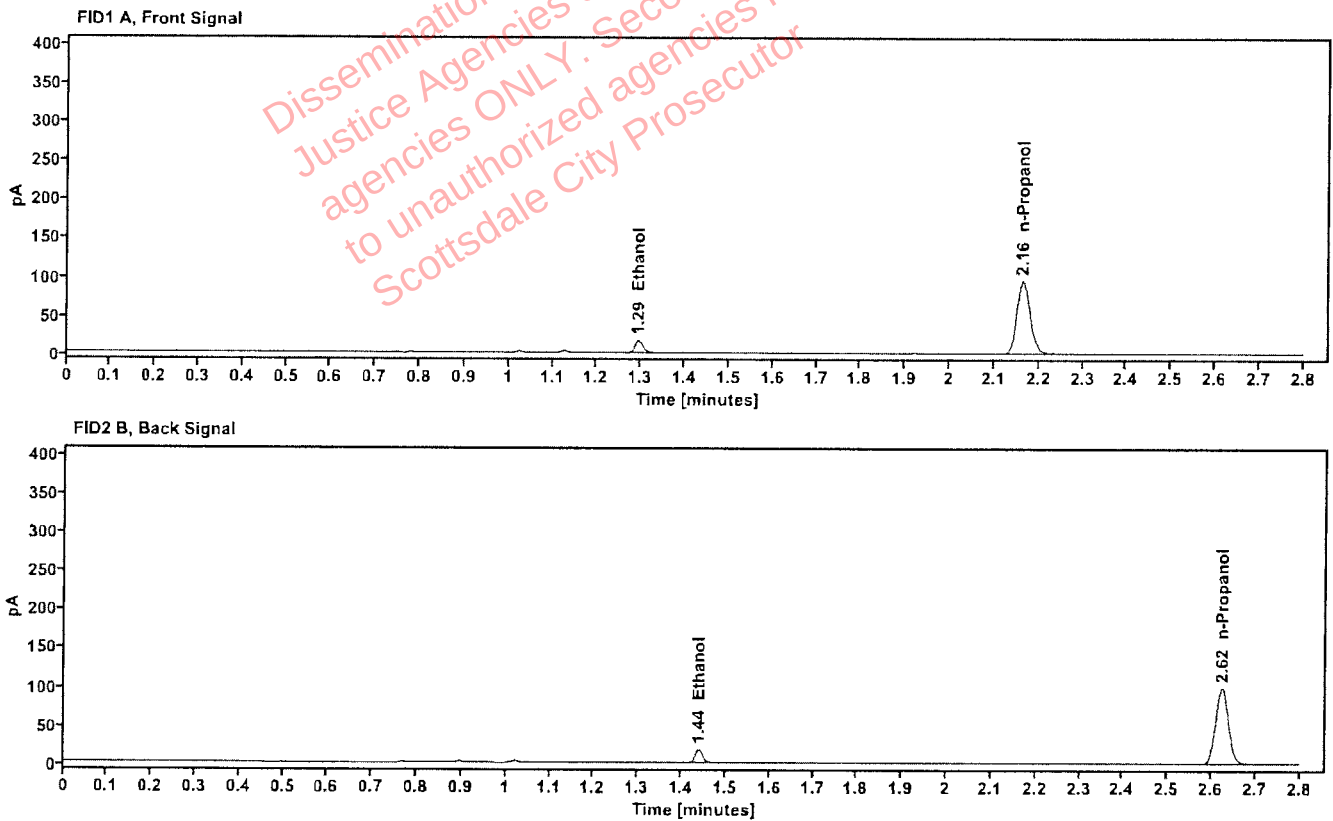


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0287           | 1.294      | 20.524    |
| n-Propanol         | ----             | 2.164      | 196.371   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 21.475    |
| n-Propanol | 2.624      | 203.682   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618621-2                     | Vial:           | 49                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:25:35 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\49.D | Analyst:        | Abbott               |

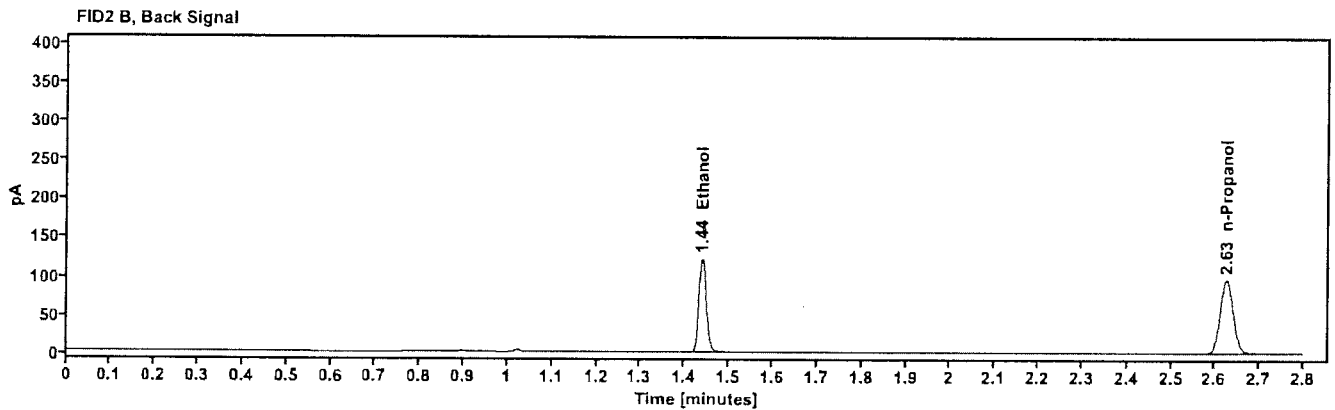
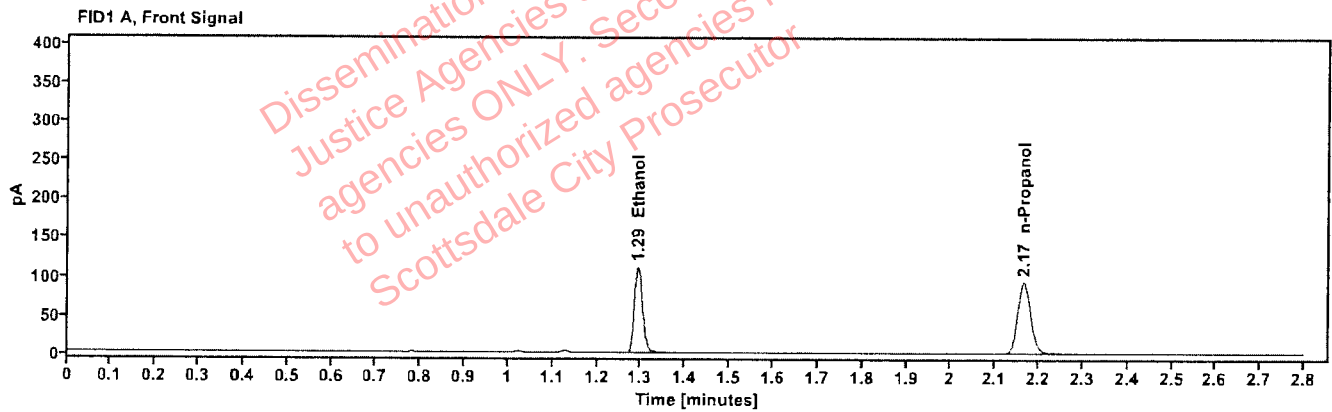


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1999           | 1.293      | 141.632   |
| n-Propanol         | ---              | 2.165      | 191.254   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 146.333   |
| n-Propanol | 2.626      | 197.801   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618621-2                     | Vial:           | 50                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:29:23 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\50.D | Analyst:        | Abbott               |

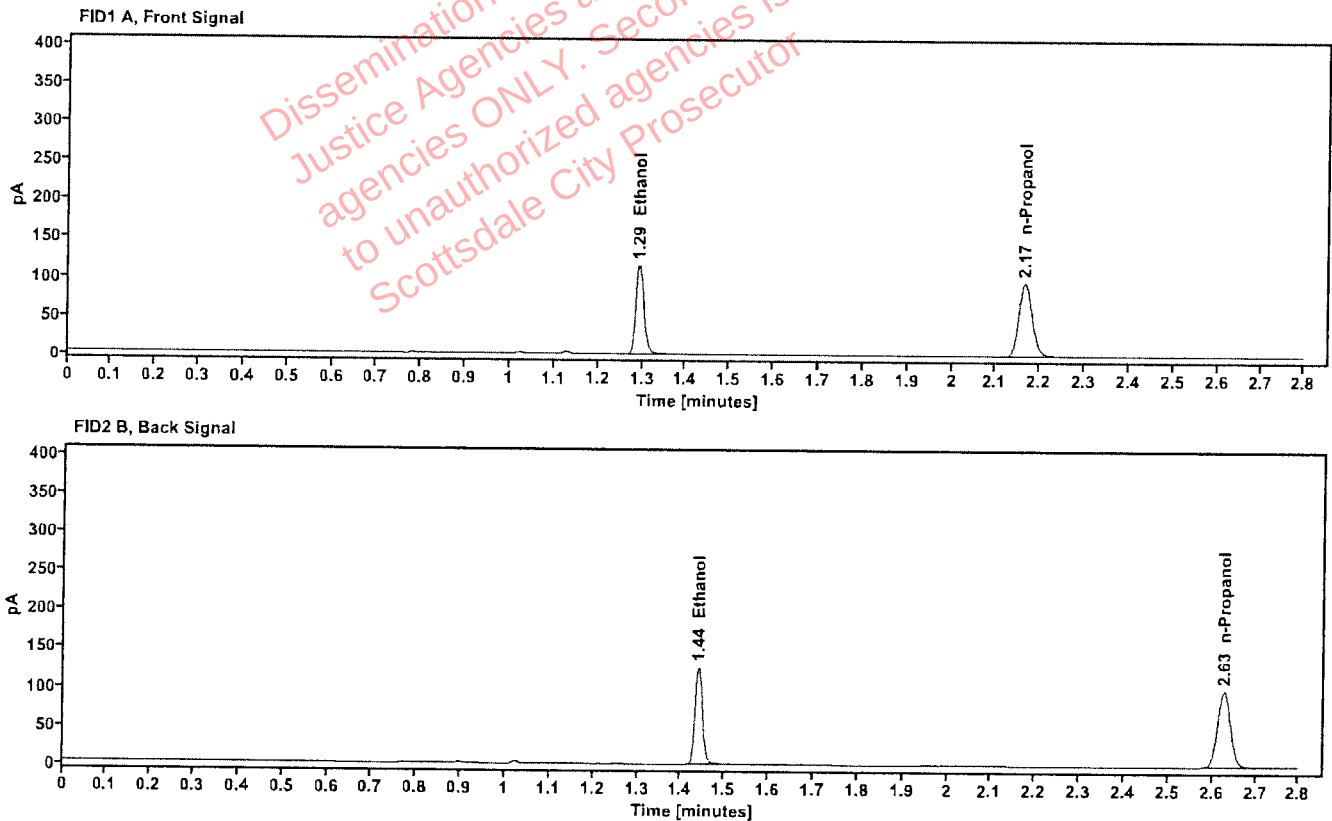


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2015           | 1.293      | 147,412   |
| n-Propanol         | —                | 2.166      | 197,510   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 152,272   |
| n-Propanol | 2.627      | 204,153   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618610-1                     | Vial:           | 51                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:33:11 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\51.D | Analyst:        | Abbott               |

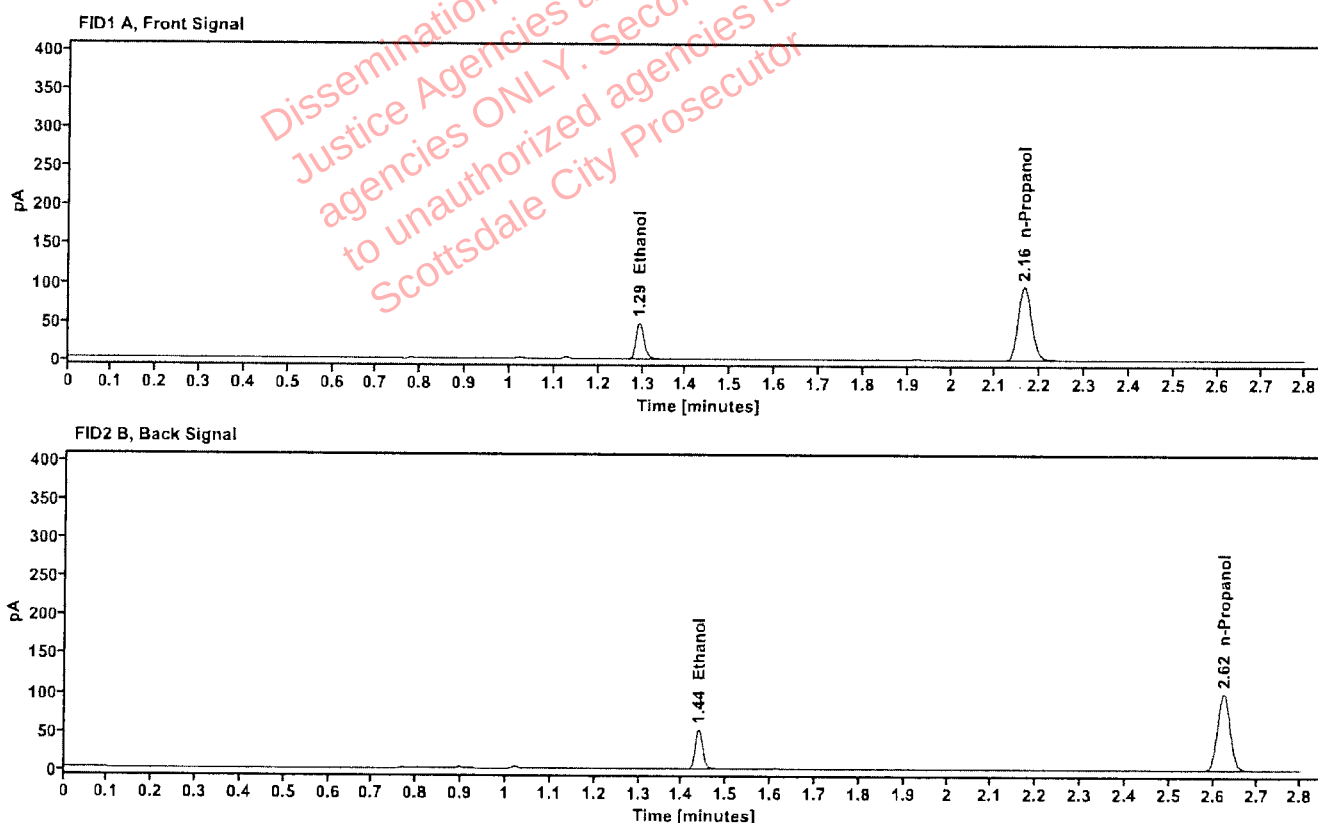


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0821           | 1.292      | 60.557    |
| n-Propanol         | ---              | 2.164      | 199.973   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.437      | 62.546    |
| n-Propanol | 2.624      | 207.260   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618610-1                     | Vial:           | 52                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:36:58 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\52.D | Analyst:        | Abbott               |

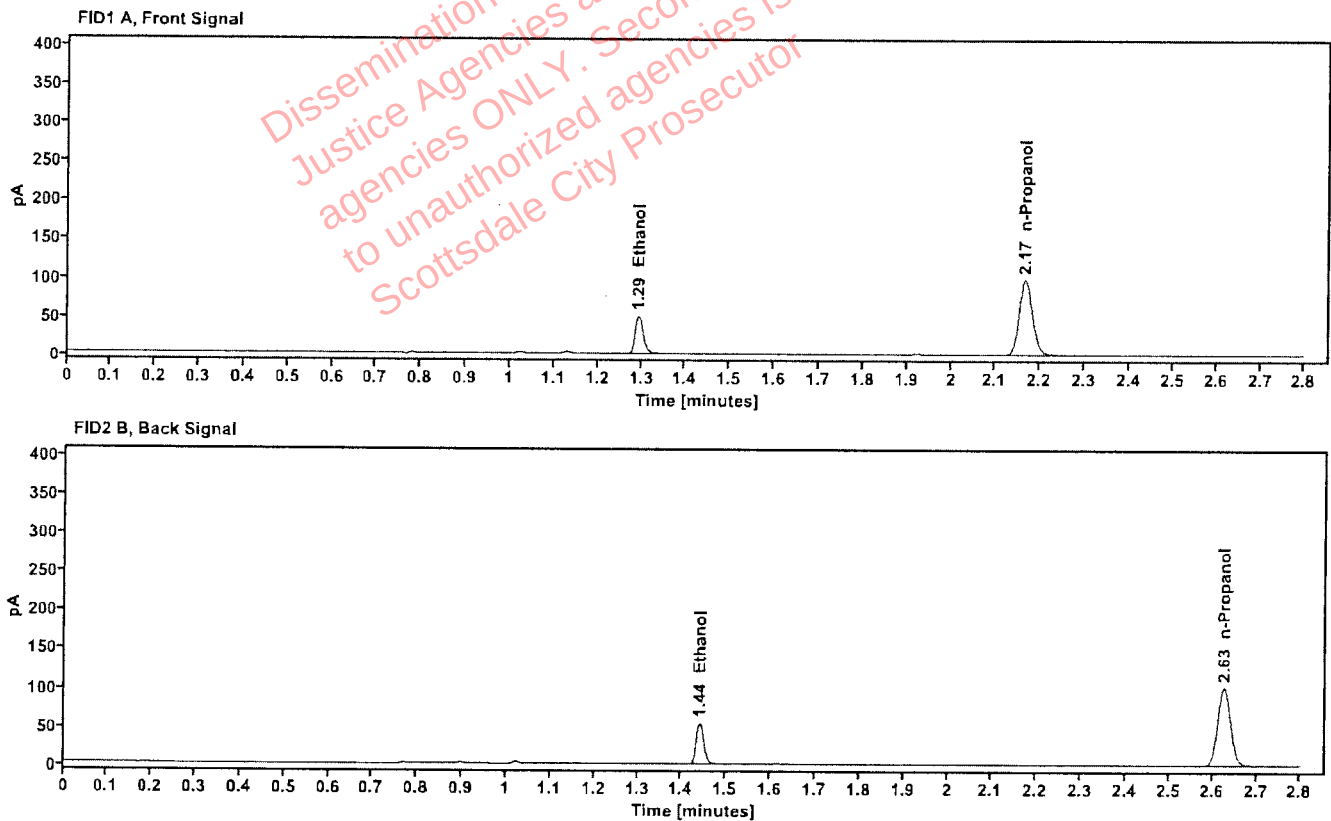


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0842           | 1.294      | 63.003    |
| n-Propanol         | —                | 2.166      | 202.729   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 64.981    |
| n-Propanol | 2.627      | 210.077   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618625-2                     | Vial:           | 54                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:44:34 PM |
| Data file:   | C:\Chem32\11Data\22Sep26\54.D | Analyst:        | Abbott               |

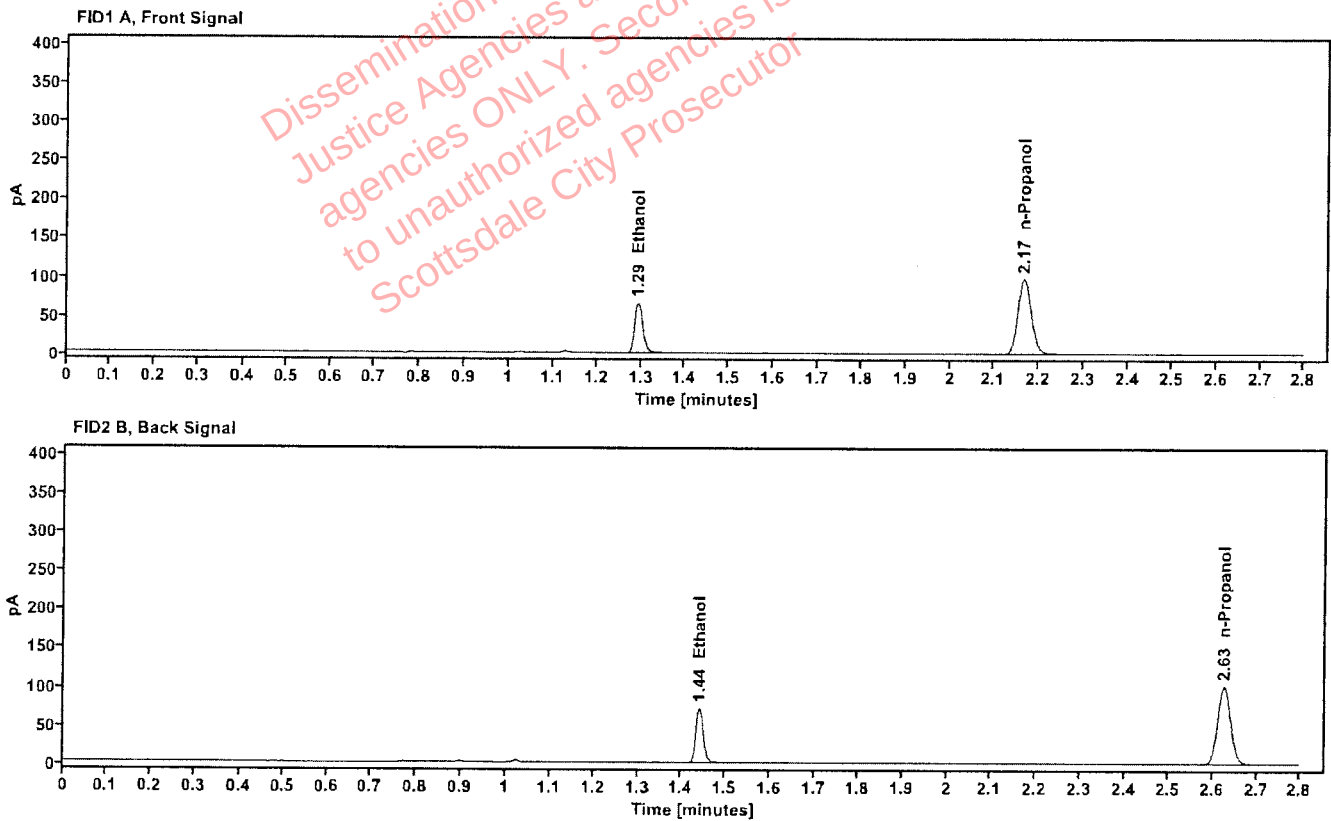


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1119           | 1.294      | 83.300    |
| n-Propanol         | ---              | 2.166      | 201.523   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 85.888    |
| n-Propanol | 2.627      | 208.617   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618625-2                      | Vial:           | 55                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 6:49:34 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\55.D | Analyst:        | Abbott               |

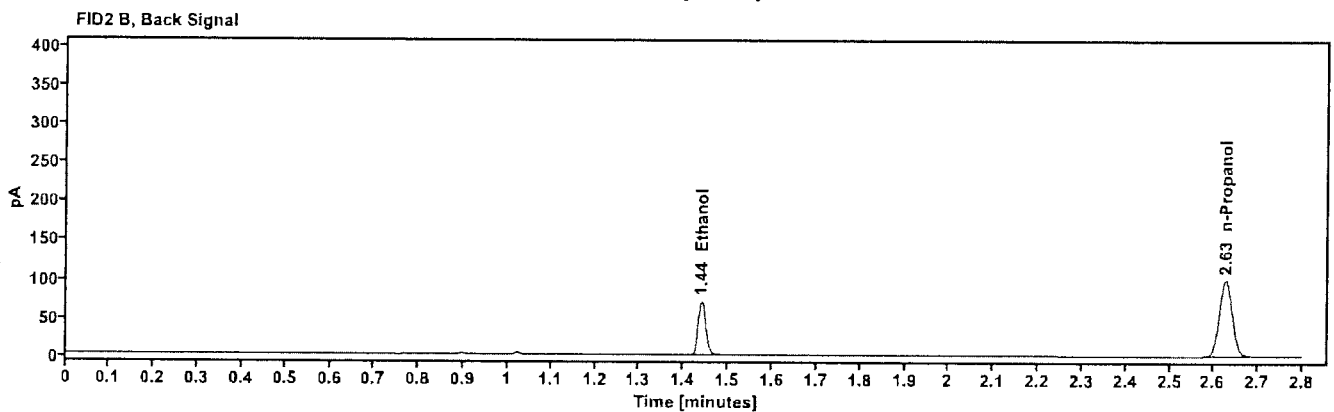
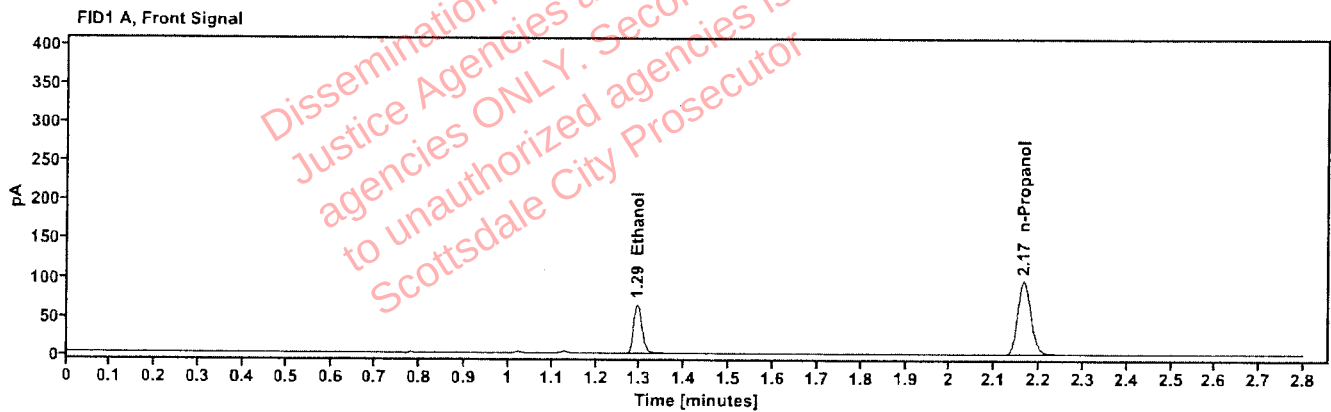


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.1116           | 1.294      | 81.713    |
| n-Propanol         | ---              | 2.166      | 198.116   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 84.309    |
| n-Propanol | 2.627      | 205.136   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618602-2                     | Vial:           | 56                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 6:53:22 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\56.D | Analyst:        | Abbott               |

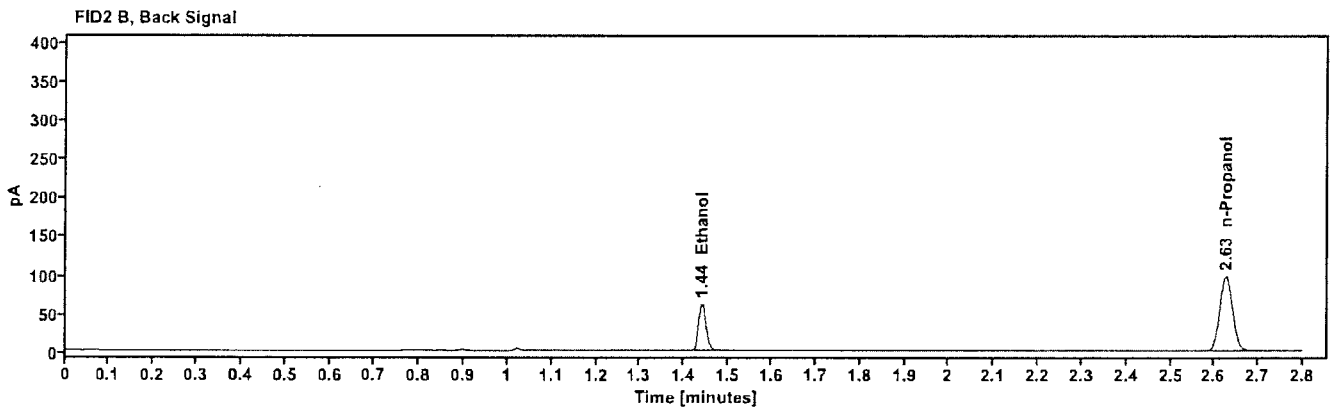
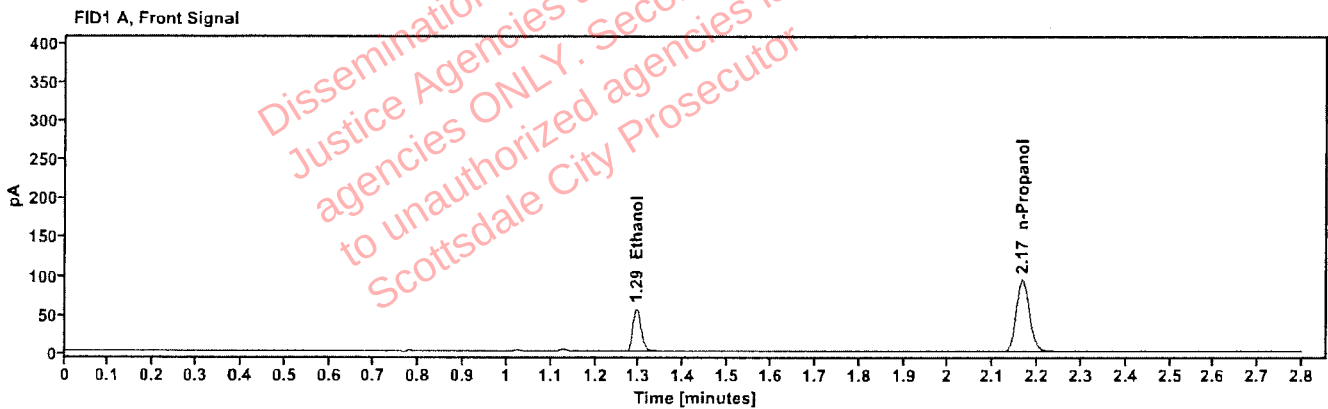


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0996           | 1.294      | 71.217    |
| n-Propanol         | -----            | 2.166      | 193.548   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 73.457    |
| n-Propanol | 2.627      | 200.508   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2618602-2                      | Vial:           | 57                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 6:57:11 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\57.D | Analyst:        | Abbott               |

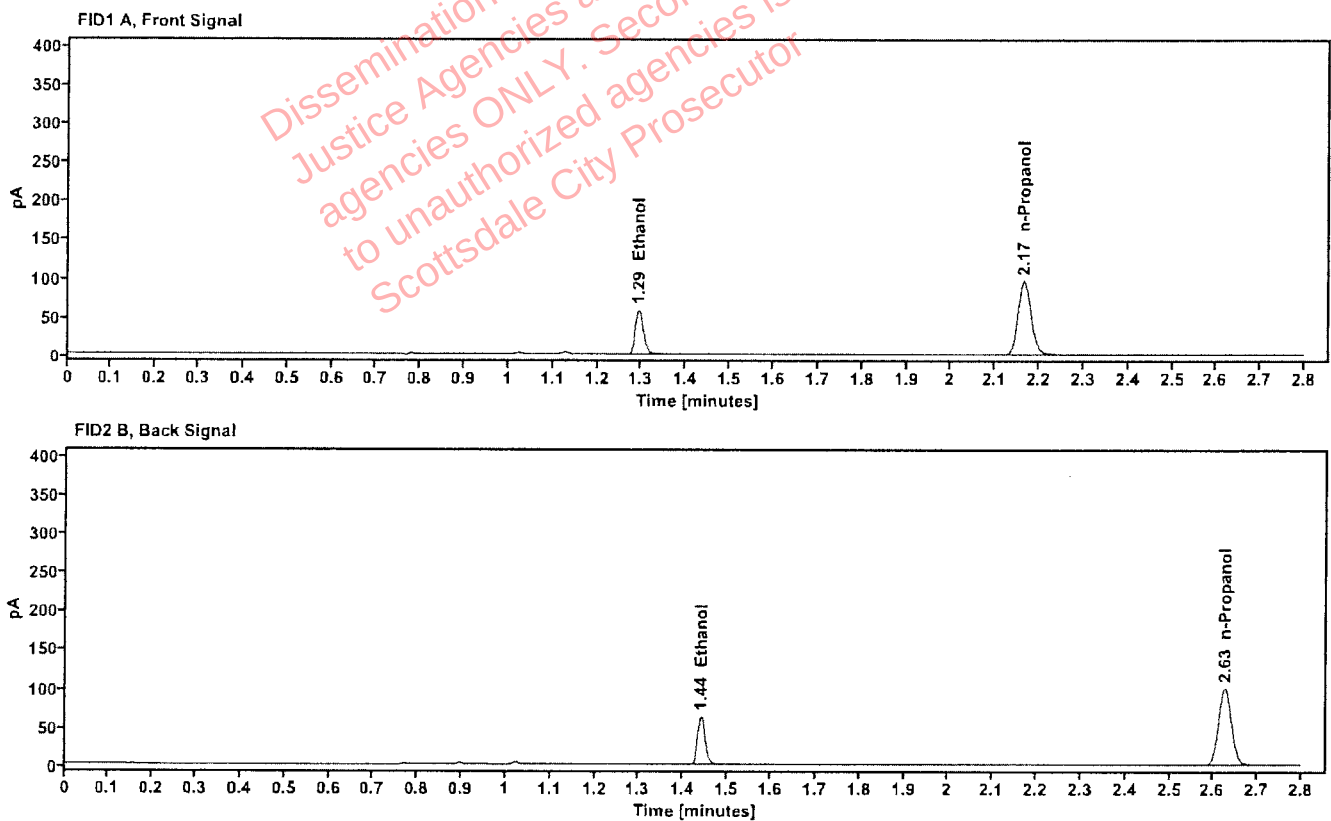


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.0997           | 1.294      | 72.959    |
| n-Propanol         | -----            | 2.166      | 198.213   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.440      | 75.265    |
| n-Propanol | 2.627      | 205.306   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2609690-2                     | Vial:           | 58                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:00:59 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\58.D | Analyst:        | Abbott               |

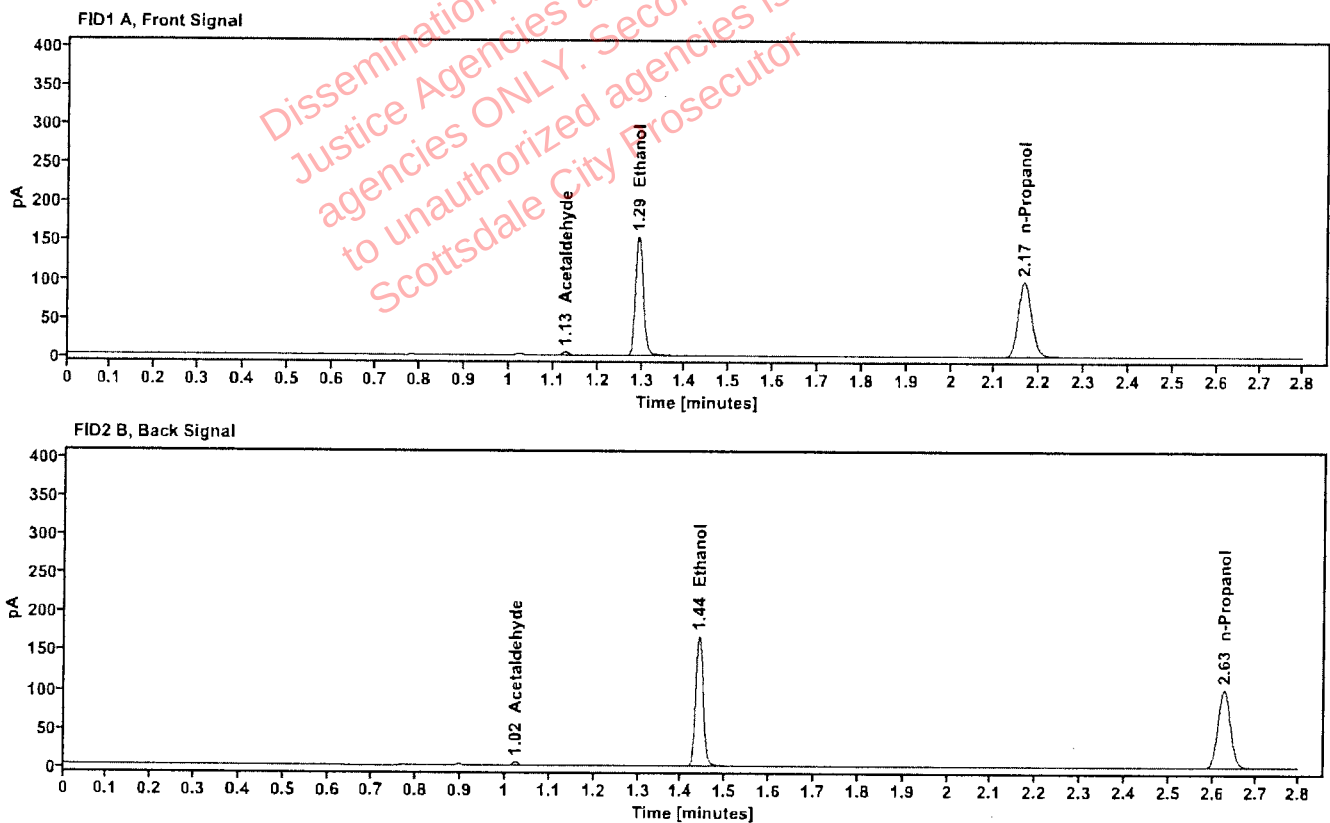


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| Acetaldehyde       | -----            | 1.127      | 4.370     |
| quantified Ethanol | 0.2621           | 1.293      | 196.351   |
| n-Propanol         | -----            | 2.166      | 202.145   |

Table 2: FID 2 B (column DB-ALC2)

| Compound     | Time (min) | Peak Area |
|--------------|------------|-----------|
| Acetaldehyde | 1.023      | 4.533     |
| Ethanol      | 1.440      | 203.379   |
| n-Propanol   | 2.627      | 209.429   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2609690-2                     | Vial:           | 59                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:04:47 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\59.D | Analyst:        | Abbott               |

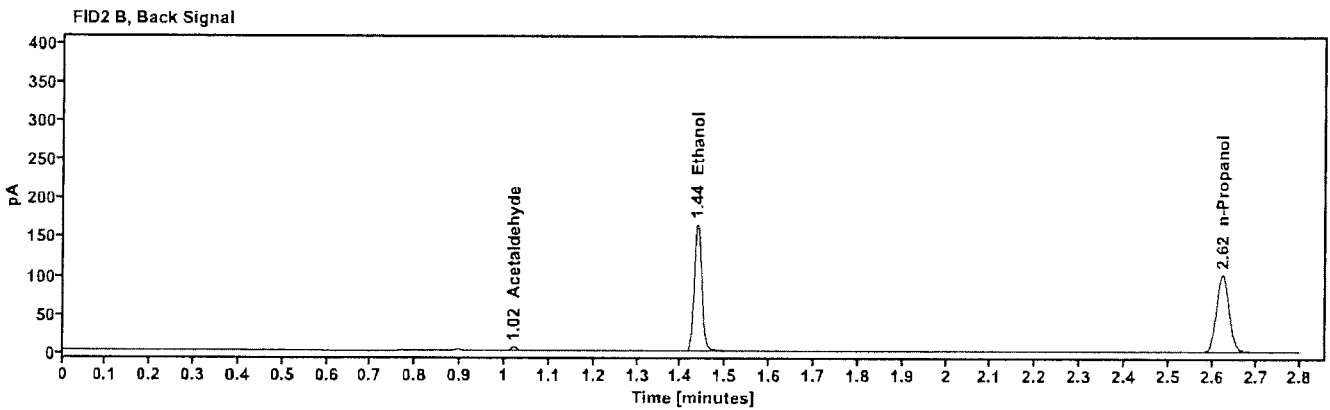
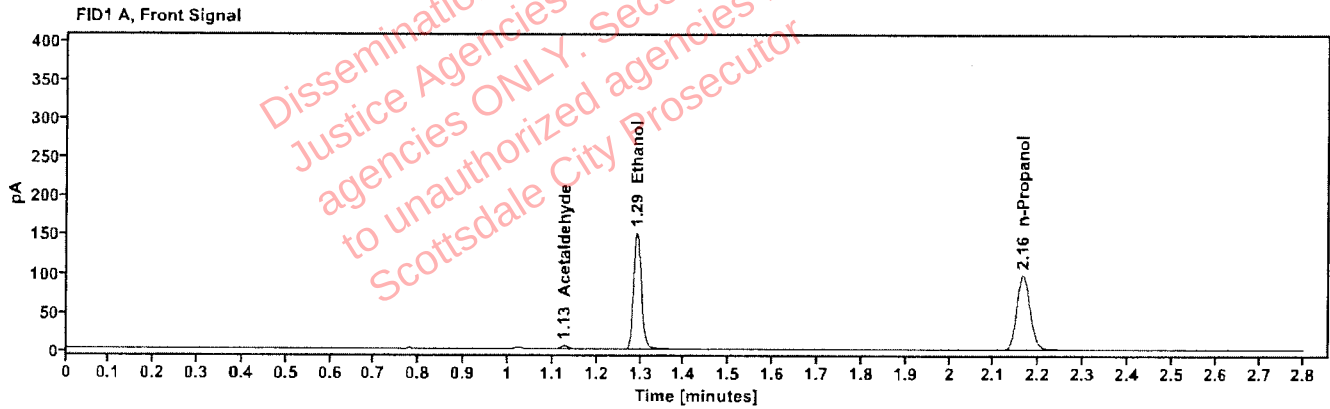


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| Acetaldehyde       | ---              | 1.126      | 4.425     |
| quantified Ethanol | 0.2598           | 1.291      | 193.686   |
| n-Propanol         | ---              | 2.164      | 201.148   |

Table 2: FID 2 B (column DB-ALC2)

| Compound     | Time (min) | Peak Area |
|--------------|------------|-----------|
| Acetaldehyde | 1.021      | 4.672     |
| Ethanol      | 1.437      | 200.602   |
| n-Propanol   | 2.624      | 208.361   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                                |                 |                      |
|--------------|--------------------------------|-----------------|----------------------|
| Sample:      | 2605726-2                      | Vial:           | 60                   |
| Description: |                                | LIMS ID:        |                      |
| Method:      | ethanol quant.M                | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007          | Injection date: | 9/22/2026 7:08:35 PM |
| Data file:   | C:\Chem32\11\Data\22Sep26\60.D | Analyst:        | Abbott               |

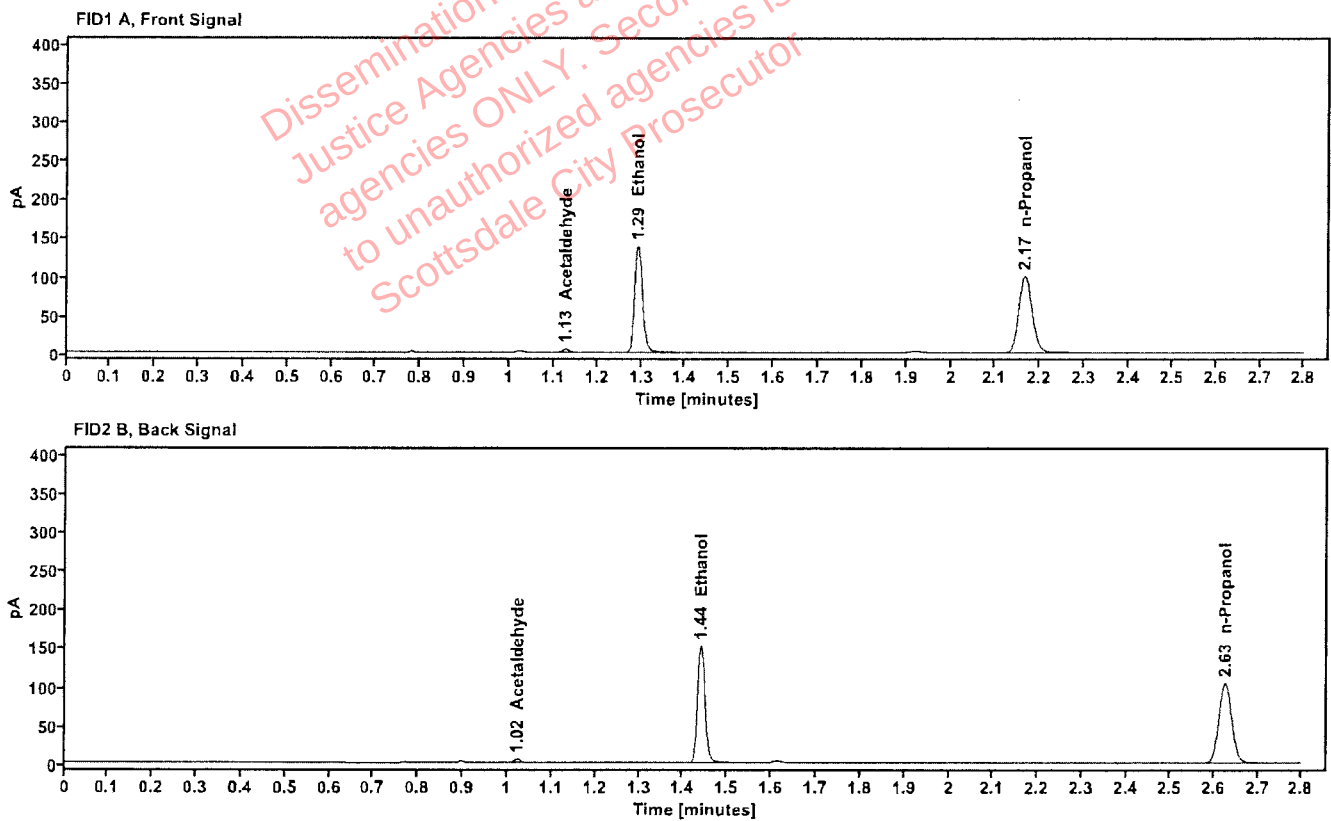


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| Acetaldehyde       | ---              | 1.127      | 4.435     |
| quantified Ethanol | 0.2311           | 1.293      | 177.192   |
| n-Propanol         | ---              | 2.166      | 206.879   |

Table 2: FID 2 B (column DB-ALC2)

| Compound     | Time (min) | Peak Area |
|--------------|------------|-----------|
| Acetaldehyde | 1.023      | 4.588     |
| Ethanol      | 1.440      | 183.067   |
| n-Propanol   | 2.627      | 214.244   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2605726-2                     | Vial:           | 61                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:13:34 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\61.D | Analyst:        | Abbott               |

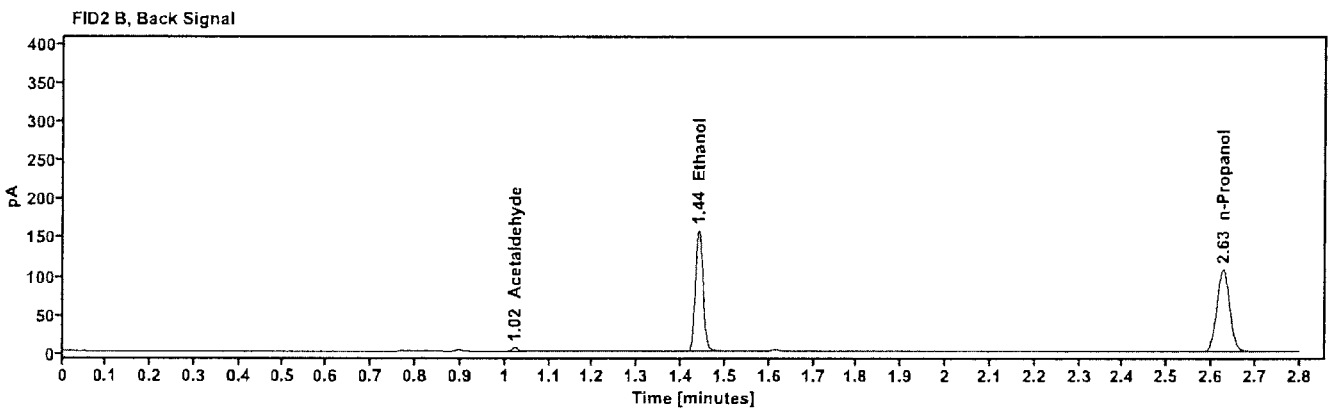
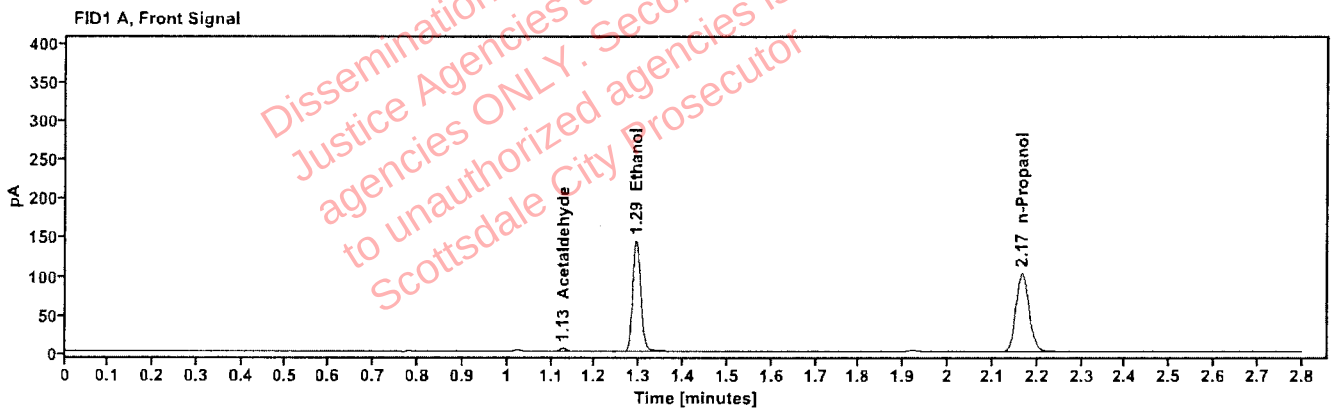


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| Acetaldehyde       | -----            | 1.127      | 4.338     |
| quantified Ethanol | 0.2331           | 1.293      | 183.451   |
| n-Propanol         | -----            | 2.166      | 212.401   |

Table 2: FID 2 B (column DB-ALC2)

| Compound     | Time (min) | Peak Area |
|--------------|------------|-----------|
| Acetaldehyde | 1.023      | 4.507     |
| Ethanol      | 1.440      | 189.587   |
| n-Propanol   | 2.627      | 219.938   |

# Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618613-2                     | Vial:           | 62                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:17:23 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\62.D | Analyst:        | Abbott               |

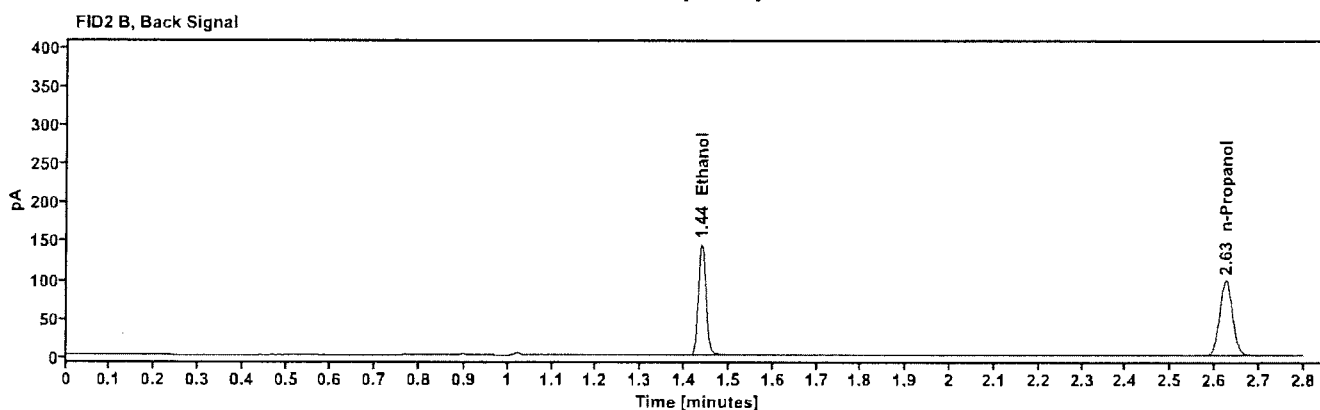
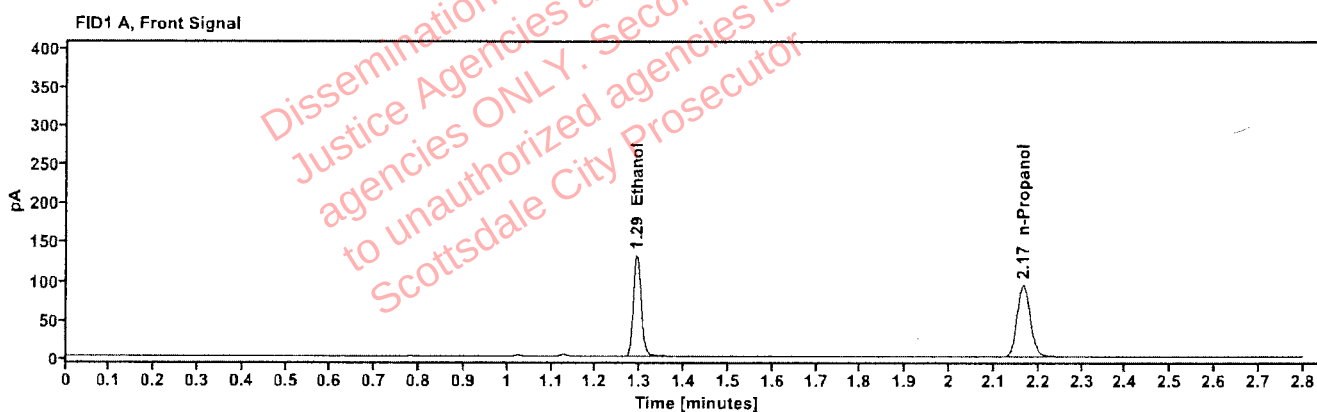


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2333           | 1.292      | 167.655   |
| n-Propanol         | -----            | 2.165      | 193.888   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 173.694   |
| n-Propanol | 2.626      | 200.606   |

## Scottsdale Police Department Crime Lab Volatiles Analysis

|              |                               |                 |                      |
|--------------|-------------------------------|-----------------|----------------------|
| Sample:      | 2618613-2                     | Vial:           | 63                   |
| Description: |                               | LIMS ID:        |                      |
| Method:      | ethanol quant.M               | Sequence:       | 22Sep26              |
| Instrument:  | US2320A024 US23180007         | Injection date: | 9/22/2026 7:21:10 PM |
| Data file:   | C:\Chem32\1\Data\22Sep26\63.D | Analyst:        | Abbott               |

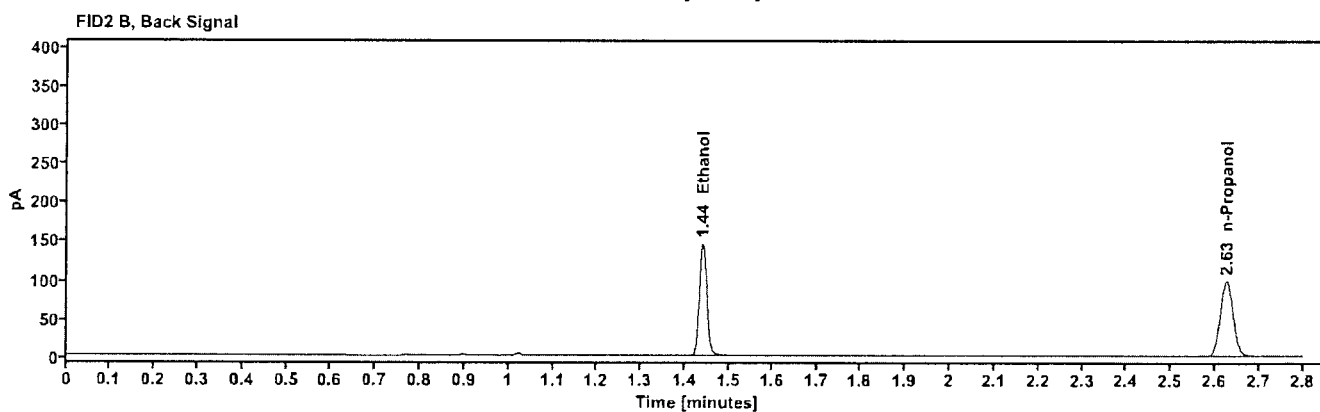
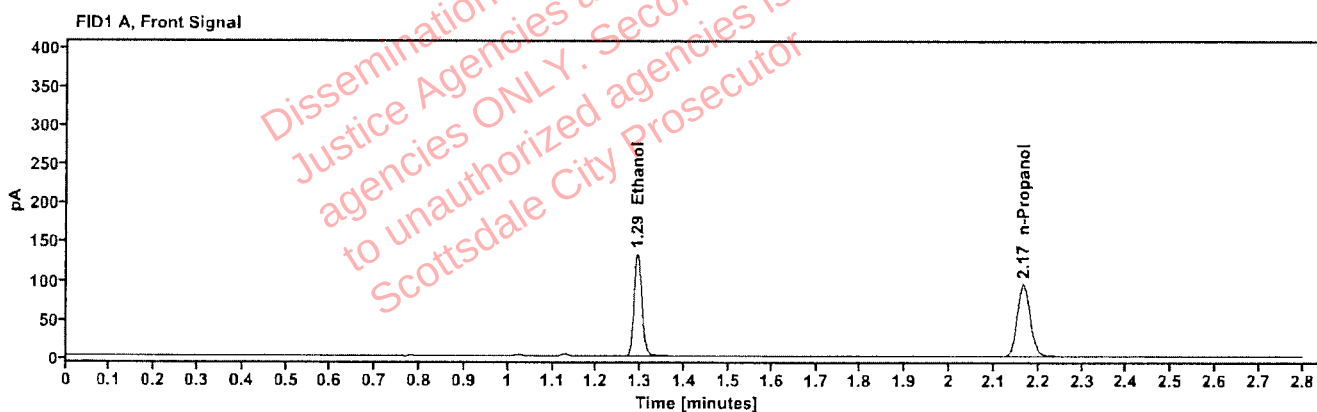


Table 1: FID 1 A (column DB-ALC1)

| Compound           | Amount (g/100mL) | Time (min) | Peak Area |
|--------------------|------------------|------------|-----------|
| quantified Ethanol | 0.2339           | 1.292      | 168.967   |
| n-Propanol         | ---              | 2.165      | 194.906   |

Table 2: FID 2 B (column DB-ALC2)

| Compound   | Time (min) | Peak Area |
|------------|------------|-----------|
| Ethanol    | 1.439      | 174.739   |
| n-Propanol | 2.626      | 201.626   |